

# RAYSTATION 10A

DICOM Conformance Statement



RayStation

10A

### *Declaration of conformity*



Complies with 93/42/EEC Medical Device Directive as amended by M1 to M5. A copy of the corresponding Declaration of Conformity is available on request.

### *Copyright*

This document contains proprietary information that is protected by copyright. No part of this document may be photocopied, reproduced or translated to another language without prior written consent of RaySearch Laboratories AB (publ).

All Rights Reserved. ©2020 RaySearch Laboratories AB (publ).

### *Trademarks*

RayStation, RayBiology, RayCare, RayStore, RayCloud, RayPlan, RaySearch Laboratories and the RaySearch Laboratories logotype are trademarks of RaySearch Laboratories AB (publ).

Third party trademarks as used herein are the property of their respective owners, which are not affiliated with RaySearch Laboratories AB (publ).

RaySearch Laboratories AB (publ) including its subsidiaries is hereafter referred to as RaySearch.

# TABLE OF CONTENTS

TABLE OF CONTENTS	1
OVERVIEW	2
INTRODUCTION	3
REVISION HISTORY	3
REMARKS	3
NETWORKING AND MEDIA INTERCHANGE	4
IMPLEMENTATION MODEL	4
TRANSFER SYNTAX	4
IMPORTING DATA	4
EXPORTING DATA	5
SECURITY	7
SUPPORT OF CHARACTER SETS	8
KNOWN ISSUES	9
IMPORT OF RESPIRATORY PHASE GATING DUTY CYCLE	9
ANNEXES	10
IMPORT IOD DEFINITIONS	10
EXPORT IOD DEFINITIONS	50
DATA DICTIONARY OF PRIVATE ATTRIBUTES	90

# 1 OVERVIEW

The following DICOM services for receiving and sending imaging and radiotherapy data are supported.

SOP Class Name	SOP Class UID	Provider of Service (SCP)	User of Service (SCU)
Transfer and Query/Retrieve			
CT Image Storage	1.2.840.10008.5.1.4.1.1.2	Yes	Yes
MR Image Storage	1.2.840.10008.5.1.4.1.1.4	Yes	Yes
PET Image Storage	1.2.840.10008.5.1.4.1.1.128	Yes	Yes
Spatial Registration Storage	1.2.840.10008.5.1.4.1.1.66.1	Yes	Yes
RT Structure Set Storage	1.2.840.10008.5.1.4.1.1.481.3	Yes	Yes
RT Plan Storage	1.2.840.10008.5.1.4.1.1.481.5	Yes	Yes
RT Ion Plan Storage	1.2.840.10008.5.1.4.1.1.481.8	Yes	Yes
RT Dose Storage	1.2.840.10008.5.1.4.1.1.481.2	Yes	Yes
RT Image Storage	1.2.840.10008.5.1.4.1.1.481.1	No	Yes
Query/Retrieve			
Patient Root Query/Retrieve Information Model – FIND	1.2.840.10008.5.1.4.1.2.1.1	No	Yes
Study Root Query/Retrieve Information Model – FIND	1.2.840.10008.5.1.4.1.2.2.1	No	Yes
Patient Root Query/Retrieve Information Model – MOVE	1.2.840.10008.5.1.4.1.2.1.2	No	Yes

## 2 INTRODUCTION

### 2.1 REVISION HISTORY

Date	Version	Comment
2020-04-29	1.0	First version.

### 2.2 REMARKS

The scope of this DICOM Conformance Statement is to facilitate integration between RayStation and other applications. The document should be read and understood in conjunction with the DICOM Standard. DICOM adherence by itself does not guarantee interoperability. However, this document facilitate a first hand comparison for interoperability between different applications supporting compatible DICOM functionality.

This document is written for users that need to understand how RayStation will integrate into their healthcare facility. This includes both those responsible for overall imaging network policy and architecture, as well as integrators who need to have a detailed understanding of the DICOM features of the product. This document contains some basic DICOM definitions so that any reader may understand how this product implements DICOM features. However, integrators are expected to fully understand all the DICOM terminology, how the tables in this document relate to the product's functionality, and how that functionality integrates with other devices that support compatible DICOM features.

This document is not supposed to replace validation with other systems to ensure proper exchange of intended information. The user should be aware of the following important issues:

- The comparison of different Conformance Statements is just the first step towards assessing interconnectivity and interoperability between the product and other DICOM conformant equipment.
- Test procedures should be defined and executed to validate the required level of interoperability with specific compatible DICOM equipment, as established by the healthcare facility.

When using RayStation together with other software, the DICOM conformance statements must be compared and relevant validation tests run. The DICOM standard by itself does not guarantee interoperability. The Conformance Statement does, however, facilitate a first-level comparison for interoperability between different applications supporting compatible DICOM functionality. RaySearch is also active within the IHE-RO. Contact RaySearch for more info regarding adherence to IHE-RO profiles.


## 3 NETWORKING AND MEDIA INTERCHANGE

### 3.1 IMPLEMENTATION MODEL

#### 3.1.1 Application Data Flow Diagram

RayStation is a treatment planning system, i.e. a software program for planning and analysis of radiation therapy plans. Typically, a treatment plan is created by importing patient images obtained from a CT scanner, importing and/or defining regions of interest (possibly with the help of MR and PET-images), deciding on a treatment setup and objectives, optimizing the treatment parameters, comparing rival plans to find the best compromise, computing the clinical dose distribution, approving the plan and finally exporting the data.

The diagram below describes how RayStation interacts with the DICOM world. RayStation consists of two DICOM Application Entities, the client and the standalone Storage SCP service.

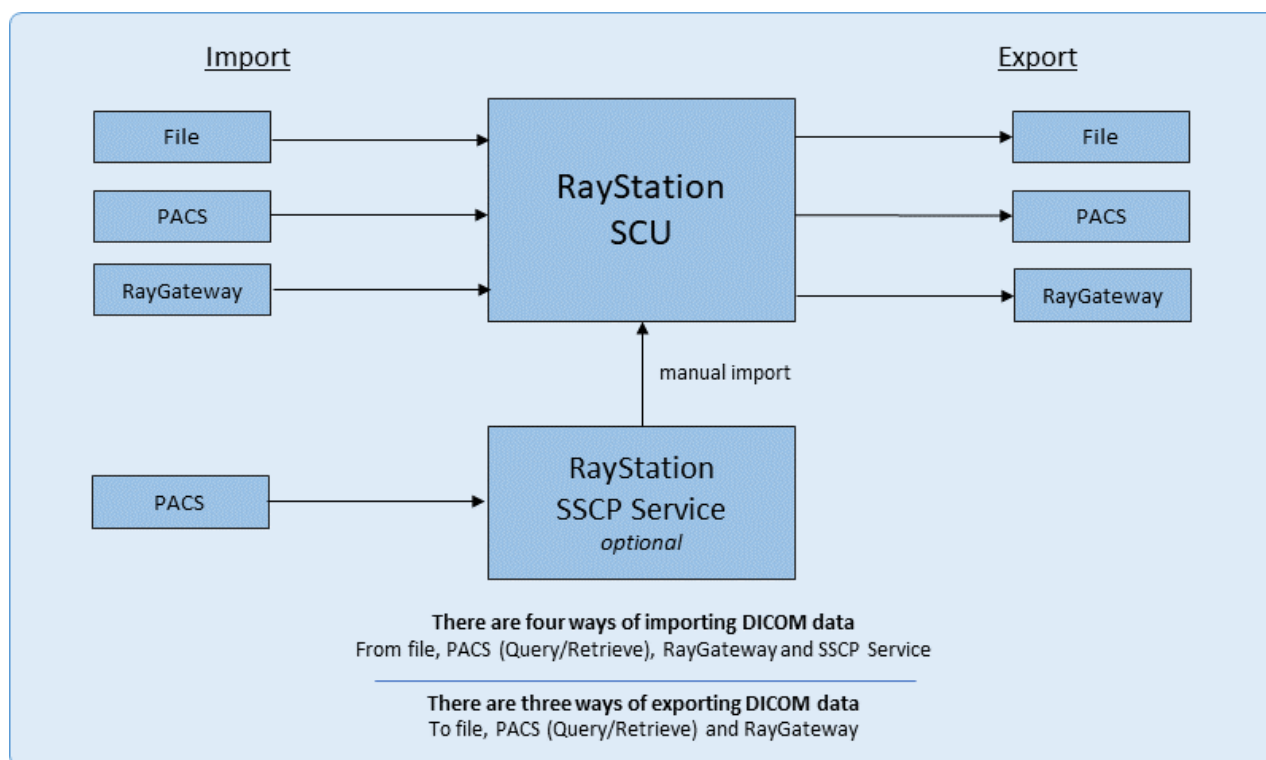

#### 3.1.2 Functional Definition of AEs

##### 3.1.2.1 Functional Definition of "RayStation Application Entity"

RayStation can act both as an SCU and SCP. When RayStation performs import using Query/Retrieve an internal Storage SCP temporarily is invoked automatically when querying the PACS. AE settings for RayStation are done in the separate ClinicSettings application.

##### 3.1.2.2 Functional Definition of "RayStation Storage SCP Service"

For transfers when the RayStation client is not running the RayStation Storage SCP Service is provided. This is a standalone Storage SCP service accepting incoming DICOM transfers and saving the files to a network share. The user then manually imports the data to the RayStation database using a version of the file-based import in the RayStation client.

### 3.2 TRANSFER SYNTAX

RayStation import supports Implicit VR Little Endian Transfer Syntax (1.2.840.10008.1.2) and Explicit VR Little Endian Transfer Syntax (1.2.840.10008.1.2.1). The Transfer Syntax used for export can be configured in the Clinic Settings application. Additional Transfer Syntaxes might be supported for import, such as compressed images. Please contact RaySearch if more information is required.

### 3.3 IMPORTING DATA

### 3.3.1 Supported Sources

There are four ways to import DICOM data into RayStation: File, PACS, RayGateway and Storage SCP. Imported CT, MR and PET datasets are stored without modification. All RT modalities are converted into internal representations and the original data will NOT be available for export.

Applying filters permanently modifies the imported data but usually does not regenerate the UID values.

#### 3.3.1.1 File

Using File source will scan directories and filter files according to the selected file mask. DICOM datasets will be presented according to a Patient/Study/Series model similar to Query/Retrieve.

Sub-directories is scanned when expanded. The complete directory structure can be scanned by checking the top directory. This will recursively scan all directories and can take a long time. Cancellation of this process is not supported in this version of RayStation.

#### 3.3.1.2 PACS

Using PACS source will use the Query/Retrieve services using DIMSE-C C-FIND and C-MOVE. PACS endpoints must be configured in the Clinic Settings applications. Universal matching for patients is not allowed (when the search term is empty or contains wildcards only, e.g. "\*", "\*\*", or "?").

Only a single patient can be selected. Patient, Study and Series nodes can be expanded to list more details, which initiates C-FIND requests. Complete patient data can be retrieved by checking the patient node. It is not possible to cancel ongoing queries.

Once the Import button is clicked, RayStation will C-MOVE requests for all selected series. Cancellation of this process is not supported in this version of RayStation.

#### 3.3.1.3 RayGateway

Using RayGateway source will behave differently depending on the endpoint system. Please contact RaySearch if more information is required.

#### 3.3.1.4 Storage SCP

If "Storage SCP Directory" has been configured in the Clinic Settings application the Storage SCP source will be available. The directory is scanned for session folders are listed in the same way as a directory structure. The user also has the option of automatically deleting the files from the DICOM Inbox after successful import.

The RayStation Storage SCP service is designed to listen to incoming DICOM datasets and writes them to the "DICOM Inbox". All incoming DICOM datasets are stored in session folders grouped by Remote AE Title, Date and Patient ID.

If the service is unable to save the incoming files (e.g. not enough storage space), the network association is aborted. As the DICOM Inbox is only regarded as a temporary storage location, the service is setup to delete the folders in the DICOM Inbox after a defined number of days.

## 3.4 EXPORTING DATA

### 3.4.1 Supported Targets

RayStation supports DICOM export to file, PACS or RayGateway. Exported image sets will be identical to the imported data with the exception that patient information is replaced to match the RayStation patient (as shown in Patient Data Management).

### 3.4.2 Generated RayStation UID

RayStation generates UID with the prefix 1.2.752.243.1.1. No assumptions shall be made on the format and RaySearch reserves the right to change this at any time.

### 3.4.3 Checksums

Checksums can be enabled in the Clinic Settings application. If this feature is enabled exported RT Plans will contain the following private attributes carrying the checksum details:

- RaySearch Checksum Algorithm Version {4001,1060}
- RaySearch Checksum Data {4001,1060}

Contact RaySearch for a specification of the algorithm to verify this checksum.

### 3.4.4 Anonymization

RayStation supports the Basic Profile de-identifier. In addition the following options are supported:

- Retain Longitudinal Temporal Information with Full Dates Option
- Retain Device Identity Option
- Retain Institution Identity Option
- Retain UIDs Option
- Retain Safe Private Option

See Table E.1-1. Application Level Confidentiality Profile Attributes in DICOM PS3.15 2018e.

Anonymization can be used when exporting DICOM and when creating patient backups. Both these features use the same underlying mechanism for cleaning DICOM datasets. Note the following limitations or deviations from the Basic Profile and supported options:

- Patient's Name (0010,0010) is replaced by a user-supplied value.
- Patient ID (0010,0020) is replaced by a user-supplied value.
- Patient's Birth Date (0010,0030) is replaced with the current date.
- Patient's Sex (0010,0040) is replaced with "0".
- Private attributes (gggg.eeee) where gggg is odd, are completely removed if the Basic Profile is used. *This is according to the standard definition of the Basic Profile but explicitly noted since this can significantly strip a DICOM dataset from information.*
- All Sequences in DICOM PS3.15 2018e Table E.1-1 attributes with the action "D" (replace with a non-zero length value that may be a dummy value and consistent with the VR) will be removed.
- If the Retain Device Identity Option is selected, private attributes containing device information is also retained.
- If the Retain Safe Private Option is selected, private attributes containing no patient information is retained.



## 4 SECURITY

No Security Profiles are supported.

## 5 SUPPORT OF CHARACTER SETS

RayStation supports the Default Character Repertoire, ISO\_IR 100 (Latin-1) and ISO\_IR 192 (Unicode). The exported character set is ISO\_IR 100 by default and can be changed in the Clinic Settings application.

Importing datasets with other character sets might lead to undefined behavior. Validation of the character set is done when importing or exporting data and warnings will be presented to the user. Usage of non-latin characters must be validated carefully by end users.

## 6 KNOWN ISSUES

### 6.1 IMPORT OF RESPIRATORY PHASE GATING DUTY CYCLE

There is a known issue when importing plan with the IMPAC private attribute Respiratory Phase Gating Duty Cycle (300B,1020). The attribute tag conflicts with another attribute definition and will therefore not be imported properly. Two possible work arounds:

1. Exporting a plan with Respiratory Phase Gating Duty Cycle first, the attribute definition will be registered in the DICOM dictionary and import will be possible.
2. Import the plan twice. The attribute definition will be registered in the DICOM dictionary and import will be successful in the second attempt.

Both work arounds must be done within the same RayStation session.

## 7 ANNEXES

### 7.1 IMPORT IOD DEFINITIONS

#### 7.1.1 CT Image IOD

Imported images are stored internally. Filters modifies data permanently.

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	Yes
	Clinical Trial Study Module	Yes
Series	General Series Module	Yes
Frame of Reference	Frame of Reference Module	Yes
Equipment	General Equipment Module	Yes
Image	General Image Module	No
	Image Plane Module	Yes
	Image Pixel Module	Yes
	Contrast/Bolus Module	No
	CT Image Module	Yes
	Multi-energy CT Image Module	No
	SOP Common Module	Yes

##### 7.1.1.1 Patient Module

Patient attributes are read and stored when importing a new patient.

Patient ID and the first component (Last Name) of Patient's Name is the unique identifier. Differences in Patient's Sex and Patient's Birth Date are checked.

Attribute name	Tag	Vr	Type	Comment
Patient's Name	{0010,0010}	PN	2	
Patient ID	{0010,0020}	LO	2	
Patient's Birth Date	{0010,0030}	DA	2	
Patient's Sex	{0010,0040}	CS	2	Empty value defaults to Unknown gender. Supported values: <ul style="list-style-type: none"> <li>• M - Male</li> <li>• F - Female</li> <li>• 0 - Other</li> </ul>
Patient Identity Removed	{0012,0062}	CS	3	Supported values: YES, NO.
De-identification Method	{0012,0063}	LO	1C	

##### 7.1.1.2 General Study Module

Attribute name	Tag	Vr	Type	Comment
Study Instance UID	(0020,000D)	UI	1	Stored when importing a new case but not used.
Study Date	(0008,0020)	DA	2	
Study Time	(0008,0030)	TM	2	
Referring Physician's Name	(0008,0090)	PN	2	Stored as Case Physician when importing a new case.
Study ID	(0020,0010)	SH	2	
Accession Number	(0008,0050)	SH	2	
Study Description	(0008,1030)	LO	3	

### 7.1.1.3 Clinical Trial Study Module

Attribute name	Tag	Vr	Type	Comment
Clinical Trial Time Point ID	(0012,0050)	LO	2	
Consent for Clinical Trial Use Sequence	(0012,0083)	SQ	3	Only used for display of contract information. Mainly for internal use.
>Distribution Type	(0012,0084)	CS	1C	
>Clinical Trial Protocol ID	(0012,0020)	LO	1C	
>Consent for Distribution Flag	(0012,0085)	CS	1	Supported values: YES, NO, WITHDRAWN.
>Clinical Trial Subject ID	(0012,0040)	LO	3	
>RaySearch Private Creator	(4001,0010)	LO	3	RAYSEARCHLABS 2.0
>Consent Expiry Date	(4001,1067)	DA	3	Date when consent expires.

### 7.1.1.4 General Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	(0008,0060)	CS	1	Supported value: CT.
Series Instance UID	(0020,000E)	UI	1	Stored to verify references from RT Structure Sets.
Series Number	(0020,0011)	IS	2	
Series Date	(0008,0021)	DA	3	
Series Time	(0008,0031)	TM	3	
Protocol Name	(0018,1030)	LO	3	Used to identify Imaging System during import.
Series Description	(0008,103E)	LO	3	
Patient Position	(0018,5100)	CS	2C	Supported values: HFS, HFP, FFS, FFP, HFDR, HFDL, FFDR, FFDL.

### 7.1.1.5 Frame of Reference Module

Attribute name	Tag	Vr	Type	Comment
Frame of Reference UID	(0020,0052)	UI	1	Stored internally to be used if referenced from other dataset.

### 7.1.1.6 General Equipment Module

Attribute name	Tag	Vr	Type	Comment
Station Name	(0008,1010)	SH	3	Used to identify the correct CT-to-Density table.
Pixel Padding Value	(0028,0120)	US or SS	1C	If pixel padding is present within the region used for dose computation, the resulting dose may be wrong.

#### 7.1.1.7 Image Plane Module

Attribute name	Tag	Vr	Type	Comment
Pixel Spacing	(0028,0030)	DS	1	Cannot be greater than 5 mm.
Image Orientation (Patient)	(0020,0037)	DS	1	Image orientation must be transversal.
Image Position (Patient)	(0020,0032)	DS	1	Converted to internal image corner. Used to find slice direction, which must be constant throughout the series and along the Z-axis.

#### 7.1.1.8 Image Pixel Module

Attribute name	Tag	Vr	Type	Comment
Samples per Pixel	(0028,0002)	US	1	Supported value: 1.
Photometric Interpretation	(0028,0004)	CS	1	Supported value: MONOCHROME2.
Rows	(0028,0010)	US	1	
Columns	(0028,0011)	US	1	
Bits Allocated	(0028,0100)	US	1	Supported value: 16.
Bits Stored	(0028,0101)	US	1	
High Bit	(0028,0102)	US	1	
Pixel Representation	(0028,0103)	US	1	Supported values: 0, 1.
Pixel Data	(7FE0,0010)	OB or OW	1C	
Pixel Padding Range Limit	(0028,0121)	US or SS	1C	

#### 7.1.1.9 CT Image Module

Attribute name	Tag	Vr	Type	Comment
Image Type	(0008,0008)	CS	1	
Samples per Pixel	(0028,0002)	US	1	Supported value: 1.
Photometric Interpretation	(0028,0004)	CS	1	Supported value: MONOCHROME2.
Bits Allocated	(0028,0100)	US	1	Supported value: 16.
Bits Stored	(0028,0101)	US	1	Must be greater than 0 and less than or equal to Bits Allocated. Supported values: 12, 13, 14, 15, 16.
High Bit	(0028,0102)	US	1	Must be equal to "Bits Stored - 1".
Rescale Intercept	(0028,1052)	DS	1	
Rescale Slope	(0028,1053)	DS	1	

Rescale Type	{0028,1054}	LO	1C	Supported value: HU.
--------------	-------------	----	----	----------------------

### 7.1.1.10 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	{0008,0016}	UI	1	Supported value: 1.2.840.10008.5.1.4.1.1.2.
SOP Instance UID	{0008,0018}	UI	1	Stored internally to be used if referenced from other dataset.
Specific Character Set	{0008,0005}	CS	1C	Supported values: ISO_IR 100, ISO_IR 192.
Instance Creation Date	{0008,0012}	DA	3	
Instance Creation Time	{0008,0013}	TM	3	
Contributing Equipment Sequence	{0018,A001}	SQ	3	Only used for display of contract information. Mainly for internal use.
>Purpose of Reference Code Sequence	{0040,A170}	SQ	1	
>>Code Value	{0008,0100}	SH	1C	Supported value: <ul style="list-style-type: none"> <li>109103 - Modifying Equipment</li> </ul>
>>Code Meaning	{0008,0104}	LO	1	Supported value: Modifying Equipment.
>Manufacturer	{0008,0070}	LO	1	
>Manufacturer's Model Name	{0008,1090}	LO	3	
>Software Versions	{0018,1020}	LO	3	
>Contribution DateTime	{0018,A002}	DT	3	
>Contribution Description	{0018,A003}	ST	3	

### 7.1.2 MR Image IOD

Imported images are stored internally. Filters modifies data permanently.

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	Yes
	Clinical Trial Study Module	Yes
Series	General Series Module	Yes
Frame of Reference	Frame of Reference Module	Yes
Equipment	General Equipment Module	Yes
Image	General Image Module	No
	Image Plane Module	Yes
	Image Pixel Module	Yes
	Contrast/Bolus Module	No
	MR Image Module	Yes

	SOP Common Module	Yes
--	-------------------	-----

### 7.1.2.1 Patient Module

Patient attributes are read and stored when importing a new patient.

Patient ID and the first component (Last Name) of Patient's Name is the unique identifier. Differences in Patient's Sex and Patient's Birth Date are checked.

Attribute name	Tag	Vr	Type	Comment
Patient's Name	(0010,0010)	PN	2	
Patient ID	(0010,0020)	LO	2	
Patient's Birth Date	(0010,0030)	DA	2	
Patient's Sex	(0010,0040)	CS	2	Empty value defaults to Unknown gender. Supported values: <ul style="list-style-type: none"> <li>• M - Male</li> <li>• F - Female</li> <li>• O - Other</li> </ul>
Patient Identity Removed	(0012,0062)	CS	3	Supported values: YES, NO.
De-identification Method	(0012,0063)	LO	1C	

### 7.1.2.2 General Study Module

Attribute name	Tag	Vr	Type	Comment
Study Instance UID	(0020,0000)	UI	1	Stored when importing a new case but not used.
Study Date	(0008,0020)	DA	2	
Study Time	(0008,0030)	TM	2	
Referring Physician's Name	(0008,0090)	PN	2	Stored as Case Physician when importing a new case.
Study ID	(0020,0010)	SH	2	
Accession Number	(0008,0050)	SH	2	

### 7.1.2.3 Clinical Trial Study Module

Attribute name	Tag	Vr	Type	Comment
Clinical Trial Time Point ID	(0012,0050)	LO	2	
Consent for Clinical Trial Use Sequence	(0012,0083)	SQ	3	Only used for display of contract information. Mainly for internal use.
>Distribution Type	(0012,0084)	CS	1C	
>Clinical Trial Protocol ID	(0012,0020)	LO	1C	
>Consent for Distribution Flag	(0012,0085)	CS	1	Supported values: YES, NO, WITHDRAWN.
>Clinical Trial Subject ID	(0012,0040)	LO	3	
>RaySearch Private Creator	(4001,0010)	LO	3	RAYSEARCHLABS 2.0
>Consent Expiry Date	(4001,1067)	DA	3	Date when consent expires.

### 7.1.2.4 General Series Module



Attribute name	Tag	Vr	Type	Comment
Modality	(0008,0060)	CS	1	Supported value: MR.
Series Instance UID	(0020,000E)	UI	1	Stored to verify references from RT Structure Sets.
Series Number	(0020,0011)	IS	2	
Series Date	(0008,0021)	DA	3	
Series Time	(0008,0031)	TM	3	
Protocol Name	(0018,1030)	LO	3	Used to identify Imaging System during import.
Series Description	(0008,103E)	LO	3	
Patient Position	(0018,5100)	CS	2C	Supported values: HFS, HFP, FFS, FFP.

#### 7.1.2.5 Frame of Reference Module

Attribute name	Tag	Vr	Type	Comment
Frame of Reference UID	(0020,0052)	UI	1	Stored internally to be used if referenced from other dataset.

#### 7.1.2.6 General Equipment Module

Attribute name	Tag	Vr	Type	Comment
Pixel Padding Value	(0028,0120)	US or SS	1C	If pixel padding is present within the region used for dose computation, the resulting dose may be wrong.

#### 7.1.2.7 Image Plane Module

Attribute name	Tag	Vr	Type	Comment
Pixel Spacing	(0028,0030)	DS	1	Cannot be greater than 5 mm.
Image Orientation (Patient)	(0020,0037)	DS	1	All image orientations are supported, as long as row and column directions are orthogonal.
Image Position (Patient)	(0020,0032)	DS	1	Converted to internal image corner. Used to find slice direction, which must be constant throughout the series.

#### 7.1.2.8 Image Pixel Module

Attribute name	Tag	Vr	Type	Comment
Samples per Pixel	(0028,0002)	US	1	Supported value: 1.
Photometric Interpretation	(0028,0004)	CS	1	Supported value: MONOCHROME2.
Rows	(0028,0010)	US	1	
Columns	(0028,0011)	US	1	
Bits Allocated	(0028,0100)	US	1	Supported value: 16.
Bits Stored	(0028,0101)	US	1	Must be equal to or less than Bits Allocated.
High Bit	(0028,0102)	US	1	
Pixel Representation	(0028,0103)	US	1	Supported values: 0, 1.

Pixel Data	{7FE0,0010}	OB or OW	1C	
Pixel Padding Range Limit	{0028,0121}	US or SS	1C	

### 7.1.2.9 MR Image Module

Attribute name	Tag	Vr	Type	Comment
Image Type	{0008,0008}	CS	1	
Samples per Pixel	{0028,0002}	US	1	Supported value: 1.
Photometric Interpretation	{0028,0004}	CS	1	Supported value: MONOCHROME2.
Bits Allocated	{0028,0100}	US	1	Supported value: 16.
Bits Stored	{0028,0101}	US	1	
High Bit	{0028,0102}	US	1	
Scanning Sequence	{0018,0020}	CS	1	
Sequence Variant	{0018,0021}	CS	1	

### 7.1.2.10 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	{0008,0016}	UI	1	Supported value: 1.2.840.10008.5.1.4.1.1.4.
SOP Instance UID	{0008,0018}	UI	1	Stored internally to be used if referenced from other dataset.
Specific Character Set	{0008,0005}	CS	1C	Supported values: ISO_IR 100, ISO_IR 192.
Instance Creation Date	{0008,0012}	DA	3	
Instance Creation Time	{0008,0013}	TM	3	
Contributing Equipment Sequence	{0018,A001}	SQ	3	Only used for display of contract information. Mainly for internal use.
>Purpose of Reference Code Sequence	{0040,A170}	SQ	1	
>>Code Value	{0008,0100}	SH	1C	Supported value: <ul style="list-style-type: none"> <li>109103 - Modifying Equipment</li> </ul>
>>Coding Scheme Designator	{0008,0102}	SH	1C	Supported value: DCM.
>>Code Meaning	{0008,0104}	LO	1	Supported value: Modifying Equipment.
>Manufacturer	{0008,0070}	LO	1	
>Manufacturer's Model Name	{0008,1090}	LO	3	
>Software Versions	{0018,1020}	LO	3	
>Contribution DateTime	{0018,A002}	DT	3	
>Contribution Description	{0018,A003}	ST	3	

### 7.1.3 RT Dose IOD

Imported RT Dose objects are converted into internal representation. Exported doses will therefore not be identical to the imported object.

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	No
Series	RT Series Module	Yes
Frame of Reference	Frame of Reference Module	Yes
Equipment	General Equipment Module	Yes
Dose	General Image Module	No
	Image Plane Module	Yes
	Image Pixel Module	Yes
	Multi-frame Module	No
	RT Dose Module	Yes
	Structure Set Module	No
	ROI Contour Module	No
	RT Dose ROI Module	No
	SOP Common Module	Yes
Frame Extraction Module	No	

### 7.1.3.1 Patient Module

Patient attributes only used to match patient.

Attribute name	Tag	Vr	Type	Comment
Patient's Name	{0010,0010}	PN	2	
Patient ID	{0010,0020}	LO	2	
Patient's Birth Date	{0010,0030}	DA	2	
Patient's Sex	{0010,0040}	CS	2	Supported values: <ul style="list-style-type: none"> <li>• M - Male</li> <li>• F - Female</li> <li>• 0 - Other</li> </ul>

### 7.1.3.2 RT Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	{0008,0060}	CS	1	Supported value: RTDOSE.
Series Instance UID	{0020,000E}	UI	1	Only read for display during import. Not stored.

### 7.1.3.3 Frame of Reference Module

Attribute name	Tag	Vr	Type	Comment
Frame of Reference UID	{0020,0052}	UI	1	Used to match referenced image set. Not stored.

### 7.1.3.4 General Equipment Module

Attribute name	Tag	Vr	Type	Comment
Pixel Padding Value	{0028,0120}	US or SS	1C	If pixel padding is present within the region used for dose computation, the resulting dose may be wrong.

### 7.1.3.5 Image Plane Module

Attribute name	Tag	Vr	Type	Comment
Pixel Spacing	{0028,0030}	DS	1	
Image Orientation (Patient)	{0020,0037}	DS	1	All non-oblique image orientations are supported (i.e. row and column directions must be along an axis).
Image Position (Patient)	{0020,0032}	DS	1	Note that this is the center of the first pixel, i.e. the corner of the image offset by half a pixel.

### 7.1.3.6 Image Pixel Module

Attribute name	Tag	Vr	Type	Comment
Samples per Pixel	{0028,0002}	US	1	Supported value: 1.
Photometric Interpretation	{0028,0004}	CS	1	Supported value: MONOCHROME2.
Rows	{0028,0010}	US	1	
Columns	{0028,0011}	US	1	
Bits Allocated	{0028,0100}	US	1	Supported values: 16, 32.
Bits Stored	{0028,0101}	US	1	Must be equal to Bits Allocated.
High Bit	{0028,0102}	US	1	HighBit must be BitsStored-1.
Pixel Representation	{0028,0103}	US	1	Supported value: 0.
Pixel Data	{7FE0,0010}	OB or OW	1C	
Pixel Padding Range Limit	{0028,0121}	US or SS	1C	

### 7.1.3.7 RT Dose Module

Attribute name	Tag	Vr	Type	Comment
Samples per Pixel	{0028,0002}	US	1C	Supported value: 1.
Photometric Interpretation	{0028,0004}	CS	1C	Supported value: MONOCHROME2.
Bits Allocated	{0028,0100}	US	1C	Supported values: 16, 32.
Bits Stored	{0028,0101}	US	1C	Bits Stored must be equal to Bits Allocated.
High Bit	{0028,0102}	US	1C	High Bit must be Bits Stored - 1.
Pixel Representation	{0028,0103}	US	1C	Supported value: 0.
Dose Units	{3004,0002}	CS	1	Supported value: GY.
Dose Type	{3004,0004}	CS	1	Supported values: PHYSICAL, EFFECTIVE.

Dose Summation Type	(3004,000A)	CS	1	<p>Note: Summation type BEAM is interpreted as the dose contribution from this beam over the entire treatment course, not per fraction.</p> <p>Supported values:</p> <ul style="list-style-type: none"> <li>• PLAN</li> <li>• BEAM</li> <li>• FRACTION</li> <li>• TOTALHOMO - Interpreted as PLAN</li> <li>• TOTALHETERO - Interpreted as PLAN</li> <li>• EVALUATION - Interpreted as PLAN</li> <li>• BNCT_GAMMA - BNCT Gamma Plan Dose Component</li> <li>• BNCT_BORON - BNCT Boron Plan Dose Component</li> <li>• BNCT_NITROGEN - BNCT Nitrogen Plan Dose Component</li> <li>• BNCT_HYDROGEN - BNCT Hydrogen Plan Dose Component</li> <li>• BNCT_NORMAL_RBE - BNCT Normal RBE Plan Dose Component</li> <li>• BNCT_TUMOR_RBE - BNCT Tumor RBE Plan Dose Component</li> <li>• MULTI_PLAN</li> </ul>
Referenced RT Plan Sequence	(300C,0002)	SQ	1C	Is required to contain only one item. For doses with summation type BEAM the referenced plan must be found and all beam doses must be imported together. For doses with summation type PLAN a dummy plan will be created if the referenced plan couldn't be found but an examination with a matching Frame of Reference UID is found.
>Referenced SOP Class UID	(0008,1150)	UI	1	
>Referenced SOP Instance UID	(0008,1155)	UI	1	Plan UID, used to establish references.
>Referenced Fraction Group Sequence	(300C,0020)	SQ	1C	
>>Referenced Fraction Group Number	(300C,0022)	IS	1	
>>Referenced Beam Sequence	(300C,0004)	SQ	1C	
>>>Referenced Beam Number	(300C,0006)	IS	1	Identify a beam within a plan
Grid Frame Offset Vector	(3004,000C)	DS	1C	If the first element of the Grid Frame Offset Vector is zero, the vector should be interpreted as relative positive offsets in the direction of the cross product of the row and column directions. This means that if Image Orientation (Patient) (0020,0037) is for example (1,0,0,0,-1,0), the cross product will be negative and hence the elements in the Grid Frame Offset Vector should be read/written with opposite sign. It should be noted that some manufacturers are non-compliant regarding this and do not change the sign of the elements in the Grid Frame Offset Vector and could cause the dose grid to not be positioned correctly around the patient.
Dose Grid Scaling	(3004,000E)	DS	1C	
Tissue Heterogeneity Correction	(3004,0014)	CS	3	Supported values: WATER, IMAGE, ROI_OVERRIDE.
RaySearch Private Creator	(4001,0010)	LO	3	RAYSEARCHLABS 2.0

ROI Index Pixel Data	{4001,1016}	OB	3	
----------------------	-------------	----	---	--

#### 7.1.3.8 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	{0008,0016}	UI	1	Supported value: 1.2.840.10008.5.1.4.1.1.481.2.
SOP Instance UID	{0008,0018}	UI	1	Stored internally to references from other objects.
Specific Character Set	{0008,0005}	CS	1C	Supported values: ISO_IR 100, ISO_IR 192.
Instance Creation Date	{0008,0012}	DA	3	Stored as last save date.
Instance Creation Time	{0008,0013}	TM	3	Stored as last save time.

#### 7.1.4 RT Structure Set IOD

Imported RT Structure Sets are converted into internal representation. Exported structure sets will therefore not be identical to the imported object.

RT Structure Sets must referenced a single image series.

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	No
Series	RT Series Module	Yes
Equipment	General Equipment Module	Yes
Frame of Reference	Frame of Reference Module	Yes
Structure Set	Structure Set Module	Yes
	ROI Contour Module	Yes
	RT ROI Observations Module	Yes
	Approval Module	Yes
	SOP Common Module	Yes

##### 7.1.4.1 Patient Module

Patient attributes only used to match patient.

Attribute name	Tag	Vr	Type	Comment
Patient's Name	{0010,0010}	PN	2	
Patient ID	{0010,0020}	LO	2	
Patient's Birth Date	{0010,0030}	DA	2	
Patient's Sex	{0010,0040}	CS	2	Supported values: <ul style="list-style-type: none"> <li>• M - Male</li> <li>• F - Female</li> <li>• O - Other</li> </ul>

##### 7.1.4.2 RT Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	(0008,0060)	CS	1	Supported value: RTSTRUCT.
Series Instance UID	(0020,000E)	UI	1	Only read for display during import. Not stored.

#### 7.1.4.3 General Equipment Module

Attribute name	Tag	Vr	Type	Comment
Manufacturer	(0008,0070)	LO	2	Used to handle ROI Physical Properties from TomoTherapy Hi-Art and Nucletron during import. Not stored.

#### 7.1.4.4 Frame of Reference Module

Attribute name	Tag	Vr	Type	Comment
Frame of Reference UID	(0020,0052)	UI	1	

#### 7.1.4.5 Structure Set Module

Attribute name	Tag	Vr	Type	Comment
Structure Set Label	(3006,0002)	SH	1	Not stored.
Structure Set Date	(3006,0008)	DA	2	Stored as last save date.
Structure Set Time	(3006,0009)	TM	2	Stored as last save time.
Referenced Frame of Reference Sequence	(3006,0010)	SQ	3	May contain one or more items. A single referenced examinations must exist. Contours on unknown examinations are ignored.
>Frame of Reference UID	(0020,0052)	UI	1	Used to match referenced image set.
>RT Referenced Study Sequence	(3006,0012)	SQ	3	Must contain exactly one referenced study.
>>Referenced SOP Class UID	(0008,1150)	UI	1	
>>Referenced SOP Instance UID	(0008,1155)	UI	1	
>>RT Referenced Series Sequence	(3006,0014)	SQ	1	Must contain exactly one referenced series.
>>>Series Instance UID	(0020,000E)	UI	1	
>>>Contour Image Sequence	(3006,0016)	SQ	1	
>>>>Referenced SOP Class UID	(0008,1150)	UI	1	
>>>>Referenced SOP Instance UID	(0008,1155)	UI	1	Referenced image has to be imported to patient (prior or at the same time as the RTSTRUCT dataset is imported).
Structure Set ROI Sequence	(3006,0020)	SQ	1	Number of elements must match number of elements in ROI Contour Sequence (3006,0039). At least one sequence must be defined.
>ROI Number	(3006,0022)	IS	1	Stored but only used internally.
>Referenced Frame of Reference UID	(3006,0024)	UI	1	
>ROI Name	(3006,0026)	LO	2	ROI:s are shared internally between multiple structure sets. If a ROI with the same name is already contoured on the referenced examination these contours will be skipped.

## 7.1.4.6 ROI Contour Module

Attribute name	Tag	Vr	Type	Comment
ROI Contour Sequence	{3006,0039}	SQ	1	
>Referenced ROI Number	{3006,0084}	IS	1	
>ROI Display Color	{3006,002A}	IS	3	Used for visualization.
>Contour Sequence	{3006,0040}	SQ	3	
>>Attached Contours	{3006,0049}	IS	3	Information about attached contours are not imported. Must not have a value.
>>Contour Image Sequence	{3006,0016}	SQ	3	
>>>Referenced SOP Class UID	{0008,1150}	UI	1	
>>>Referenced SOP Instance UID	{0008,1155}	UI	1	Must be found in the referenced Examination.
>>Contour Geometric Type	{3006,0042}	CS	1	
>>Number of Contour Points	{3006,0046}	IS	1	
>>Contour Data	{3006,0050}	DS	1	Converted to internal contour representation. All contour vertices must be located on the same slice.

## 7.1.4.7 RT ROI Observations Module

Attribute name	Tag	Vr	Type	Comment
RT ROI Observations Sequence	{3006,0080}	SQ	1	
>Observation Number	{3006,0082}	IS	1	
>Referenced ROI Number	{3006,0084}	IS	1	
>ROI Observation Label	{3006,0085}	SH	3	ROI Observation Label is "Localization Poi" and ROI Interpreted Type is "MARKER" the ROI is automatically converted to a localization point.
>RT ROI Interpreted Type	{3006,00A4}	CS	2	All types in the DICOM standard are supported and converted to internal equivalent. Furthermore if ROI Observation Label has value "Localization Poi" and ROI Interpreted Type is "MARKER", then the ROI is automatically converted to a localization point. Supported values: ACQ_ISOCENTER, AVOIDANCE, BOLUS, BRACHY_ACCESSORY, BRACHY_CHANNEL, BRACHY_CHNL_SHLD, BRACHY_SRC_APP, CAVITY, CONTRAST_AGENT, CONTROL, CTV, DOSE_REGION, EXTERNAL, FIXATION, GTV, INIT_LASER_ISO, INIT_MATCH_ISO, IRRAD_VOLUME, ISOCENTER, MARKER, ORGAN, PTV, REGISTRATION, SUPPORT, TREATED_VOLUME, NONE.
>Material ID	{300A,00E1}	SH	3	Imported as material name. Note that material names must be unique.



>ROI Physical Properties Sequence	{3006,00B0}	SQ	3	<p>The values REL_MASS_DENSITY, REL_ELEC_DENSITY, EFFECTIVE_Z, EFF_Z_PER_A and ELEM_FRACTION are supported. The physical properties expressed in REL_MASS_DENSITY, REL_ELEC_DENSITY and EFF_Z_PER_A are connected as <math>REL\_ELEC\_DENSITY = REL\_MASS\_DENSITY * EFF\_Z\_PER\_A / EFF\_Z\_PER\_A\ water</math>. Any combination of two of these three properties are regarded as giving information also about the third. If all four are available REL_ELEC_DENSITY is disregarded in the import. If ELEM_FRACTION and REL_MASS_DENSITY are available, other given properties are not imported but are computed from these two. If ELEM_FRACTION and REL_ELEC_DENSITY and EFF_Z_PER_A are available but not REL_MASS_DENSITY, then REL_MASS_DENSITY is computed using the relation above and properties other than the resulting REL_MASS_DENSITY and ELEM_FRACTION are not imported but are computed from these two. If neither ELEM_FRACTION nor EFFECTIVE_Z are available RayStation must make assumptions about the material properties. RayStation then assumes the material to be water with REL_MASS_DENSITY set to the REL_ELEC_DENSITY from the RT Struct file if it is available or can be computed. This assumption leads to photon attenuations which are appropriate for materials with EFF_Z_PER_A similar to water. For high-Z materials, this can lead to significant dose deviations. If REL_ELEC_DENSITY cannot be established, REL_MASS_DENSITY is used.</p>
>>ROI Physical Property	{3006,00B2}	CS	1	<p>Relative mass and electron densities must be between 0.0 - 22.7. MEAN_EXCI_ENERGY depicts the Mean Excitation Energy for a given material. Supported values:</p> <ul style="list-style-type: none"> <li>● REL_MASS_DENSITY</li> <li>● REL_ELEC_DENSITY</li> <li>● EFFECTIVE_Z</li> <li>● EFF_Z_PER_A</li> <li>● ELEM_FRACTION</li> <li>● MEAN_EXCI_ENERGY</li> <li>● CBE_BORON - Used for BNCT dose computation</li> <li>● RBE_HYDROGEN - Used for BNCT dose computation</li> <li>● RBE_NITROGEN - Used for BNCT dose computation</li> <li>● RBE_GAMMA - Used for BNCT dose computation</li> <li>● BORON_CONC_RATIO - Boron Concentration Ratio. Used for BNCT dose computation</li> <li>● BNCT_DOSE_LIMIT - Used for BNCT dose computation</li> </ul>
>>ROI Elemental Composition Sequence	{3006,00B6}	SQ	1C	
>>>ROI Elemental Composition Atomic Number	{3006,00B7}	US	1	
>>>ROI Elemental Composition Atomic Mass Fraction	{3006,00B8}	FL	1	
>>ROI Physical Property Value	{3006,00B4}	DS	1	
>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0

>Tissue Name	{4001,1010}	ST	3	RaySearch Private Tag. Contains the tissue name given to this Region of Interest.
>RBE Cell Type Name	{4001,1014}	ST	3	RaySearch Private Tag. Contains the RBE Cell Type Name given to this Region of Interest.

#### 7.1.4.8 Approval Module

Attribute name	Tag	Vr	Type	Comment
Approval Status	{300E,0002}	CS	1	Used to set the approval status of the imported structure set. Supported values: <ul style="list-style-type: none"> <li>• APPROVED - Structure set is approved.</li> <li>• UNAPPROVED - Structure set is not approved.</li> <li>• REJECTED - Same rules as for UNAPPROVED.</li> </ul>
Review Date	{300E,0004}	DA	2C	Stored if Approval Status is APPROVED.
Review Time	{300E,0005}	TM	2C	Stored if Approval Status is APPROVED.
Reviewer Name	{300E,0008}	PN	2C	Stored if Approval Status is APPROVED.

#### 7.1.4.9 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	{0008,0016}	UI	1	Supported value: 1.2.840.10008.5.1.4.1.1.481.3.
SOP Instance UID	{0008,0018}	UI	1	Stored internally to find references from other objects.
Specific Character Set	{0008,0005}	CS	1C	Supported values: ISO_IR 100, ISO_IR 192.
Instance Creation Date	{0008,0012}	DA	3	Stored as last save date if Structure Set Date {3006,0008} is missing.
Instance Creation Time	{0008,0013}	TM	3	Stored as last save time if Structure Set Time {3006,00109} is missing.

#### 7.1.5 RT Plan IOD

Imported RT Plans are converted into internal representation. Exported plans will therefore not be identical to the imported object.

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	No
Series	RT Series Module	Yes
Frame of Reference	Frame of Reference Module	Yes
Equipment	General Equipment Module	Yes
Plan	RT General Plan Module	Yes
	RT Prescription Module	Yes
	RT Patient Setup Module	Yes
	RT Fraction Scheme Module	Yes

RT Beams Module	Yes
RT Brachy Application Setups Module	No
Approval Module	Yes
SOP Common Module	Yes

### 7.1.5.1 Patient Module

Patient attributes only used to match patient.

Attribute name	Tag	Vr	Type	Comment
Patient's Name	{0010,0010}	PN	2	
Patient ID	{0010,0020}	LO	2	
Patient's Birth Date	{0010,0030}	DA	2	
Patient's Sex	{0010,0040}	CS	2	Supported values: <ul style="list-style-type: none"> <li>• M - Male</li> <li>• F - Female</li> <li>• 0 - Other</li> </ul>

### 7.1.5.2 RT Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	{0008,0060}	CS	1	Supported value: RTPLAN.
Series Instance UID	{0020,000E}	UI	1	Only read for display during import. Not stored.

### 7.1.5.3 Frame of Reference Module

Attribute name	Tag	Vr	Type	Comment
Frame of Reference UID	{0020,0052}	UI	1	Used to match referenced image set. Not stored.

### 7.1.5.4 General Equipment Module

Attribute name	Tag	Vr	Type	Comment
Manufacturer	{0008,0070}	LO	2	Used to handle special RT Plan cases. Not stored.
Manufacturer's Model Name	{0008,1090}	LO	3	Used to handle special RT Plan cases. Not stored.

### 7.1.5.5 RT General Plan Module

Attribute name	Tag	Vr	Type	Comment
RT Plan Label	{300A,0002}	SH	1	Stored as Beam Set Name.
RT Plan Name	{300A,0003}	LO	3	Stored as Treatment Plan Name.
RT Plan Description	{300A,0004}	ST	3	Stored as Beam Set Description.
RT Plan Date	{300A,0006}	DA	2	Stored as last save date.
RT Plan Time	{300A,0007}	TM	2	Stored as last save time.

RT Plan Geometry	{300A,000C}	CS	1	The plan is only required to reference structure set if value is PATIENT. Supported values: PATIENT, TREATMENT_DEVICE.
Referenced Structure Set Sequence	{300C,0060}	SQ	1C	If populated, the referenced Structure Set must be found. Not required.
>Referenced SOP Class UID	{0008,1150}	UI	1	Supported value: 1.2.840.10008.5.1.4.1.1.481.3.
>Referenced SOP Instance UID	{0008,1155}	UI	1	

#### 7.1.5.6 RT Prescription Module

Attribute name	Tag	Vr	Type	Comment
Dose Reference Sequence	{300A,0010}	SQ	3	Primary Prescription supported only. Evaluation functions are read from Brainlab iPlan RT Plan.
>Dose Reference Number	{300A,0012}	IS	1	Not stored.
>Dose Reference Structure Type	{300A,0014}	CS	1	Supported values: POINT, VOLUME, SITE.
>Dose Reference Description	{300A,0016}	LO	3	Stored as Site Description if Dose Reference Structure Type {300A,0014} is SITE.
>Referenced ROI Number	{3006,0084}	IS	1C	Used to find ROI when Dose Reference Structure Type {300A,0014} is VOLUME.
>Nominal Prior Dose	{300A,001A}	DS	3	Not supported. Supported value: 0.
>Dose Reference Type	{300A,0020}	CS	1	Only TARGET is supported for Primary Prescription. ORGAN_AT_RISK is supported when importing evaluation functions from Brainlab plans only. Supported values: TARGET, ORGAN_AT_RISK.
>Constraint Weight	{300A,0021}	DS	3	Used when importing evaluation functions from Brainlab plans only. Not stored.
>Target Minimum Dose	{300A,0025}	DS	3	
>Target Prescription Dose	{300A,0026}	DS	3	Used when importing Primary Prescription.
>Target Maximum Dose	{300A,0027}	DS	3	
>Target Underdose Volume Fraction	{300A,0028}	DS	3	Used when importing Primary Prescription.
>Organ at Risk Full-volume Dose	{300A,002A}	DS	3	Stored as evaluation functions from Brainlab plans only.
>Organ at Risk Limit Dose	{300A,002B}	DS	3	Stored as evaluation functions from Brainlab plans only.
>Organ at Risk Maximum Dose	{300A,002C}	DS	3	Stored as evaluation functions from Brainlab plans only.
>Organ at Risk Overdose Volume Fraction	{300A,002D}	DS	3	Stored as evaluation functions from Brainlab plans only.

#### 7.1.5.7 RT Patient Setup Module

Attribute name	Tag	Vr	Type	Comment
Patient Setup Sequence	{300A,0180}	SQ	1	Must contain at least one item.
>Patient Setup Number	{300A,0182}	IS	1	Not stored.

>Patient Position	{0018,5100}	CS	1C	Must be compatible with the planning image set patient position. Supported values: HFS, HFP, FFS, FFP, HFDR, HFDL, FFDR, FFDL.
>Motion Synchronization Sequence	{300A,0410}	SQ	3	
>>Respiratory Motion Compensation Technique	{0018,9170}	CS	1	Supported values: <ul style="list-style-type: none"> <li>• NONE</li> <li>• GATING</li> <li>• TRACKING - Static target tracking</li> <li>• REALTIME - Realtime tracking</li> </ul>
>>Respiratory Signal Source	{0018,9171}	CS	1	Supported values: <ul style="list-style-type: none"> <li>• NONE</li> <li>• EXTERNAL_MARKER</li> <li>• SPIROMETER</li> <li>• BELT</li> <li>• INTERNAL_MARKER</li> <li>• FIDUCIAL - Tracking of fiducial markers</li> <li>• TRACKING_ROI - Tracking of volume</li> </ul>
>>Respiratory Motion Compensation Technique Description	{0018,9185}	ST	3	
>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>Tracking ROI Sequence	{4001,1062}	SQ	3	Contains referenced ROI used for tracking
>>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>>Referenced ROI Number	{4001,1051}	IS	1C	
>Imaging Angles	{4001,1063}	DS	3	Selected imaging angles for Tomo Synchrony
medPhoton Private Creator	{30BB,0010}	LO	3	medPhoton 1.0
Patient Setup ID	{30BB,1000}	SH	3	Stored in Patient Setup as IDCAS parameters.
Imaging Protocol ID	{30BB,1001}	SH	3	Stored in Patient Setup as IDCAS parameters.

#### 7.1.5.8 RT Fraction Scheme Module

Attribute name	Tag	Vr	Type	Comment
Fraction Group Sequence	{300A,0070}	SQ	1	Each Fraction Group is imported as a separate Beam Set. Empty sequence is handled as special zero-fraction case. Dose cannot be imported to a plan with multiple Beam Sets.
>Fraction Group Number	{300A,0071}	IS	1	Used to identify fraction groups. Stored but only used internally.
>Number of Fractions Planned	{300A,0078}	IS	2	Used to create fraction schedule. If this value is 0 it is interpreted as a one-fraction plan.
>Number of Fraction Pattern Digits Per Day	{300A,0079}	IS	3	Not stored.
>Repeat Fraction Cycle Length	{300A,007A}	IS	3	Not stored.
>Fraction Pattern	{300A,007B}	LT	3	Not stored.

>Number of Beams	(300A,0080)	IS	1	
>Referenced Beam Sequence	(300C,0004)	SQ	1C	There must be at least one beam in the fraction group.
>>Referenced Beam Number	(300C,0006)	IS	1	
>>Beam Dose	(300A,0084)	DS	3	Stored as Beam Dose.
>>Beam Meterset	(300A,0086)	DS	3	Stored as Beam MU. Must be defined in MU.
>>Beam Dose Point Depth	(300A,0088)	FL	3	Stored as Physical Depth.
>>Beam Dose Point Equivalent Depth	(300A,0089)	FL	3	Stored as Radiological Depth.
>>Beam Dose Specification Point	(300A,0082)	DS	3	Stored as Dose Specification Point.
>Number of Brachy Application Setups	(300A,00A0)	IS	1	Not supported. Supported value: 0.

### 7.1.5.9 RT Beams Module

Attribute name	Tag	Vr	Type	Comment
Beam Sequence	(300A,00B0)	SQ	1	
>Beam Number	(300A,00C0)	IS	1	Stored as Beam Number.
>Beam Name	(300A,00C2)	LO	3	Stored as Beam Name.
>Beam Description	(300A,00C3)	ST	3	Stored as Beam Description.
>Beam Type	(300A,00C4)	CS	1	Used to identify plan technique. Beam Type STATIC is converted to SMLC plan. Beam Type DYNAMIC which have Gantry Rotation Direction != NONE are converted to Arc plans Otherwise plans is considered to be DMLC plan. Supported values: STATIC, DYNAMIC.
>Radiation Type	(300A,00C6)	CS	2	Supported values: PHOTON, ELECTRON, NEUTRON.
>Primary Fluence Mode Sequence	(3002,0050)	SQ	3	
>>Fluence Mode	(3002,0051)	CS	1	
>>Fluence Mode ID	(3002,0052)	SH	1C	
>Treatment Machine Name	(300A,00B2)	SH	2	Stored as Machine Reference. Used to find treatment machine.
>Primary Dosimeter Unit	(300A,00B3)	CS	3	Used for precondition check only. Supported values: MU, MINUTE.
>Source-Axis Distance	(300A,00B4)	DS	3	Used to calculate Block isocenter to tray distance.
>Beam Limiting Device Sequence	(300A,00B6)	SQ	1	
>>RT Beam Limiting Device Type	(300A,00B8)	CS	1	Supported values: X, Y, ASYMX, ASYMY, MLCX, MLCY.
>>Number of Leaf/Jaw Pairs	(300A,00BC)	IS	1	
>>Leaf Position Boundaries	(300A,00BE)	DS	2C	
>Referenced Patient Setup Number	(300C,006A)	IS	3	If this attribute is not set, the patient setup sequence must contain one patient setup. If multiple patient setups are defined in the patient setup sequence, this attribute must be set. All referenced patient setups must use the same patient position.

>Treatment Delivery Type	{300A,00CE}	CS	3	If the Treatment Delivery type is null, beam is assumed to be of type TREATMENT. Supported values: TREATMENT, SETUP.
>Number of Wedges	{300A,00D0}	IS	1	
>Wedge Sequence	{300A,00D1}	SQ	1C	
>>Wedge Number	{300A,00D2}	IS	1	Not stored.
>>Wedge Type	{300A,00D3}	CS	2	Converted to internal Wedge Type. Supported values: STANDARD, DYNAMIC, MOTORIZED.
>>Wedge ID	{300A,00D4}	SH	3	Stored as Wedge ID.
>>Wedge Angle	{300A,00D5}	IS	2	Stored as Wedge Angle.
>>Wedge Orientation	{300A,00D8}	DS	2	Stored as Wedge Orientation. Supported values: 0.0, 90.0, 180.0, 270.0.
>Number of Compensators	{300A,00E0}	IS	1	
>Compensator Sequence	{300A,00E3}	SQ	1C	Compensator is not supported. Varian fluences stored as compensators are ignored.
>>Compensator Number	{300A,00E4}	IS	1C	
>>Compensator Type	{300A,00EE}	CS	3	
>>Compensator Rows	{300A,00E7}	IS	1	
>>Compensator Columns	{300A,00E8}	IS	1	
>>Compensator Pixel Spacing	{300A,00E9}	DS	1	
>>Compensator Position	{300A,00EA}	DS	1	
>Number of Boli	{300A,00ED}	IS	1	
>Referenced Bolus Sequence	{300C,00B0}	SQ	1C	Referenced boli must exist.
>>Referenced ROI Number	{3006,0084}	IS	1	
>>Bolus ID	{300A,00DC}	SH	3	
>Number of Blocks	{300A,00F0}	IS	1	Electron plans support 0 to 1 blocks. Photon plans support 0 to 2 blocks.
>Total Block Tray Factor	{300A,00F2}	DS	3	Must be set to between 0.5 and 1.0 for photon plans if block exist in the plan. Not used for electron plans
>Block Sequence	{300A,00F4}	SQ	1C	
>>Block Tray ID	{300A,00F5}	SH	3	Stored as Block Tray ID.
>>Source to Block Tray Distance	{300A,00F6}	DS	2	Used to calculate Block isocenter to tray distance.
>>Block Type	{300A,00F8}	CS	1	Stored as Block Type. Supported values: SHIELDING, APERTURE.
>>Block Divergence	{300A,00FA}	CS	2	Stored as Block Divergence. Electron blocks will convert PRESENT to ABSENT because of model restrictions. Supported values: PRESENT, ABSENT.
>>Block Mounting Position	{300A,00FB}	CS	3	Stored as Block Mounting Position. Electron blocks will convert PATIENT_SIDE to SOURCE_SIDE because of model restrictions. Supported values: PATIENT_SIDE, SOURCE_SIDE.

>>Block Number	{300A,00FC}	IS	1	Stored as Block Number.
>>Block Name	{300A,00FE}	LO	3	Stored as Block Name.
>>Material ID	{300A,00E1}	SH	2	Stored as Block Material Name.
>>Block Thickness	{300A,0100}	DS	2C	Not supported for Photon blocks. Stored as Block Thickness for Electron blocks.
>>Block Transmission	{300A,0102}	DS	2C	Stored for Photon blocks only. Must be between 0.0 - 0.2.
>>Block Number of Points	{300A,0104}	IS	2	
>>Block Data	{300A,0106}	DS	2	Stored as Block Data.
>Applicator Sequence	{300A,0107}	SQ	3	Stored as Electron applicator or Photon cone depending on Radiation Type.
>>Applicator ID	{300A,0108}	SH	1	Stored as Applicator ID for Electron applicator. Must be "Apex 2.5" for Photon cone.
>>Applicator Type	{300A,0109}	CS	1	Stored as Applicator Type. Supported values: ELECTRON_SQUARE, ELECTRON_RECT, ELECTRON_CIRC, ELECTRON_SHORT, ELECTRON_OPEN, PHOTON_CIRC, STEREOTACTIC.
>>Applicator Geometry Sequence	{300A,0431}	SQ	3	Only used for photon cones.
>>>Applicator Aperture Shape	{300A,0432}	CS	1	Not stored. Supported value: SYM_CIRCULAR.
>>>Applicator Opening	{300A,0433}	FL	1C	Stored as Applicator Opening (converted from MM to CM).
>>Applicator Description	{300A,010A}	LO	3	
>>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>>Applicator Thickness	{4001,1054}	DS	3	Physical thickness of applicator (in mm) parallel to radiation beam axis.
>Final Cumulative Meterset Weight	{300A,010E}	DS	1C	
>Number of Control Points	{300A,0110}	IS	1	Must be greater than or equal to 2. Must be exactly 2 if Radiation Type {300A, 00C6} is ELECTRON.
>Control Point Sequence	{300A,0111}	SQ	1	
>>Control Point Index	{300A,0112}	IS	1	
>>Cumulative Meterset Weight	{300A,0134}	DS	2	
>>Nominal Beam Energy	{300A,0114}	DS	3	Required for the first control point and is assumed to be constant throughout beam. Nominal Beam Energy will be mapped against nominal energy of a machine beam quality when assigning imported beams to a machine. The nominal energy can differ from the energy used internally for dose calculation.
>>Dose Rate Set	{300A,0115}	DS	3	
>>Wedge Position Sequence	{300A,0116}	SQ	1C	
>>>Referenced Wedge Number	{300C,00C0}	IS	1	
>>>Wedge Position	{300A,0118}	CS	1	Supported values: IN, OUT.
>>Beam Limiting Device Position Sequence	{300A,011A}	SQ	1C	Leaf/jaw positions are required to be constant within a control point pair for an SMLC plan.



>>>RT Beam Limiting Device Type	{300A,00B8}	CS	1	Supported values: X, Y, ASYMX, ASYMY, MLCX, MLCY.
>>>Leaf/Jaw Positions	{300A,011C}	DS	1	
>>Gantry Angle	{300A,011E}	DS	1C	Required for the first control point and required to not change for SMLC and DMLC. See Beam Type {300A,00C4} for more information.
>>Gantry Rotation Direction	{300A,011F}	CS	1C	Required for the first control point and required to not be "NONE" for SMLC and DMLC. See Beam Type {300A,00C4} for more information. Supported values: NONE, CW, CC.
>>Gantry Pitch Angle	{300A,014A}	FL	3	Must be null or 0.
>>Gantry Pitch Rotation Direction	{300A,014C}	CS	3	Supported values: NONE, CW, CC.
>>Beam Limiting Device Angle	{300A,0120}	DS	1C	Required for the first control point, must be empty for subsequent control points or have the same value as in the first control point.
>>Beam Limiting Device Rotation Direction	{300A,0121}	CS	1C	Required for the first control point. Supported values: NONE, CW, CC.
>>Patient Support Angle	{300A,0122}	DS	1C	Required for the first control point, must be empty for subsequent control points or have the same value as in the first control point.
>>Patient Support Rotation Direction	{300A,0123}	CS	1C	Required for the first control point. Supported values: NONE, CW, CC.
>>Table Top Eccentric Axis Distance	{300A,0124}	DS	3	
>>Table Top Eccentric Angle	{300A,0125}	DS	1C	Required for the first control point. Supported value: 0.0.
>>Table Top Eccentric Rotation Direction	{300A,0126}	CS	1C	Required for the first control point. Supported values: NONE, CW, CC.
>>Table Top Pitch Angle	{300A,0140}	FL	1C	
>>Table Top Pitch Rotation Direction	{300A,0142}	CS	1C	Supported values: NONE, CW, CC.
>>Table Top Roll Angle	{300A,0144}	FL	1C	
>>Table Top Roll Rotation Direction	{300A,0146}	CS	1C	Supported values: NONE, CW, CC.
>>Table Top Vertical Position	{300A,0128}	DS	2C	
>>Table Top Longitudinal Position	{300A,0129}	DS	2C	
>>Table Top Lateral Position	{300A,012A}	DS	2C	
>>Isocenter Position	{300A,012C}	DS	2C	Required for the first control point. Must be constant throughout beam.
>>TomoTherapy HA Private Creator	{300D,0010}	LO	3	TOMO_HA_01
>>Tomo Projection Sinogram Data	{300D,10A7}	DS	3	Projection sinogram values for all 64 leaves of projection starting at the current control point. Empty if all 64 leaves have a sinogram value of zero. Empty for last control point.
>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>Treatment Machine Commission Time	{4001,1001}	DT	3	The commission time of the treatment machine. Used together with Treatment Machine Name to determine the correct machine.
>Tomo IDMS Beam ID	{4001,1021}	LO	3	

>Tomo IDMS Machine ID	{4001,1020}	LO	3	
>Tomo Beam Isocenter	{4001,1023}	DS	3	
>Tomo machine revision	{4001,1025}	IS	3	
>Tomo beam revision	{4001,1026}	IS	3	
>Tomo intended back jaw position	{4001,1027}	DS	3	
>Tomo intended front jaw position	{4001,1028}	DS	3	
>Tomo Projection Time	{4001,1053}	DS	3	Projection time for Tomo Helical and Tomo Direct plans. Used for import of legacy plans.
>Planned Blood Boron Concentration	{4001,1015}	DS	3	
>Tomo Projection Time FL	{4001,1068}	FL	3	Projection time for Tomo Helical and Tomo Direct plans.
>TomoTherapy HA Private Creator	{300D,0010}	LO	3	TOMO_HA_01
>Tomo Gantry Period	{300D,1040}	DS	3	Used for import of TomoHelical plans exported from RayStation 6 and 7.

#### 7.1.5.10 Approval Module

Attribute name	Tag	Vr	Type	Comment
Approval Status	{300E,0002}	CS	1	Used to set the approval status of the imported structure set. Supported values: <ul style="list-style-type: none"> <li>APPROVED - Plan is approved.</li> <li>UNAPPROVED - Plan is not approved.</li> <li>REJECTED - Same rules as for UNAPPROVED.</li> </ul>
Review Date	{300E,0004}	DA	2C	Stored if Approval Status is APPROVED.
Review Time	{300E,0005}	TM	2C	Stored if Approval Status is APPROVED.
Reviewer Name	{300E,0008}	PN	2C	Stored if Approval Status is APPROVED.

#### 7.1.5.11 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	{0008,0016}	UI	1	Supported value: 1.2.840.10008.5.1.4.1.1.481.5.
SOP Instance UID	{0008,0018}	UI	1	Stored internally to find references from other objects.
Specific Character Set	{0008,0005}	CS	1C	Supported values: ISO_IR 100, ISO_IR 192.
Instance Creation Date	{0008,0012}	DA	3	Stored if last save date if RT Plan Date {300A, 0006} is missing.
Instance Creation Time	{0008,0013}	TM	3	Stored if last save time if RT Plan Time {300A, 0007} is missing.
RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
Tomo Plan SOP Instance UID	{4001,102A}	UI	3	The SOP Instance UID that shall be used by RayGateway for creating plan instances in IDMS based on this RT Plan.

## 7.1.6 PET Image IOD

Imported images are stored internally. Filters modifies data permanently.

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	Yes
	Patient Study Module	Yes
	Clinical Trial Study Module	Yes
Series	General Series Module	Yes
	PET Series Module	Yes
	PET Isotope Module	Yes
	NM/PET Patient Orientation Module	Yes
Frame of Reference	Frame of Reference Module	Yes
	Synchronization Module	No
Equipment	General Equipment Module	Yes
Image	General Image Module	No
	Image Plane Module	Yes
	Image Pixel Module	Yes
	PET Image Module	Yes
	SOP Common Module	Yes

### 7.1.6.1 Patient Module

Patient attributes are read and stored when importing a new patient.

Patient ID and the first component (Last Name) of Patient's Name is the unique identifier. Differences in Patient's Sex and Patient's Birth Date are checked.

Attribute name	Tag	Vr	Type	Comment
Patient's Name	{0010,0010}	PN	2	
Patient ID	{0010,0020}	LO	2	
Patient's Birth Date	{0010,0030}	DA	2	
Patient's Sex	{0010,0040}	CS	2	Empty value defaults to Unknown gender. Supported values: <ul style="list-style-type: none"> <li>• M - Male</li> <li>• F - Female</li> <li>• 0 - Other</li> </ul>
Patient Identity Removed	{0012,0062}	CS	3	Supported value: YES.
De-identification Method	{0012,0063}	LO	1C	

### 7.1.6.2 General Study Module

Attribute name	Tag	Vr	Type	Comment
----------------	-----	----	------	---------

Study Instance UID	{0020,000D}	UI	1	Stored when importing a new case but not used.
Study Date	{0008,0020}	DA	2	
Study Time	{0008,0030}	TM	2	
Referring Physician's Name	{0008,0090}	PN	2	Stored as Case Physician when importing a new case.
Study ID	{0020,0010}	SH	2	
Accession Number	{0008,0050}	SH	2	

### 7.1.6.3 Patient Study Module

Attribute name	Tag	Vr	Type	Comment
Patient's Weight	{0010,1030}	DS	3	Used for SUV calculation.

### 7.1.6.4 Clinical Trial Study Module

Attribute name	Tag	Vr	Type	Comment
Clinical Trial Time Point ID	{0012,0050}	LO	2	
Consent for Clinical Trial Use Sequence	{0012,0083}	SQ	3	Only used for display of contract information. Mainly for internal use.
>Distribution Type	{0012,0084}	CS	1C	
>Clinical Trial Protocol ID	{0012,0020}	LO	1C	
>Consent for Distribution Flag	{0012,0085}	CS	1	Supported values: YES, NO, WITHDRAWN.
>Clinical Trial Subject ID	{0012,0040}	LO	3	
>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>Consent Expiry Date	{4001,1067}	DA	3	Date when consent expires.

### 7.1.6.5 General Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	{0008,0060}	CS	1	Supported value: PT.
Series Instance UID	{0020,000E}	UI	1	Stored to verify references from RT Structure Sets.
Series Number	{0020,0011}	IS	2	
Series Date	{0008,0021}	DA	3	
Series Time	{0008,0031}	TM	3	
Protocol Name	{0018,1030}	LO	3	Used to identify the protocol during import.
Series Description	{0008,103E}	LO	3	Only used for display.
Patient Position	{0018,5100}	CS	2C	Supported values: HFS, HFP, FFS, FFP.

### 7.1.6.6 PET Series Module

Attribute name	Tag	Vr	Type	Comment
Series Date	{0008,0021}	DA	1	
Series Time	{0008,0031}	TM	1	

Units	{0054,1001}	CS	1	Supported values: CNTS, NONE, CM2, PCNT, CPS, BQML, MGINML, UMOLMINML, MLMING, MLG, 1CM, UMOLML, PROPCNTS, PROPCPS, MLMINML, MLML, GML, STDDEV.
Counts Source	{0054,1002}	CS	1	
Series Type	{0054,1000}	CS	1	
Number of Slices	{0054,0081}	US	1	
Corrected Image	{0028,0051}	CS	2	Must contain DECYATTN in order to compute SUV.
Decay Correction	{0054,1102}	CS	1	Supported values: <ul style="list-style-type: none"> <li>• NONE</li> <li>• START - Required to compute SUV</li> <li>• ADMIN</li> </ul>

### 7.1.6.7 PET Isotope Module

Attribute name	Tag	Vr	Type	Comment
Radiopharmaceutical Information Sequence	{0054,0016}	SQ	2	
>Radiopharmaceutical Start Time	{0018,1072}	TM	3	
>Radionuclide Total Dose	{0018,1074}	DS	3	
>Radionuclide Half Life	{0018,1075}	DS	3	

### 7.1.6.8 NM/PET Patient Orientation Module

Attribute name	Tag	Vr	Type	Comment
Patient Orientation Code Sequence	{0054,0410}	SQ	2	
>Code Value	{0008,0100}	SH	1C	Supported value: <ul style="list-style-type: none"> <li>• F-10450 - recumbent</li> </ul>
>Coding Scheme Designator	{0008,0102}	SH	1C	Supported value: 99SDM.
>Code Meaning	{0008,0104}	LO	1	Supported value: recumbent.
>Patient Orientation Modifier Code Sequence	{0054,0412}	SQ	2C	
>>Code Value	{0008,0100}	SH	1C	Supported values: <ul style="list-style-type: none"> <li>• F-10340 - supine</li> <li>• F-10310 - prone</li> </ul>
>>Coding Scheme Designator	{0008,0102}	SH	1C	Supported value: 99SDM.
>>Code Meaning	{0008,0104}	LO	1	Supported values: supine, prone.
Patient Gantry Relationship Code Sequence	{0054,0414}	SQ	2	

>Code Value	{0008,0100}	SH	1C	Supported values: <ul style="list-style-type: none"> <li>• F-10470 - headfirst</li> <li>• F-10480 - feet-first</li> <li>• G-5190 - headfirst</li> <li>• G-5191 - feet-first</li> </ul>
>Coding Scheme Designator	{0008,0102}	SH	1C	Supported value: 99SDM.
>Code Meaning	{0008,0104}	LO	1	Supported values: headfirst, feet-first.

#### 7.1.6.9 Frame of Reference Module

Attribute name	Tag	Vr	Type	Comment
Frame of Reference UID	{0020,0052}	UI	1	Stored internally to be used if referenced from other dataset.

#### 7.1.6.10 General Equipment Module

Attribute name	Tag	Vr	Type	Comment
Pixel Padding Value	{0028,0120}	US or SS	1C	If pixel padding is present within the region used for dose computation, the resulting dose may be wrong.

#### 7.1.6.11 Image Plane Module

Attribute name	Tag	Vr	Type	Comment
Pixel Spacing	{0028,0030}	DS	1	Cannot be greater than 5 mm.
Image Orientation (Patient)	{0020,0037}	DS	1	Image orientation must be transversal.
Image Position (Patient)	{0020,0032}	DS	1	Converted to internal image corner. Used to find slice direction, which must be constant throughout the series and along the Z-axis.

#### 7.1.6.12 Image Pixel Module

Attribute name	Tag	Vr	Type	Comment
Samples per Pixel	{0028,0002}	US	1	Supported value: 1.
Photometric Interpretation	{0028,0004}	CS	1	Supported value: MONOCHROME2.
Rows	{0028,0010}	US	1	
Columns	{0028,0011}	US	1	
Bits Allocated	{0028,0100}	US	1	Supported value: 16.
Bits Stored	{0028,0101}	US	1	Must be equal to or less than Bits Allocated.
High Bit	{0028,0102}	US	1	
Pixel Representation	{0028,0103}	US	1	Supported values: 0, 1.
Pixel Data	{7FE0,0010}	OB or OW	1C	
Pixel Padding Range Limit	{0028,0121}	US or SS	1C	

## 7.1.6.13 PET Image Module

Attribute name	Tag	Vr	Type	Comment
Image Type	{0008,0008}	CS	1	
Samples per Pixel	{0028,0002}	US	1	Supported value: 1.
Photometric Interpretation	{0028,0004}	CS	1	Supported value: MONOCHROME2.
Bits Allocated	{0028,0100}	US	1	Supported value: 16.
Bits Stored	{0028,0101}	US	1	Must be greater than 0 and less than or equal to Bits Allocated.
High Bit	{0028,0102}	US	1	HighBit must be equal to Bits Stored - 1.
Rescale Intercept	{0028,1052}	DS	1	Supported value: 0.0.
Rescale Slope	{0028,1053}	DS	1	
Frame Reference Time	{0054,1300}	DS	1	
Image Index	{0054,1330}	US	1	
Acquisition Date	{0008,0022}	DA	2	
Acquisition Time	{0008,0032}	TM	2	
Actual Frame Duration	{0018,1242}	IS	2	
Decay Factor	{0054,1321}	DS	1C	
Philips Private Creator	{7053,0010}	LO	3	Philips PET Private Group
Philips SUV Scale Factor	{7053,1000}	DS	3	
Philips Activity Concentration Scale Factor	{7053,1009}	DS	3	
GE Private Creator	{0009,0010}	LO	3	GEMS_PETD_01
GE PET Scan Time	{0009,1000}	DT	3	

## 7.1.6.14 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	{0008,0016}	UI	1	Supported value: 1.2.840.10008.5.1.4.1.1.128.
SOP Instance UID	{0008,0018}	UI	1	Stored internally to be used if referenced from other dataset.
Specific Character Set	{0008,0005}	CS	1C	Supported values: ISO_IR 100, ISO_IR 192.
Instance Creation Date	{0008,0012}	DA	3	
Instance Creation Time	{0008,0013}	TM	3	
Contributing Equipment Sequence	{0018,A001}	SQ	3	Only used for display of contract information. Mainly for internal use.
>Purpose of Reference Code Sequence	{0040,A170}	SQ	1	
>>Code Value	{0008,0100}	SH	1C	Supported value: <ul style="list-style-type: none"> <li>109103 - Modifying Equipment</li> </ul>

>>Coding Scheme Designator	{0008,0102}	SH	1C	Supported value: DCM.
>>Code Meaning	{0008,0104}	LO	1	Supported value: Modifying Equipment.
>Manufacturer	{0008,0070}	LO	1	
>Manufacturer's Model Name	{0008,1090}	LO	3	
>Software Versions	{0018,1020}	LO	3	
>Contribution DateTime	{0018,A002}	DT	3	
>Contribution Description	{0018,A003}	ST	3	

### 7.1.7 Spatial Registration IOD

Imported Spatial Registrations are always converted into internal representation. Exported registrations will therefore not be identical to the imported object.

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	Yes
Series	General Series Module	Yes
	Spatial Registration Series Module	Yes
Frame of Reference	Frame of Reference Module	Yes
Equipment	General Equipment Module	Yes
Spatial Registration	Spatial Registration Module	Yes
	Common Instance Reference Module	Yes
	SOP Common Module	Yes

#### 7.1.7.1 Patient Module

Patient attributes only used to match patient.

Attribute name	Tag	Vr	Type	Comment
Patient's Name	{0010,0010}	PN	2	Only used for matching patient.
Patient ID	{0010,0020}	LO	2	Only used for matching patient.
Patient's Birth Date	{0010,0030}	DA	2	Only used for matching patient.
Patient's Sex	{0010,0040}	CS	2	Only used for matching patient. Supported values: <ul style="list-style-type: none"> <li>• M - Male</li> <li>• F - Female</li> <li>• 0 - Other</li> </ul>

#### 7.1.7.2 General Study Module

Attribute name	Tag	Vr	Type	Comment
Study Instance UID	{0020,0000}	UI	1	



## 7.1.7.3 General Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	{0008,0060}	CS	1	Supported value: REG.
Series Instance UID	{0020,000E}	UI	1	Only read for display during import. Not stored.
Series Number	{0020,0011}	IS	2	
Series Date	{0008,0021}	DA	3	
Series Time	{0008,0031}	TM	3	
Protocol Name	{0018,1030}	LO	3	
Series Description	{0008,103E}	LO	3	
Patient Position	{0018,5100}	CS	2C	

## 7.1.7.4 Spatial Registration Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	{0008,0060}	CS	1	Supported value: REG.

## 7.1.7.5 Frame of Reference Module

Attribute name	Tag	Vr	Type	Comment
Frame of Reference UID	{0020,0052}	UI	1	Used to match referenced image set.

## 7.1.7.6 General Equipment Module

Attribute name	Tag	Vr	Type	Comment
----------------	-----	----	------	---------

## 7.1.7.7 Spatial Registration Module

Attribute name	Tag	Vr	Type	Comment
Content Date	{0008,0023}	DA	1	
Content Time	{0008,0033}	TM	1	
Instance Number	{0020,0013}	IS	1	
Content Label	{0070,0080}	CS	1	
Content Description	{0070,0081}	LO	2	
Content Creator's Name	{0070,0084}	PN	2	
Registration Sequence	{0070,0308}	SQ	1	
>Frame of Reference UID	{0020,0052}	UI	1C	Used to match referenced image sets.
>Referenced Image Sequence	{0008,1140}	SQ	1C	
>>Referenced SOP Class UID	{0008,1150}	UI	1	
>>Referenced SOP Instance UID	{0008,1155}	UI	1	
>Matrix Registration Sequence	{0070,0309}	SQ	1	
>>Matrix Sequence	{0070,030A}	SQ	1	

>>>Frame of Reference Transformation Matrix	(3006,00C6)	DS	1	
>>>Frame of Reference Transformation Matrix Type	(0070,030C)	CS	1	Supported value: RIGID.

#### 7.1.7.8 Common Instance Reference Module

Attribute name	Tag	Vr	Type	Comment
Referenced Series Sequence	(0008,1115)	SQ	1C	
>Series Instance UID	(0020,000E)	UI	1	
>Referenced Instance Sequence	(0008,114A)	SQ	1	
>>Referenced SOP Class UID	(0008,1150)	UI	1	
>>Referenced SOP Instance UID	(0008,1155)	UI	1	
Studies Containing Other Referenced Instances Sequence	(0008,1200)	SQ	1C	If referenced images are located in a Study different from the REG object.
>Study Instance UID	(0020,000D)	UI	1	
>Referenced Series Sequence	(0008,1115)	SQ	1	
>>Series Instance UID	(0020,000E)	UI	1	
>>Referenced Instance Sequence	(0008,114A)	SQ	1	
>>>Referenced SOP Class UID	(0008,1150)	UI	1	
>>>Referenced SOP Instance UID	(0008,1155)	UI	1	

#### 7.1.7.9 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	(0008,0016)	UI	1	Supported value: 1.2.840.10008.5.1.4.1.1.66.1.
SOP Instance UID	(0008,0018)	UI	1	Stored internally to references from other objects.
Specific Character Set	(0008,0005)	CS	1C	Supported values: ISO_IR 100, ISO_IR 192.
Instance Creation Date	(0008,0012)	DA	3	Stored as last save date.
Instance Creation Time	(0008,0013)	TM	3	Stored as last save time.

#### 7.1.8 RT Ion Plan IOD

Imported RT Ion Plans are converted into internal representation. Exported plans will therefore not be identical to the imported object.

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	No
Series	RT Series Module	Yes
Frame of Reference	Frame of Reference Module	Yes
Equipment	General Equipment Module	Yes
Plan	RT General Plan Module	Yes

RT Prescription Module	Yes
RT Patient Setup Module	Yes
RT Fraction Scheme Module	Yes
RT Ion Beams Module	Yes
Approval Module	Yes
SOP Common Module	Yes

### 7.1.8.1 Patient Module

Patient attributes only used to match patient.

Attribute name	Tag	Vr	Type	Comment
Patient's Name	(0010,0010)	PN	2	
Patient ID	(0010,0020)	LO	2	
Patient's Birth Date	(0010,0030)	DA	2	
Patient's Sex	(0010,0040)	CS	2	Supported values: <ul style="list-style-type: none"> <li>• M - Male</li> <li>• F - Female</li> <li>• O - Other</li> </ul>

### 7.1.8.2 RT Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	(0008,0060)	CS	1	Supported value: RTPLAN.
Series Instance UID	(0020,000E)	UI	1	Only read for display during import. Not stored.

### 7.1.8.3 Frame of Reference Module

Attribute name	Tag	Vr	Type	Comment
Frame of Reference UID	(0020,0052)	UI	1	Used to match referenced image set. Not stored.

### 7.1.8.4 General Equipment Module

Attribute name	Tag	Vr	Type	Comment
Manufacturer	(0008,0070)	LO	2	Used to handle special RT Plan cases. Not stored.
Manufacturer's Model Name	(0008,1090)	LO	3	Used to handle special RT Plan cases. Not stored.

### 7.1.8.5 RT General Plan Module

Attribute name	Tag	Vr	Type	Comment
RT Plan Label	(300A,0002)	SH	1	Stored as Beam Set Name.
RT Plan Name	(300A,0003)	LO	3	Stored as Treatment Plan Name.
RT Plan Description	(300A,0004)	ST	3	Stored as Beam Set Description.
RT Plan Date	(300A,0006)	DA	2	Stored as last save date.

RT Plan Time	{300A,0007}	TM	2	Stored as last save time.
Plan Intent	{300A,000A}	CS	3	
RT Plan Geometry	{300A,000C}	CS	1	The plan is only required to reference structure set if value is PATIENT. Supported values: PATIENT, TREATMENT_DEVICE.
Referenced Structure Set Sequence	{300C,0060}	SQ	1C	
>Referenced SOP Class UID	{0008,1150}	UI	1	Supported value: 1.2.840.10008.5.1.4.1.1.481.3.
>Referenced SOP Instance UID	{0008,1155}	UI	1	Referenced Structure Set must exist or be imported together with plan.

#### 7.1.8.6 RT Prescription Module

Attribute name	Tag	Vr	Type	Comment
Dose Reference Sequence	{300A,0010}	SQ	3	Primary Prescription supported only.
>Dose Reference Number	{300A,0012}	IS	1	Not stored.
>Dose Reference Structure Type	{300A,0014}	CS	1	Supported values: POINT, VOLUME, SITE.
>Dose Reference Description	{300A,0016}	LO	3	Stored as Site Description if Dose Reference Structure Type {300A,0014} is SITE.
>Referenced ROI Number	{3006,0084}	IS	1C	Used to find referenced ROI if Dose Reference Structure Type {300A,0014} is VOLUME.
>Nominal Prior Dose	{300A,001A}	DS	3	Not supported. Supported value: 0.
>Dose Reference Type	{300A,0020}	CS	1	Supported value: TARGET.
>Target Prescription Dose	{300A,0026}	DS	3	Used when importing Primary Prescription. The value is interpreted as an RBE-weighted prescription dose.
>Target Underdose Volume Fraction	{300A,0028}	DS	3	Used when importing Primary Prescription.
>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>Target Prescription Effective Dose	{4001,1011}	DS	3	RaySearch Private Tag. Prescribed dose to Dose Reference if Dose Reference Type {300A,0020} is TARGET. The dose has been corrected for biological effect using user-defined modeling technique.

#### 7.1.8.7 RT Patient Setup Module

Attribute name	Tag	Vr	Type	Comment
Patient Setup Sequence	{300A,0180}	SQ	1	Must contain at least one item.
>Patient Setup Number	{300A,0182}	IS	1	Not stored.
>Patient Position	{0018,5100}	CS	1C	Must be compatible with the planning image set patient position. Supported values: HFS, HFP, FFS, FFP, HFDR, HFDL, FFDR, FFDL, SITTING.

>Motion Synchronization Sequence	{300A,0410}	SQ	3	Stored in Patient Setup Motion Synchronization Technique. Only the following combinations of Respiratory Motion Compensation Technique and Respiratory Signal Source are supported: [NONE, NONE], [GATING, EXTERNAL_MARKER], [GATING, SPIROMETER] and [GATING, BELT]. All other combinations will be disregarded.
>>Respiratory Motion Compensation Technique	{0018,9170}	CS	1	Supported values: NONE, GATING.
>>Respiratory Signal Source	{0018,9171}	CS	1	Supported values: NONE, EXTERNAL_MARKER, SPIROMETER, BELT, INTERNAL_MARKER.
>>Respiratory Interval Time	{0020,9254}	FD	3	Measured interval time in ms from maximum respiration peak to the next peak for the respiratory cycle in which this frame occurs.
>>IMPAC Private Creator	{300B,0010}	LO	3	IMPAC
>>Respiratory Phase Gating Duty Cycle	{300B,1020}	DS	3	Start and end phases for gating window for respiratory motion, given as percentages of the respiratory cycle.
medPhoton Private Creator	{30BB,0010}	LO	3	medPhoton 1.0
Patient Setup ID	{30BB,1000}	SH	3	Stored in Patient Setup as IDCAS parameters.
Imaging Protocol ID	{30BB,1001}	SH	3	Stored in Patient Setup as IDCAS parameters.

#### 7.1.8.8 RT Fraction Scheme Module

Attribute name	Tag	Vr	Type	Comment
Fraction Group Sequence	{300A,0070}	SQ	1	Each Fraction Group is imported as a separate Beam Set within the same Treatment Plan. Empty sequence is handled as special zero-fraction case. NOTE: Dose cannot be imported to a plan with multiple Beam Sets.
>Fraction Group Number	{300A,0071}	IS	1	Used to identify fraction groups. Stored but only used internally.
>Number of Fractions Planned	{300A,0078}	IS	2	Used to create fraction schedule. If this value is 0 it is interpreted as a one-fraction plan.
>Number of Fraction Pattern Digits Per Day	{300A,0079}	IS	3	Not stored.
>Repeat Fraction Cycle Length	{300A,007A}	IS	3	Not stored.
>Fraction Pattern	{300A,007B}	LT	3	Not stored.
>Number of Beams	{300A,0080}	IS	1	
>Referenced Beam Sequence	{300C,0004}	SQ	1C	There must be at least one beam in the fraction group.
>>Referenced Beam Number	{300C,0006}	IS	1	
>>Beam Dose	{300A,0084}	DS	3	Stored as Beam Dose.
>>Beam Dose Type	{300A,0090}	CS	1C	Supported values: PHYSICAL, EFFECTIVE.
>>Beam Meterset	{300A,0086}	DS	3	Stored as Beam MU. Must be defined in MU.
>>Beam Dose Specification Point	{300A,0082}	DS	3	Stored as Dose Specification Point.
>Number of Brachy Application Setups	{300A,00A0}	IS	1	Not supported. Supported value: 0.

## 7.1.8.9 RT Ion Beams Module

Attribute name	Tag	Vr	Type	Comment
Ion Beam Sequence	(300A,03A2)	SQ	1	
>Beam Number	(300A,00C0)	IS	1	Stored as Beam Number.
>Beam Name	(300A,00C2)	LO	1	Stored as Beam Name.
>Beam Description	(300A,00C3)	ST	3	Stored as Beam Description.
>Beam Type	(300A,00C4)	CS	1	Supported value: STATIC.
>Radiation Type	(300A,00C6)	CS	1	Supported values: PROTON, ION.
>Radiation Mass Number	(300A,0302)	IS	1C	Checked if Radiation Type (300A,00C6) is ION. Not stored. Supported values: 4, 12.
>Radiation Atomic Number	(300A,0304)	IS	1C	Checked if Radiation Type (300A,00C6) is ION. Not stored. Supported values: 2, 6.
>Radiation Charge State	(300A,0306)	SS	1C	Checked if Radiation Type (300A,00C6) is ION. Not stored. Supported values: 2, 6.
>Scan Mode	(300A,0308)	CS	1	If value is NONE, the private tag IBA Scattered Mode (300D, 1002) depicts if the beam delivery type is Single Scattering or Double Scattering. Supported values: NONE, UNIFORM, MODULATED, MODULATED_SPEC, LINE, WOBBLING.
>Modulated Scan Mode Type	(300A,0309)	CS	1C	Supported values: STATIONARY, LEAPING.
>Treatment Machine Name	(300A,00B2)	SH	2	Stored as Machine Reference. Used for finding treatment machine.
>Primary Dosimeter Unit	(300A,00B3)	CS	1	Supported values: MU, NP.
>Virtual Source-Axis Distances	(300A,030A)	FL	1	Equal values in X and Y if Beam Limiting Device Angle (300A, 0120) is not 0
>Ion Beam Limiting Device Sequence	(300A,03A4)	SQ	3	
>>RT Beam Limiting Device Type	(300A,00B8)	CS	1	Supported value: MLCX.
>>Isocenter to Beam Limiting Device Distance	(300A,00BB)	FL	2	
>>Number of Leaf/Jaw Pairs	(300A,00BC)	IS	1	
>>Leaf Position Boundaries	(300A,00BE)	DS	1C	
>Referenced Patient Setup Number	(300C,006A)	IS	3	
>Treatment Delivery Type	(300A,00CE)	CS	1	If the Treatment Delivery type is null, beam is assumed to be of type TREATMENT. Supported value: TREATMENT.
>Number of Wedges	(300A,00D0)	IS	1	Supported value: <ul style="list-style-type: none"> <li>0 - Import of wedges for ion plans is not supported.</li> </ul>
>Number of Compensators	(300A,00E0)	IS	1	Supported values: 0, 1.

>Ion Range Compensator Sequence	{300A,02EA}	SQ	1C	
>>Compensator Number	{300A,00E4}	IS	1	
>>Material ID	{300A,00E1}	SH	2	
>>Compensator ID	{300A,00E5}	SH	3	
>>Isocenter to Compensator Tray Distance	{300A,02E4}	FL	1C	
>>Compensator Divergence	{300A,02E0}	CS	1	Supported values: ABSENT, PRESENT.
>>Compensator Mounting Position	{300A,02E1}	CS	1	Supported values: PATIENT_SIDE, SOURCE_SIDE, DOUBLE_SIDED.
>>Compensator Rows	{300A,00E7}	IS	1	
>>Compensator Columns	{300A,00E8}	IS	1	
>>Compensator Pixel Spacing	{300A,00E9}	DS	1	
>>Compensator Position	{300A,00EA}	DS	1	
>>Compensator Column Offset	{300A,02E5}	FL	1C	
>>Compensator Thickness Data	{300A,00EC}	DS	1	
>>Compensator Milling Tool Diameter	{300A,02E8}	FL	3	
>Number of Boli	{300A,00ED}	IS	1	Supported value: 0.
>Number of Blocks	{300A,00F0}	IS	1	Supported values: 0, 1.
>Ion Block Sequence	{300A,03A6}	SQ	1C	
>>Isocenter to Block Tray Distance	{300A,00F7}	FL	1	
>>Block Type	{300A,00F8}	CS	1	Supported values: SHIELDING, APERTURE.
>>Block Divergence	{300A,00FA}	CS	1	Supported values: PRESENT, ABSENT.
>>Block Mounting Position	{300A,00FB}	CS	1	Supported values: PATIENT_SIDE, SOURCE_SIDE.
>>Block Number	{300A,00FC}	IS	1	
>>Block Name	{300A,00FE}	LO	3	
>>Material ID	{300A,00E1}	SH	2	
>>Block Thickness	{300A,0100}	DS	1	
>>Block Number of Points	{300A,0104}	IS	1	
>>Block Data	{300A,0106}	DS	1	
>>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>>Collimation Mode	{4001,1057}	CS	3	Private tag. Used to signal that the ion beam is collimated using Multi-Layer Aperture (MLA). Supported value: DYNAMIC.
>Snout Sequence	{300A,030C}	SQ	3	Maximum one Snout is supported.
>>Snout ID	{300A,030F}	SH	1	
>Number of Range Shifters	{300A,0312}	IS	1	Supported values: 0, 1.

>Range Shifter Sequence	{300A,0314}	SQ	1C	Supported values: <ul style="list-style-type: none"> <li>• 0</li> <li>• 1 - Scan Mode {300A,0308} must be MODULATED</li> </ul>
>>Range Shifter Number	{300A,0316}	IS	1	
>>Range Shifter ID	{300A,0318}	SH	1	Must be 'HYPER_RS' for Mevion Hyperscan.
>>Range Shifter Type	{300A,0320}	CS	1	Supported values: ANALOG, BINARY.
>Number of Lateral Spreading Devices	{300A,0330}	IS	1	0-2 Scatterers and 0 or 2 Magnets Supported values: 0, 1, 2, 3, 4.
>Lateral Spreading Device Sequence	{300A,0332}	SQ	1C	Only used for Wobbling and Mitsubishi US plans
>>Lateral Spreading Device Number	{300A,0334}	IS	1	
>>Lateral Spreading Device ID	{300A,0336}	SH	1	
>>Lateral Spreading Device Type	{300A,0338}	CS	1	Supported values: SCATTERER, MAGNET.
>Number of Range Modulators	{300A,0340}	IS	1	Supported values: 0, 1, 2.
>Range Modulator Sequence	{300A,0342}	SQ	1C	Not used for IBA and Mevion plans, since private SOBPs Width and Depth tags are used to communicate nozzle settings.
>>Range Modulator Number	{300A,0344}	IS	1	
>>Range Modulator ID	{300A,0346}	SH	1	
>>Range Modulator Type	{300A,0348}	CS	1	Supported values: FIXED, WHL_FIXEDWEIGHTS, WHL_MODWEIGHTS.
>>Range Modulator Description	{300A,034A}	LO	3	
>Patient Support Type	{300A,0350}	CS	1	Supported values: TABLE, CHAIR.
>Final Cumulative Meterset Weight	{300A,010E}	DS	1C	
>Number of Control Points	{300A,0110}	IS	1	Must be 2 for Wobbling, Double Scattering, and Uniform Scanning. Must be equal or greater than 2 for Pencil Beam Scanning.
>Ion Control Point Sequence	{300A,03A8}	SQ	1	
>>Control Point Index	{300A,0112}	IS	1	
>>Cumulative Meterset Weight	{300A,0134}	DS	2	
>>Nominal Beam Energy	{300A,0114}	DS	1C	
>>Meterset Rate	{300A,035A}	FL	3	
>>Ion Wedge Position Sequence	{300A,03AC}	SQ	1C	
>>>Referenced Wedge Number	{300C,00C0}	IS	1	
>>>Wedge Position	{300A,0118}	CS	1	Supported values: IN, OUT.
>>>Wedge Thin Edge Position	{300A,00DB}	FL	1C	
>>Range Shifter Settings Sequence	{300A,0360}	SQ	1C	Only read for Mevion Hyperscan.
>>>Referenced Range Shifter Number	{300C,0100}	IS	1	



>>>Range Shifter Setting	(300A,0362)	LO	1	For Mevion Hyperscan, the imported value must consist of a string '101011' representing energy selector plates starting with the most upstream plate (furthest away from isocenter), where each '1' and '0' indicates whether the specific plate is used or not.
>>Lateral Spreading Device Settings Sequence	(300A,0370)	SQ	1C	
>>>Referenced Lateral Spreading Device Number	(300C,0102)	IS	1	
>>>Lateral Spreading Device Setting	(300A,0372)	LO	1	
>>>Isocenter to Lateral Spreading Device Distance	(300A,0374)	FL	3	
>>>Lateral Spreading Device Water Equivalent Thickness	(300A,033C)	FL	3	
>>Range Modulator Settings Sequence	(300A,0380)	SQ	1C	
>>>Referenced Range Modulator Number	(300C,0104)	IS	1	
>>>Range Modulator Gating Start Water Equivalent Thickness	(300A,0386)	FL	3	
>>>Range Modulator Gating Stop Water Equivalent Thickness	(300A,0388)	FL	3	
>>Beam Limiting Device Position Sequence	(300A,011A)	SQ	1C	
>>>RT Beam Limiting Device Type	(300A,00B8)	CS	1	Supported value: MLCX.
>>>Leaf/Jaw Positions	(300A,011C)	DS	1	
>>Gantry Angle	(300A,011E)	DS	1C	
>>Gantry Rotation Direction	(300A,011F)	CS	1C	Supported value: NONE.
>>Gantry Pitch Angle	(300A,014A)	FL	2C	
>>Gantry Pitch Rotation Direction	(300A,014C)	CS	2C	Supported value: NONE.
>>Beam Limiting Device Angle	(300A,0120)	DS	1C	
>>Beam Limiting Device Rotation Direction	(300A,0121)	CS	1C	Supported value: NONE.
>>Scan Spot Tune ID	(300A,0390)	SH	1C	
>>Number of Scan Spot Positions	(300A,0392)	IS	1C	
>>Scan Spot Position Map	(300A,0394)	FL	1C	
>>Scan Spot Meterset Weights	(300A,0396)	FL	1C	
>>Scanning Spot Size	(300A,0398)	FL	3	
>>Number of Paintings	(300A,039A)	IS	1C	
>>Patient Support Angle	(300A,0122)	DS	1C	Must be constant within beam.
>>Patient Support Rotation Direction	(300A,0123)	CS	1C	Supported value: NONE.
>>Table Top Pitch Angle	(300A,0140)	FL	2C	

>>Table Top Pitch Rotation Direction	{300A,0142}	CS	2C	Supported value: NONE.
>>Table Top Roll Angle	{300A,0144}	FL	2C	
>>Table Top Roll Rotation Direction	{300A,0146}	CS	2C	Supported value: NONE.
>>Table Top Vertical Position	{300A,0128}	DS	2C	Supported value: 0.0.
>>Table Top Longitudinal Position	{300A,0129}	DS	2C	Supported value: 0.0.
>>Table Top Lateral Position	{300A,012A}	DS	2C	Supported value: 0.0.
>>Snout Position	{300A,030D}	FL	2C	
>>Isocenter Position	{300A,012C}	DS	2C	
>>IMPAC Private Creator	{300B,0010}	LO	3	IMPAC
>>Line Spot Tune ID	{300B,1090}	SH	3	Private tag. Used for Sumitomo Line Scanning.
>>Number of Line Scan Spot Positions	{300B,1092}	IS	3	Private tag. Used for Sumitomo Line Scanning.
>>Line Scan Position Map	{300B,1094}	FL	3	Private tag. Used for Sumitomo Line Scanning.
>>Line Scan Meterset Weights	{300B,1096}	FL	3	Private tag. Used for Sumitomo Line Scanning.
>>Line Scanning Spot Size	{300B,1098}	FL	3	Private tag. Used for Sumitomo Line Scanning.
>>Number of Line Scan Spot Paintings	{300B,109A}	IS	3	Private tag. Used for Sumitomo Line Scanning.
>>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>>Spill Length	{4001,1005}	FL	3	RaySearch Private tag. The synchrotron spill length in seconds.
>>Degradation	{4001,1006}	FL	3	RaySearch Private tag. The degradation applied by the degrader before the synchrotron.
>>Particles Per Spill	{4001,1007}	FL	3	RaySearch Private tag. Number of particles delivered during a spill for synchrotrons.
>>CNAO Nominal Beam Energy	{4001,1008}	DS	3	RaySearch Private tag. Nominal Beam Energy in MeV/nucleon. This value is copied to "Nominal Beam Energy" if not existing.
>>Reference Beam Data ID	{4001,1030}	SH	3	RaySearch Private tag. Selected Reference Beam Data ID. Only imported for Mitsubishi Uniform Scanning
>>Number Of Aperture Contour Points	{4001,1055}	IS	3	Private tag. Used for Multi-Layer Aperture (MLA). Number of contour points in the control point specific aperture contour.
>>Aperture Contour Point Map	{4001,1056}	FL	3	Private tag. Used for Multi-Layer Aperture (MLA). Contour point positions defined at isocenter.
>IMPAC Private Creator	{300B,0010}	LO	3	IMPAC
>Maximum Collimated Field Diameter	{300B,1002}	FL	3	Private tag. The maximum diameter [in mm] of a circle, centered about the beam axis, which contains the collimated field.
>Nominal SOBP Width	{300B,100E}	FL	3	Private tag. Distance [in mm] between maximal Water-Equivalent distance to distal border of target and minimal Water-Equivalent Distance to proximal border of target. Required for passive plans. Supported values: [0,320]

>Planned Distal Target Distance	{300B,1004}	FL	3	Private tag. Maximal Water-Equivalent distance (in mm) to distal border of target. Required for passive plans. Supported values: [35,320].
>IBA Private Creator	{300D,0010}	LO	3	IBA
>IBA Scattered Mode	{300D,1002}	SH	3	Private tag. Used to depict the scattering mode. Only used when Scan Mode {300A, 0308} is NONE. Supported values: SINGLE, DOUBLE.
>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>Treatment Machine Commission Time	{4001,1001}	DT	3	RaySearch Private tag. The commission time of the treatment machine. Used together with Treatment Machine Name to determine the correct machine.
>RBE Model Name	{4001,1002}	ST	3	RaySearch Private tag. The RBE Model Name. Used to determine the correct RBE Model.
>RBE Model Commission Time	{4001,1003}	DT	3	RaySearch Private tag. The commission time of the RBE Model. Used together with RBE Model Name to determine the correct RBE Model.
>Block Milling Tool Diameter	{4001,1004}	FL	3	RaySearch Private tag. The block milling tool diameter in mm.
>Internal Treatment Machine Name	{4001,1012}	SH	3	RaySearch Private tag. The internal treatment machine name. This value will differ from Treatment Machine Name {300A,00B2} if a treatment machine name alias have been specified on the ion beam quality.
>Treatment Machine RBE Mode	{4001,1033}	CS	3	Private tag. Indicates RBE mode of treatment machine. Supported values: RBE_INCLUDED, RBE_NOT_INCLUDED.

#### 7.1.8.10 Approval Module

Attribute name	Tag	Vr	Type	Comment
Approval Status	{300E,0002}	CS	1	Used to set the approval status of the imported structure set. Supported values: <ul style="list-style-type: none"> <li>APPROVED - Plan is approved.</li> <li>UNAPPROVED - Plan is not approved.</li> <li>REJECTED - Same rules as for UNAPPROVED.</li> </ul>
Review Date	{300E,0004}	DA	2C	Stored if Approval Status is APPROVED.
Review Time	{300E,0005}	TM	2C	Stored if Approval Status is APPROVED.
Reviewer Name	{300E,0008}	PN	2C	Stored if Approval Status is APPROVED.

#### 7.1.8.11 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	{0008,0016}	UI	1	Supported value: 1.2.840.10008.5.1.4.1.1.481.8.
SOP Instance UID	{0008,0018}	UI	1	Stored internally to find references from other objects.
Specific Character Set	{0008,0005}	CS	1C	Supported values: ISO_IR 100, ISO_IR 192.
Instance Creation Date	{0008,0012}	DA	3	Stored if last save date if RT Plan Date {300A, 0006} is missing.

Instance Creation Time	{0008,0013}	TM	3	Stored if last save time if RT Plan Time {300A, 0007} is missing.
------------------------	-------------	----	---	---

## 7.2 EXPORT IOD DEFINITIONS

### 7.2.1 CT Image IOD

Patient demographics is replaced in exported images to match patient. The following applies to exported 4DCT projections.

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	Yes
	Clinical Trial Study Module	Yes
Series	General Series Module	Yes
Frame of Reference	Frame of Reference Module	Yes
Equipment	General Equipment Module	Yes
Image	General Image Module	No
	Image Plane Module	Yes
	Image Pixel Module	Yes
	Contrast/Bolus Module	No
	CT Image Module	Yes
	Multi-energy CT Image Module	No
	SOP Common Module	Yes

#### 7.2.1.1 Patient Module

Attribute name	Tag	Vr	Type	Comment
Patient's Name	{0010,0010}	PN	2	
Patient ID	{0010,0020}	LO	2	
Patient's Birth Date	{0010,0030}	DA	2	
Patient's Sex	{0010,0040}	CS	2	Possible values: <ul style="list-style-type: none"> <li>• M - Male</li> <li>• F - Female</li> <li>• O - Other</li> </ul>
Other Patient Names	{0010,1001}	PN	3	Populated with the original Patient's Name if different from patient.
Other Patient IDs	{0010,1000}	LO	3	Populated with the original Patient ID if different from patient.

#### 7.2.1.2 General Study Module

Attribute name	Tag	Vr	Type	Comment
Study Instance UID	{0020,0000}	UI	1	Based on the referenced image set.

Study Date	{0008,0020}	DA	2	Date when 4DCT projection was created.
Study Time	{0008,0030}	TM	2	Time when 4DCT projection was created.
Accession Number	{0008,0050}	SH	2	
Study Description	{0008,1030}	LO	3	Always: "4DCT Projection".

### 7.2.1.3 Clinical Trial Study Module

Attribute name	Tag	Vr	Type	Comment
----------------	-----	----	------	---------

### 7.2.1.4 General Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	{0008,0060}	CS	1	Always CT.
Series Instance UID	{0020,000E}	UI	1	Generated RaySearch UID.
Series Number	{0020,0011}	IS	2	
Series Date	{0008,0021}	DA	3	
Series Time	{0008,0031}	TM	3	
Series Description	{0008,103E}	LO	3	Always: "Projection method: [projection method used]".
Patient Position	{0018,5100}	CS	2C	

### 7.2.1.5 Frame of Reference Module

Attribute name	Tag	Vr	Type	Comment
Frame of Reference UID	{0020,0052}	UI	1	Based on the image sets used for projection.

### 7.2.1.6 General Equipment Module

Attribute name	Tag	Vr	Type	Comment
Manufacturer	{0008,0070}	LO	2	Always RaySearch Laboratories.
Manufacturer's Model Name	{0008,1090}	LO	3	Always RayStation.
Software Versions	{0018,1020}	LO	3	Always the current application version and the application version when object was last saved. May also contain the filter version applied when exporting the object.

### 7.2.1.7 Image Plane Module

Attribute name	Tag	Vr	Type	Comment
Pixel Spacing	{0028,0030}	DS	1	
Image Orientation (Patient)	{0020,0037}	DS	1	
Image Position (Patient)	{0020,0032}	DS	1	Note that this is the center of the first pixel, i.e. the corner of the image offset by half a pixel.

### 7.2.1.8 Image Pixel Module

Attribute name	Tag	Vr	Type	Comment
----------------	-----	----	------	---------

Samples per Pixel	{0028,0002}	US	1	Always 1.
Photometric Interpretation	{0028,0004}	CS	1	Always MONOCHROME2.
Rows	{0028,0010}	US	1	
Columns	{0028,0011}	US	1	
Bits Allocated	{0028,0100}	US	1	Always 16.
Bits Stored	{0028,0101}	US	1	Always 16.
High Bit	{0028,0102}	US	1	Always 15.
Pixel Representation	{0028,0103}	US	1	Always 1.
Pixel Data	{7FE0,0010}	OB or OW	1C	
Pixel Padding Range Limit	{0028,0121}	US or SS	1C	

### 7.2.1.9 CT Image Module

Attribute name	Tag	Vr	Type	Comment
Image Type	{0008,0008}	CS	1	
Samples per Pixel	{0028,0002}	US	1	Always 1.
Photometric Interpretation	{0028,0004}	CS	1	Always MONOCHROME2.
Bits Allocated	{0028,0100}	US	1	Always 16.
Bits Stored	{0028,0101}	US	1	Always 16.
High Bit	{0028,0102}	US	1	Always 16.
Rescale Intercept	{0028,1052}	DS	1	Always 0.0.
Rescale Slope	{0028,1053}	DS	1	Always 1.0.
Rescale Type	{0028,1054}	LO	1C	Always HU.

### 7.2.1.10 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	{0008,0016}	UI	1	Always 1.2.840.10008.5.1.4.1.1.2.
SOP Instance UID	{0008,0018}	UI	1	Generated RaySearch UID.
Specific Character Set	{0008,0005}	CS	1C	Possible values: ISO_IR 100, ISO_IR 192.
Instance Creation Date	{0008,0012}	DA	3	
Instance Creation Time	{0008,0013}	TM	3	

### 7.2.2 RT Image IOD

RT Images (DRR) are created during export. Content and annotations may not be identical to the DRR window.

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	Yes

Series	RT Series Module	Yes
Frame of Reference	Frame of Reference Module	Yes
Equipment	General Equipment Module	Yes
Image	General Image Module	Yes
	Image Pixel Module	Yes
	Contrast/Bolus Module	No
	Cine Module	No
	Multi-frame Module	No
	RT Image Module	Yes
	VOI LUT Module	Yes
	SOP Common Module	Yes
	Frame Extraction Module	No

### 7.2.2.1 Patient Module

Attribute values based on the patient demographics.

Attribute name	Tag	Vr	Type	Comment
Patient's Name	{0010,0010}	PN	2	
Patient ID	{0010,0020}	LO	2	
Patient's Birth Date	{0010,0030}	DA	2	
Patient's Sex	{0010,0040}	CS	2	Possible values: <ul style="list-style-type: none"> <li>• M - Male</li> <li>• F - Female</li> <li>• 0 - Other</li> </ul>

### 7.2.2.2 General Study Module

Attribute name	Tag	Vr	Type	Comment
Study Instance UID	{0020,000D}	UI	1	
Study Date	{0008,0020}	DA	2	
Study Time	{0008,0030}	TM	2	
Referring Physician's Name	{0008,0090}	PN	2	
Study ID	{0020,0010}	SH	2	
Accession Number	{0008,0050}	SH	2	
Study Description	{0008,1030}	LO	3	

### 7.2.2.3 RT Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	{0008,0060}	CS	1	Always RTIMAGE.

Series Instance UID	{0020,000E}	UI	1	Generated RaySearch UID.
Series Number	{0020,0011}	IS	2	Always 1.
Series Date	{0008,0021}	DA	3	
Series Time	{0008,0031}	TM	3	

#### 7.2.2.4 Frame of Reference Module

Attribute name	Tag	Vr	Type	Comment
Frame of Reference UID	{0020,0052}	UI	1	Always same as the referenced image series, except for plans with treatment setup Ocular Gaze.

#### 7.2.2.5 General Equipment Module

Attribute name	Tag	Vr	Type	Comment
Manufacturer	{0008,0070}	LO	2	Always RaySearch Laboratories.
Institution Name	{0008,0080}	LO	3	As defined in Clinic Settings.
Manufacturer's Model Name	{0008,1090}	LO	3	Always RayStation.
Software Versions	{0018,1020}	LO	3	Always the current application version and the application version when object was last saved. May also contain the filter version applied when exporting the object.

#### 7.2.2.6 General Image Module

Attribute name	Tag	Vr	Type	Comment
Instance Number	{0020,0013}	IS	2	Exported from originating Beam Number.
Patient Orientation	{0020,0020}	CS	2C	
Content Date	{0008,0023}	DA	2C	Last save date of Beam Set.
Content Time	{0008,0033}	TM	2C	Last save time of Beam Set.
Acquisition Number	{0020,0012}	IS	3	
Burned In Annotation	{0028,0301}	CS	3	Possible values: YES, NO.

#### 7.2.2.7 Image Pixel Module

Attribute name	Tag	Vr	Type	Comment
Samples per Pixel	{0028,0002}	US	1	Always 1.
Photometric Interpretation	{0028,0004}	CS	1	Always MONOCHROME2.
Rows	{0028,0010}	US	1	
Columns	{0028,0011}	US	1	
Bits Allocated	{0028,0100}	US	1	Always 16.
Bits Stored	{0028,0101}	US	1	Always 16.
High Bit	{0028,0102}	US	1	Always 15.
Pixel Representation	{0028,0103}	US	1	



Pixel Data	{7FE0,0010}	OB or OW	1C	
Pixel Padding Range Limit	{0028,0121}	US or SS	1C	

### 7.2.2.8 RT Image Module

Attribute name	Tag	Vr	Type	Comment
Samples per Pixel	{0028,0002}	US	1	Always 1.
Photometric Interpretation	{0028,0004}	CS	1	Always MONOCHROME2.
Bits Allocated	{0028,0100}	US	1	Always 16.
Bits Stored	{0028,0101}	US	1	Always 16.
High Bit	{0028,0102}	US	1	Always 15.
Pixel Representation	{0028,0103}	US	1	Always 0.
RT Image Label	{3002,0002}	SH	1	Beam Number or Beam Name (see Clinic Settings).
RT Image Name	{3002,0003}	LO	3	Beam Number or Beam Name (see Clinic Settings).
RT Image Description	{3002,0004}	ST	3	Combines Beam Name or Beam Number and Beam Description (see Clinic Settings).
Image Type	{0008,0008}	CS	1	Possible values: DERIVED, SECONDARY, DRR.
Conversion Type	{0008,0064}	CS	2	Always WSD.
RT Image Plane	{3002,000C}	CS	1	Always NORMAL.
X-Ray Image Receptor Angle	{3002,000E}	DS	2	
RT Image Orientation	{3002,0010}	DS	2C	
Image Plane Pixel Spacing	{3002,0011}	DS	2	
RT Image Position	{3002,0012}	DS	2	
Radiation Machine Name	{3002,0020}	SH	2	Name of treatment machine or setup image device.
Primary Dosimeter Unit	{300A,00B3}	CS	2	
Radiation Machine SAD	{3002,0022}	DS	2	
RT Image SID	{3002,0026}	DS	2	
Referenced RT Plan Sequence	{300C,0002}	SQ	3	
>Referenced SOP Class UID	{0008,1150}	UI	1	
>Referenced SOP Instance UID	{0008,1155}	UI	1	
Referenced Beam Number	{300C,0006}	IS	3	
Referenced Fraction Group Number	{300C,0022}	IS	3	
Gantry Angle	{300A,011E}	DS	3	
Gantry Pitch Angle	{300A,014A}	FL	3	Only used for fixed setup beams
Beam Limiting Device Angle	{300A,0120}	DS	3	
Patient Support Angle	{300A,0122}	DS	3	
Table Top Pitch Angle	{300A,0140}	FL	3	

Table Top Roll Angle	{300A,0144}	FL	3	
Isocenter Position	{300A,012C}	DS	3	
Patient Position	{0018,5100}	CS	1C	
RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
ROI Contour Sequence	{4001,1040}	SQ	3	Contours projected unto the isocenter plane.
>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>Contour Name	{4001,1041}	LO	1	
>Contour Description	{4001,1042}	ST	2	
>Contour Generation Algorithm	{4001,1043}	CS	1	Possible values: <ul style="list-style-type: none"> <li>• AUTOMATIC - Calculated ROI</li> <li>• SEMIAUTOMATIC - Contour calculated with user assistance</li> <li>• MANUAL - User-entered ROI</li> </ul>
>Contour Interpreted Type	{4001,1044}	CS	1	Possible values: AVOIDANCE, BOLUS, CAVITY, CTV, EXTERNAL, FIXATION, GTV, IRRAD_VOLUME, ORGAN, PTV, SUPPORT, TREATED_VOLUME, NONE.
>Contour Display Color	{4001,1045}	IS	2	RGB triplet color representation for ROI, specified using the range 0-255.
>Referenced Structure Set SOP Instance UID	{4001,1050}	UI	1C	
>Referenced ROI Number	{4001,1051}	IS	1C	
>Contour Data Sequence	{4001,1052}	SQ	3	
>>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>>Pixel Contour Data	{4001,1049}	OF	2	A data stream of (x,y) pairs that comprise the contours. The number of pairs shall be equal to Number of Pixel Contour Points {4001,1048}. May contain multiple contours that does not connect. Coordinates are projected onto the IEC GANTRY coordinate system and translated to image pixel indices.
>>Number of Pixel Contour Points	{4001,1048}	IS	2	
Varian Medical Systems	{3243,0010}	LO	3	Varian Medical Systems VISION 3243
Related Referenced Beam Number	{3243,1028}	IS	3	Value copied from Referenced Beam Number.

### 7.2.2.9 VOI LUT Module

Attribute name	Tag	Vr	Type	Comment
Window Center	{0028,1050}	DS	1C	
Window Width	{0028,1051}	DS	1C	

### 7.2.2.10 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	{0008,0016}	UI	1	Always 1.2.840.10008.5.1.4.1.1.481.1.

SOP Instance UID	{0008,0018}	UI	1	Generated RaySearch UID.
Specific Character Set	{0008,0005}	CS	1C	Possible values: ISO_IR 100, ISO_IR 192.
Instance Creation Date	{0008,0012}	DA	3	Last save time.
Instance Creation Time	{0008,0013}	TM	3	

### 7.2.3 RT Dose IOD

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	Yes
Series	RT Series Module	Yes
Frame of Reference	Frame of Reference Module	Yes
Equipment	General Equipment Module	Yes
Dose	General Image Module	No
	Image Plane Module	Yes
	Image Pixel Module	Yes
	Multi-frame Module	Yes
	RT Dose Module	Yes
	Structure Set Module	No
	ROI Contour Module	No
	RT Dose ROI Module	No
	SOP Common Module	Yes
	Frame Extraction Module	No

#### 7.2.3.1 Patient Module

Attribute values based on the patient demographics.

Attribute name	Tag	Vr	Type	Comment
Patient's Name	{0010,0010}	PN	2	
Patient ID	{0010,0020}	LO	2	
Patient's Birth Date	{0010,0030}	DA	2	
Patient's Sex	{0010,0040}	CS	2	Possible values: <ul style="list-style-type: none"> <li>• M - Male</li> <li>• F - Female</li> <li>• 0 - Other</li> </ul>

#### 7.2.3.2 General Study Module

Attribute values based on the referenced image set.

Attribute name	Tag	Vr	Type	Comment
----------------	-----	----	------	---------

Study Instance UID	(0020,000D)	UI	1	
Study Date	(0008,0020)	DA	2	
Study Time	(0008,0030)	TM	2	
Referring Physician's Name	(0008,0090)	PN	2	
Study ID	(0020,0010)	SH	2	
Accession Number	(0008,0050)	SH	2	
Study Description	(0008,1030)	LO	3	

### 7.2.3.3 RT Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	(0008,0060)	CS	1	Always RTDOSE.
Series Instance UID	(0020,000E)	UI	1	Generated RaySearch UID.
Series Number	(0020,0011)	IS	2	Always 1.
Series Date	(0008,0021)	DA	3	Last save date.
Series Time	(0008,0031)	TM	3	Last save time.
Series Description	(0008,103E)	LO	3	Exported from Beam Set Name.

### 7.2.3.4 Frame of Reference Module

Attribute name	Tag	Vr	Type	Comment
Frame of Reference UID	(0020,0052)	UI	1	Always same as the referenced image series, except for plans with treatment setup Ocular Gaze.

### 7.2.3.5 General Equipment Module

Attribute name	Tag	Vr	Type	Comment
Manufacturer	(0008,0070)	LO	2	Always RaySearch Laboratories.
Institution Name	(0008,0080)	LO	3	As defined in Clinic Settings.
Manufacturer's Model Name	(0008,1090)	LO	3	Always RayStation.
Software Versions	(0018,1020)	LO	3	Always the current application version and the application version when object was last saved. May also contain the filter version applied when exporting the object.

### 7.2.3.6 Image Plane Module

Attribute name	Tag	Vr	Type	Comment
Pixel Spacing	(0028,0030)	DS	1	Exported from voxel size.
Image Orientation (Patient)	(0020,0037)	DS	1	Always '1\0\0\0\1\0' for 3D doses. Plane dependant for 2D doses.
Image Position (Patient)	(0020,0032)	DS	1	Note that this is the center of the first pixel, i.e. the corner of the image offset by half a pixel.
Slice Thickness	(0018,0050)	DS	2	Only exported for 3D doses. Exported from voxel size.

## 7.2.3.7 Image Pixel Module

Attribute name	Tag	Vr	Type	Comment
Samples per Pixel	(0028,0002)	US	1	Always 1.
Photometric Interpretation	(0028,0004)	CS	1	Always MONOCHROME2.
Rows	(0028,0010)	US	1	
Columns	(0028,0011)	US	1	
Bits Allocated	(0028,0100)	US	1	Always 16.
Bits Stored	(0028,0101)	US	1	Always 16.
High Bit	(0028,0102)	US	1	Always 15.
Pixel Representation	(0028,0103)	US	1	Always 0.
Pixel Data	{7FE0,0010}	OB or OW	1C	

## 7.2.3.8 Multi-frame Module

Not written for planar (2D) doses.

Attribute name	Tag	Vr	Type	Comment
Number of Frames	(0028,0008)	IS	1	Only used for 3D doses. Always equal the number of dose slices in z-direction.
Frame Increment Pointer	(0028,0009)	AT	1	Only used for 3D doses. Always '{3004,000c}'.

## 7.2.3.9 RT Dose Module

Attribute name	Tag	Vr	Type	Comment
Samples per Pixel	(0028,0002)	US	1C	Always 1.
Photometric Interpretation	(0028,0004)	CS	1C	Always MONOCHROME2.
Bits Allocated	(0028,0100)	US	1C	Always 16.
Bits Stored	(0028,0101)	US	1C	Always 16.
High Bit	(0028,0102)	US	1C	Always 15.
Pixel Representation	(0028,0103)	US	1C	Always 0.
Content Date	(0008,0023)	DA	3	Last save date.
Content Time	(0008,0033)	TM	3	Last save time.
Dose Units	(3004,0002)	CS	1	Always GY.
Dose Type	(3004,0004)	CS	1	Possible values: PHYSICAL, EFFECTIVE.
Instance Number	(0020,0013)	IS	3	In-series counter starting from 1.
Dose Comment	(3004,0006)	LO	3	Exported from Beam Set Description.

Dose Summation Type	{3004,000A}	CS	1	Dose Summation Type FRACTION is not exported. Summation type BEAM refers to the dose contribution from this beam over the entire treatment course, not per fraction. Possible values: <ul style="list-style-type: none"> <li>• PLAN</li> <li>• BEAM</li> <li>• EVALUATION</li> <li>• BNCT_GAMMA - BNCT Gamma Plan Dose Component</li> <li>• BNCT_BORON - BNCT Boron Plan Dose Component</li> <li>• BNCT_NITROGEN - BNCT Nitrogen Plan Dose Component</li> <li>• BNCT_HYDROGEN - BNCT Hydrogen Plan Dose Component</li> <li>• BNCT_NORMAL_RBE - BNCT Normal RBE Plan Dose Component</li> <li>• BNCT_TUMOR_RBE - BNCT Tumor RBE Plan Dose Component</li> <li>• PLAN_OCULAR_GAZE - Used for plans with treatment setup Ocular Gaze.</li> <li>• BEAM_OCULAR_GAZE - Used for plans with treatment setup Ocular Gaze.</li> </ul>
Referenced RT Plan Sequence	{300C,0002}	SQ	1C	
>Referenced SOP Class UID	{0008,1150}	UI	1	Possible values: 1.2.840.10008.5.1.4.1.1.481.5, 1.2.840.10008.5.1.4.1.1.481.8.
>Referenced SOP Instance UID	{0008,1155}	UI	1	SOP Instance UID for the referenced plan.
>Referenced Fraction Group Sequence	{300C,0020}	SQ	1C	Only written if Dose Summation Type is BEAM.
>>Referenced Fraction Group Number	{300C,0022}	IS	1	Always 1.
>>Referenced Beam Sequence	{300C,0004}	SQ	1C	
>>>Referenced Beam Number	{300C,0006}	IS	1	Exported from Beam Number.
Grid Frame Offset Vector	{3004,000C}	DS	1C	Only written for 3D doses.
Dose Grid Scaling	{3004,000E}	DS	1C	Dose Grid Scaling is always maximized for highest possible resolution.
Tissue Heterogeneity Correction	{3004,0014}	CS	3	Possible values: WATER, IMAGE, ROI_OVERRIDE.
RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
ROI Index Pixel Data	{4001,1016}	OB	3	Contains ROI material index mapping used for BNCT dose computation.
ROI Index Bits Allocated	{4001,1017}	US	3	Number of bits allocated for each pixel sample

### 7.2.3.10 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	{0008,0016}	UI	1	Always 1.2.840.10008.5.1.4.1.1.481.2.
SOP Instance UID	{0008,0018}	UI	1	Generated RaySearch UID.
Specific Character Set	{0008,0005}	CS	1C	Possible values: ISO_IR 100, ISO_IR 192.

Instance Creation Date	(0008,0012)	DA	3	Last save date.
Instance Creation Time	(0008,0013)	TM	3	Last save time.

### 7.2.4 RT Structure Set IOD

All geometries are converted into contours when exporting RT Structure Sets. An imported structure set will always be converted from internal representation and a new instance will be created.

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	Yes
Series	RT Series Module	Yes
Equipment	General Equipment Module	Yes
Frame of Reference	Frame of Reference Module	Yes
Structure Set	Structure Set Module	Yes
	ROI Contour Module	Yes
	RT ROI Observations Module	Yes
	Approval Module	Yes
	SOP Common Module	Yes

#### 7.2.4.1 Patient Module

Attribute values based on the patient demographics.

Attribute name	Tag	Vr	Type	Comment
Patient's Name	(0010,0010)	PN	2	
Patient ID	(0010,0020)	LO	2	
Patient's Birth Date	(0010,0030)	DA	2	
Patient's Sex	(0010,0040)	CS	2	Possible values: <ul style="list-style-type: none"> <li>• M - Male</li> <li>• F - Female</li> <li>• O - Other</li> </ul>

#### 7.2.4.2 General Study Module

Attribute values based on the referenced image set.

Attribute name	Tag	Vr	Type	Comment
Study Instance UID	(0020,000D)	UI	1	
Study Date	(0008,0020)	DA	2	
Study Time	(0008,0030)	TM	2	
Referring Physician's Name	(0008,0090)	PN	2	
Study ID	(0020,0010)	SH	2	

Study Description	(0008,1030)	LO	3	
-------------------	-------------	----	---	--

#### 7.2.4.3 RT Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	(0008,0060)	CS	1	Always RTSTRUCT.
Series Instance UID	(0020,000E)	UI	1	Generated RaySearch UID.
Series Number	(0020,0011)	IS	2	Always 1.
Series Date	(0008,0021)	DA	3	Last saved date.
Series Time	(0008,0031)	TM	3	Last saved time.
Series Description	(0008,103E)	LO	3	Always "RS: Approved Structure Set" or "RS: Unapproved Structure Set"

#### 7.2.4.4 General Equipment Module

Attribute name	Tag	Vr	Type	Comment
Manufacturer	(0008,0070)	LO	2	Always RaySearch Laboratories.
Institution Name	(0008,0080)	LO	3	As defined in Clinic Settings.
Manufacturer's Model Name	(0008,1090)	LO	3	Always RayStation.
Software Versions	(0018,1020)	LO	3	Always the current application version and the application version when object was last saved. May also contain the filter version applied when exporting the object.

#### 7.2.4.5 Frame of Reference Module

Attribute name	Tag	Vr	Type	Comment
Frame of Reference UID	(0020,0052)	UI	1	

#### 7.2.4.6 Structure Set Module

Attribute name	Tag	Vr	Type	Comment
Structure Set Label	(3006,0002)	SH	1	Always 'RS: Approved Structure Set' or 'RS: Unapproved Structure Set'.
Instance Number	(0020,0013)	IS	3	Always 1.
Structure Set Date	(3006,0008)	DA	2	Last save date.
Structure Set Time	(3006,0009)	TM	2	Last save time.
Referenced Frame of Reference Sequence	(3006,0010)	SQ	3	
>Frame of Reference UID	(0020,0052)	UI	1	Based on the referenced image set.
>RT Referenced Study Sequence	(3006,0012)	SQ	3	Always contains one referenced study.
>>Referenced SOP Class UID	(0008,1150)	UI	1	RT Detached Study Storage UID. Always 1.2.840.10008.3.1.2.3.1.
>>Referenced SOP Instance UID	(0008,1155)	UI	1	Study UID from the image series for which the structure set is created.



>>RT Referenced Series Sequence	{3006,0014}	SQ	1	Always contains one referenced series.
>>>Series Instance UID	{0020,000E}	UI	1	Series UID from the image series for which the structure set is created.
>>>Contour Image Sequence	{3006,0016}	SQ	1	All images from the image series for which the structure set is created is added to the sequence. Images for which there is no contour defined are also added.
>>>>Referenced SOP Class UID	{0008,1150}	UI	1	
>>>>Referenced SOP Instance UID	{0008,1155}	UI	1	
Structure Set ROI Sequence	{3006,0020}	SQ	1	All ROI:s and POI:s for which there exists a contour/point on the referenced image series are added to the sequence.
>ROI Number	{3006,0022}	IS	1	
>Referenced Frame of Reference UID	{3006,0024}	UI	1	Always same as the image series frame of reference for which the structure set is created.
>ROI Name	{3006,0026}	LO	2	Exported from ROI Name.
>ROI Generation Algorithm	{3006,0036}	CS	2	Always SEMIAUTOMATIC.

#### 7.2.4.7 ROI Contour Module

Attribute name	Tag	Vr	Type	Comment
ROI Contour Sequence	{3006,0039}	SQ	1	
>Referenced ROI Number	{3006,0084}	IS	1	Exported from ROI Number.
>ROI Display Color	{3006,002A}	IS	3	Exported from ROI/POI Display Color.
>Contour Sequence	{3006,0040}	SQ	3	
>>Contour Number	{3006,0048}	IS	3	
>>Contour Image Sequence	{3006,0016}	SQ	3	May contain one or no contour image. If there is an image in the referenced image series that that is positioned on the same z-coordinate as the contour, that image will be referenced in this sequence. If no image is found, the sequence will be empty.
>>>Referenced SOP Class UID	{0008,1150}	UI	1	
>>>Referenced SOP Instance UID	{0008,1155}	UI	1	
>>Contour Geometric Type	{3006,0042}	CS	1	Possible values: POINT, CLOSED_PLANAR.
>>Number of Contour Points	{3006,0046}	IS	1	
>>Contour Data	{3006,0050}	DS	1	

#### 7.2.4.8 RT ROI Observations Module

Attribute name	Tag	Vr	Type	Comment
RT ROI Observations Sequence	{3006,0080}	SQ	1	
>Observation Number	{3006,0082}	IS	1	
>Referenced ROI Number	{3006,0084}	IS	1	

>ROI Observation Label	{3006,0085}	SH	3	Same as the name of the referenced ROI, except if the ROI is of type Localization Poi, in which case the label is set to "Localization Poi".
>RT ROI Interpreted Type	{3006,00A4}	CS	2	All DICOM standard terms are supported. Localization Poi exported as "MARKER".
>ROI Interpreter	{3006,00A6}	PN	2	
>Material ID	{300A,00E1}	SH	3	Exported from material name.
>ROI Physical Properties Sequence	{3006,00B0}	SQ	3	The values REL_MASS_DENSITY, REL_ELEC_DENSITY, EFFECTIVE_Z, EFF_Z_PER_A and ELEM_FRACTION are supported. However at least one of REL_MASS_DENSITY or REL_ELEC_DENSITY is needed to create a density override. If only one of these values exist, the physical properties of water is used and mass density is set primarily from REL_MASS_DENSITY and secondarily from REL_ELEC_DENSITY.
>>ROI Physical Property	{3006,00B2}	CS	1	MEAN_EXCI_ENERGY depicts the Mean Excitation Energy for a given material. Possible values: <ul style="list-style-type: none"> <li>• REL_MASS_DENSITY</li> <li>• REL_ELEC_DENSITY</li> <li>• EFFECTIVE_Z</li> <li>• EFF_Z_PER_A</li> <li>• ELEM_FRACTION</li> <li>• MEAN_EXCI_ENERGY</li> <li>• CBE_BORON - Used for BNCT dose computation</li> <li>• RBE_HYDROGEN - Used for BNCT dose computation</li> <li>• RBE_NITROGEN - Used for BNCT dose computation</li> <li>• RBE_GAMMA - Used for BNCT dose computation</li> <li>• BORON_CONC_RATIO - Boron Concentration Ratio. Used for BNCT dose computation</li> <li>• BNCT_DOSE_LIMIT - Used for BNCT dose computation</li> </ul>
>>ROI Elemental Composition Sequence	{3006,00B6}	SQ	1C	
>>>ROI Elemental Composition Atomic Number	{3006,00B7}	US	1	
>>>ROI Elemental Composition Atomic Mass Fraction	{3006,00B8}	FL	1	
>>ROI Physical Property Value	{3006,00B4}	DS	1	
>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>Tissue Name	{4001,1010}	ST	3	RaySearch Private Tag. Contains the tissue name given to this Region of Interest.
>RBE Cell Type Name	{4001,1014}	ST	3	RaySearch Private Tag. Contains the RBE Cell Type Name given to this Region of Interest.

### 7.2.4.9 Approval Module

Attribute name	Tag	Vr	Type	Comment
----------------	-----	----	------	---------

Approval Status	{300E,0002}	CS	1	Possible values: <ul style="list-style-type: none"> <li>APPROVED - Structure set is approved.</li> <li>UNAPPROVED - Structure set is not approved.</li> </ul>
Review Date	{300E,0004}	DA	2C	Approval date.
Review Time	{300E,0005}	TM	2C	Approval time.
Reviewer Name	{300E,0008}	PN	2C	User that approved data.

#### 7.2.4.10 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	{0008,0016}	UI	1	Always 1.2.840.10008.5.1.4.1.1.481.3.
SOP Instance UID	{0008,0018}	UI	1	Generated RaySearch UID.
Specific Character Set	{0008,0005}	CS	1C	Possible values: ISO_IR 100, ISO_IR 192.
Instance Creation Date	{0008,0012}	DA	3	Last save time.
Instance Creation Time	{0008,0013}	TM	3	Last save time.

#### 7.2.5 RT Plan IOD

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	Yes
Series	RT Series Module	Yes
Frame of Reference	Frame of Reference Module	Yes
Equipment	General Equipment Module	Yes
Plan	RT General Plan Module	Yes
	RT Prescription Module	Yes
	RT Patient Setup Module	Yes
	RT Fraction Scheme Module	Yes
	RT Beams Module	Yes
	RT Brachy Application Setups Module	No
	Approval Module	Yes
	SOP Common Module	Yes

##### 7.2.5.1 Patient Module

Attribute values based on the patient demographics.

Attribute name	Tag	Vr	Type	Comment
Patient's Name	{0010,0010}	PN	2	
Patient ID	{0010,0020}	LO	2	
Patient's Birth Date	{0010,0030}	DA	2	

Patient's Sex	(0010,0040)	CS	2	Possible values: <ul style="list-style-type: none"> <li>• M - Male</li> <li>• F - Female</li> <li>• O - Other</li> </ul>
---------------	-------------	----	---	--

### 7.2.5.2 General Study Module

Attribute values based on the referenced image set.

Attribute name	Tag	Vr	Type	Comment
Study Instance UID	(0020,000D)	UI	1	
Study Date	(0008,0020)	DA	2	
Study Time	(0008,0030)	TM	2	
Referring Physician's Name	(0008,0090)	PN	2	
Study ID	(0020,0010)	SH	2	
Accession Number	(0008,0050)	SH	2	
Study Description	(0008,1030)	LO	3	

### 7.2.5.3 RT Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	(0008,0060)	CS	1	Always RTPLAN.
Series Instance UID	(0020,000E)	UI	1	Generated RaySearch UID.
Series Number	(0020,0011)	IS	2	
Series Date	(0008,0021)	DA	3	
Series Time	(0008,0031)	TM	3	
Series Description	(0008,103E)	LO	3	Exported from Beam Set Name.
Operators' Name	(0008,1070)	PN	2	

### 7.2.5.4 Frame of Reference Module

Attribute name	Tag	Vr	Type	Comment
Frame of Reference UID	(0020,0052)	UI	1	Always same as the referenced image series.
Position Reference Indicator	(0020,1040)	LO	2	

### 7.2.5.5 General Equipment Module

Attribute name	Tag	Vr	Type	Comment
Manufacturer	(0008,0070)	LO	2	Always RaySearch Laboratories.
Institution Name	(0008,0080)	LO	3	As defined in Clinic Settings.
Manufacturer's Model Name	(0008,1090)	LO	3	Always RayStation.

Software Versions	{0018,1020}	LO	3	Always the current application version and the application version when object was last saved. May also contain the filter version applied when exporting the object.
-------------------	-------------	----	---	---

### 7.2.5.6 RT General Plan Module

Attribute name	Tag	Vr	Type	Comment
RT Plan Label	{300A,0002}	SH	1	Exported from Beam Set Name.
RT Plan Name	{300A,0003}	LO	3	Exported from Treatment Plan Name.
RT Plan Description	{300A,0004}	ST	3	Exported from Beam Set Description.
RT Plan Date	{300A,0006}	DA	2	Last save date.
RT Plan Time	{300A,0007}	TM	2	Last save time.
Treatment Protocols	{300A,0009}	LO	3	Exported as plan delivery technique. Possible values: SMLC, DMLC, TomoHelical, TomoStatic, DynamicArc, StaticArc, TomoDirect, BNCT, NCT.
Plan Intent	{300A,000A}	CS	3	Treatment plans is always CURATIVE, QA plans always VERIFICATION, Beam Commission plans always RESEARCH, and plans with recalculated dose might have RECALCULATION. Possible values: CURATIVE, VERIFICATION, MACHINE_QA, RESEARCH, RECALCULATION, BNCT_DOSE_CALC.
Treatment Sites	{300A,000B}	LO	3	
RT Plan Geometry	{300A,000C}	CS	1	Always PATIENT.
Referenced Structure Set Sequence	{300C,0060}	SQ	1C	Based on the referenced Examination Structure Set.
>Referenced SOP Class UID	{0008,1150}	UI	1	Always 1.2.840.10008.5.1.4.1.1.481.3.
>Referenced SOP Instance UID	{0008,1155}	UI	1	
Referenced RT Plan Sequence	{300C,0002}	SQ	3	
>Referenced SOP Class UID	{0008,1150}	UI	1	Always 1.2.840.10008.5.1.4.1.1.481.5.
>Referenced SOP Instance UID	{0008,1155}	UI	1	
>RT Plan Relationship	{300A,0055}	CS	1	Possible values: VERIFIED_PLAN, EQUIVALENT, CONCURRENT.
>Brainlab Private Creator	{300B,0012}	LO	3	Brainlab - ONC - Beam Parameters
>Referenced Beam List	{300B,1210}	IS	3	Only used for Brainlab plans.

### 7.2.5.7 RT Prescription Module

Attribute name	Tag	Vr	Type	Comment
Prescription Description	{300A,000E}	ST	3	Exported from RT Plan Label.
Dose Reference Sequence	{300A,0010}	SQ	3	Primary Prescription supported only.
>Dose Reference Number	{300A,0012}	IS	1	Always 1.
>Dose Reference UID	{300A,0013}	UI	3	
>Dose Reference Structure Type	{300A,0014}	CS	1	Possible values: POINT, VOLUME, SITE.

>Dose Reference Description	(300A,0016)	LO	3	Exported from Site Description or textualized prescription.
>Referenced ROI Number	(3006,0084)	IS	1C	
>Dose Reference Type	(300A,0020)	CS	1	Always TARGET.
>Constraint Weight	(300A,0021)	DS	3	
>Target Minimum Dose	(300A,0025)	DS	3	
>Target Prescription Dose	(300A,0026)	DS	3	
>Target Underdose Volume Fraction	(300A,0028)	DS	3	Written for prescriptions of type Dose At Volume, Median Dose (always 50) and Near Minimum Dose (always 2).
>Organ at Risk Maximum Dose	(300A,002C)	DS	3	
>Organ at Risk Overdose Volume Fraction	(300A,002D)	DS	3	

### 7.2.5.8 RT Patient Setup Module

Attribute name	Tag	Vr	Type	Comment
Patient Setup Sequence	(300A,0180)	SQ	1	
>Patient Setup Number	(300A,0182)	IS	1	
>Patient Position	(0018,5100)	CS	1C	Possible values: HFS, HFP, FFS, FFP, HFDR, HFDL, FFDR, FFDL.
>Referenced Setup Image Sequence	(300A,0401)	SQ	3	
>>Referenced SOP Class UID	(0008,1150)	UI	1	
>>Referenced SOP Instance UID	(0008,1155)	UI	1	
>Setup Technique	(300A,01B0)	CS	3	Possible values: ISOCENTRIC, FIXED_SSD.
>Setup Technique Description	(300A,01B2)	ST	3	
>Table Top Vertical Setup Displacement	(300A,01D2)	DS	3	Vertical displacement in IEC TABLE TOP coordinate system (in mm) relative to Localization Poi. Value is only written if the Localization Poi is defined on the planning CT. The displacement is defined as (Localization Poi - isocenter). Export of table top displacement can be disabled in the RayPhysics per machine.
>Table Top Longitudinal Setup Displacement	(300A,01D4)	DS	3	Longitudinal displacement in IEC TABLE TOP coordinate system (in mm) relative to Localization Poi. Value is only written if the Localization Poi is defined on the planning CT. The displacement is defined as (Localization Poi - isocenter). Export of table top displacement can be disabled in the RayPhysics per machine.
>Table Top Lateral Setup Displacement	(300A,01D6)	DS	3	Lateral displacement in IEC TABLE TOP coordinate system (in mm) relative to Localization Poi. Value is only written if the Localization Poi is defined on the planning CT. The displacement is defined as (Localization Poi - isocenter). Export of table top displacement can be disabled in the RayPhysics per machine.

>Motion Synchronization Sequence	{300A,0410}	SQ	3	
>>Respiratory Motion Compensation Technique	{0018,9170}	CS	1	Possible values: <ul style="list-style-type: none"> <li>• NONE</li> <li>• GATING</li> <li>• TRACKING - Static target tracking</li> <li>• REALTIME - Realtime tracking</li> </ul>
>>Respiratory Signal Source	{0018,9171}	CS	1	Possible values: <ul style="list-style-type: none"> <li>• NONE</li> <li>• EXTERNAL_MARKER</li> <li>• SPIROMETER</li> <li>• BELT</li> <li>• INTERNAL_MARKER</li> <li>• FIDUCIAL - Tracking of fiducial markers</li> <li>• TRACKING_ROI - Tracking of volume</li> </ul>
>>Respiratory Motion Compensation Technique Description	{0018,9185}	ST	3	
>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>Tomo Localization point	{4001,1022}	DS	3	
>Tracking ROI Sequence	{4001,1062}	SQ	3	Contains referenced ROI used for tracking
>>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>>Referenced ROI Number	{4001,1051}	IS	1C	
>Imaging Angles	{4001,1063}	DS	3	Selected imaging angles for Tomo Synchrony
>Tomo Couch Insertion Position Y	{4001,1066}	DS	3	
>medPhoton Private Creator	{30BB,0010}	LO	3	medPhoton 1.0
>Patient Support Angle Offset Interval	{30BB,1022}	DS	3	The offset interval for the yaw angle in degrees to use when performing collision checks. Relates to Patient Support Angle {300A, 0122}. The first value should be zero or negative and the second value should be zero or positive. Example: [0.0, 0.0], [-2.0, 2.0] or [-2.0, 3.0].
>Table Top Pitch Angle Offset Interval	{30BB,1040}	DS	3	The offset interval for the pitch angle in degrees to use when performing collision checks. Relates to Table Top Pitch Angle {300A, 0140}. The first value should be zero or negative and the second value should be zero or positive. Example: [0.0, 0.0], [-2.0, 2.0] or [-2.0, 3.0].
>Table Top Roll Angle Offset Interval	{30BB,1044}	DS	3	The offset interval for the roll angle in degrees to use when performing collision checks. Relates to Table Top Roll Angle {300A, 0144}. The first value defines the negative offset and should be zero or negative. The second value defines the positive offset and should be zero or positive. Example: [0.0, 0.0], [-2.0, 2.0] or [-2.0, 3.0].

>Isocenter Position Offset Interval	(30BB,102C)	DS	3	The offset interval for the isocenter in millimeter to use when performing collision checks. Relates to Isocenter Position (300A, 012C). The first three values define the negative offsets and should be zero or negative. The last three values define the positive offsets and should be zero or positive. Example: [0.0, 0.0, 0.0, 0.0, 0.0, 0.0], [-2.0, -3.0, 0.0, 2.0, 1.0, 0.0].
medPhoton Private Creator	(30BB,0010)	LO	3	medPhoton 1.0
Patient Setup ID	(30BB,1000)	SH	3	Exported from Patient Setup IDCAS parameters.
Imaging Protocol ID	(30BB,1001)	SH	3	Exported from Patient Setup IDCAS parameters.

### 7.2.5.9 RT Fraction Scheme Module

Attribute name	Tag	Vr	Type	Comment
Fraction Group Sequence	(300A,0070)	SQ	1	This sequence will always contain one item. Multiple fraction groups are not supported.
>Fraction Group Number	(300A,0071)	IS	1	Always 1.
>Referenced Dose Reference Sequence	(300C,0050)	SQ	3	
>>Referenced Dose Reference Number	(300C,0051)	IS	1	Always 1.
>Number of Fractions Planned	(300A,0078)	IS	2	
>Number of Fraction Pattern Digits Per Day	(300A,0079)	IS	3	
>Repeat Fraction Cycle Length	(300A,007A)	IS	3	
>Fraction Pattern	(300A,007B)	LT	3	
>Number of Beams	(300A,0080)	IS	1	
>Referenced Beam Sequence	(300C,0004)	SQ	1C	
>>Referenced Beam Number	(300C,0006)	IS	1	
>>Referenced Dose Reference UID	(300A,0083)	UI	3	
>>Beam Dose	(300A,0084)	DS	3	Exported from Beam Dose at specified point, unless the user has chosen to "Export beam dose specification point dose as fraction of prescription dose" for the treatment machine.
>>Beam Meterset	(300A,0086)	DS	3	Total MU for beam.
>>Beam Dose Point Depth	(300A,0088)	FL	3	Should be interpreted as the distance from the surface of the patient to the beam dose specification point, i.e. including support, fixture and boli.
>>Beam Dose Point Equivalent Depth	(300A,0089)	FL	3	For photons: The water equivalent depth of the beam dose specification point, including support, fixation and boli. For electrons: The integrated depth from the surface to the beam dose specification point, including support, fixation and boli, weighted with electron density relative to water.
>>Beam Dose Point SSD	(300A,008A)	FL	3	Should be interpreted as the distance to the surface of the patient, including support, fixation and boli.
>>Beam Dose Specification Point	(300A,0082)	DS	3	Exported from Dose Specification Point.



>>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>>Beam Dose Specification Point Name	{4001,1029}	ST	3	The name of the Beam Dose Specification Point.
>Number of Brachy Application Setups	{300A,00A0}	IS	1	Not supported. Always 0.

### 7.2.5.10 RT Beams Module

Attribute name	Tag	Vr	Type	Comment
Beam Sequence	{300A,00B0}	SQ	1	
>Beam Number	{300A,00C0}	IS	1	Exported from Beam Name.
>Beam Name	{300A,00C2}	LO	3	Exported from Beam Name.
>Beam Description	{300A,00C3}	ST	3	Exported from Beam Description.
>Beam Type	{300A,00C4}	CS	1	Always STATIC for SMLC plans. Always DYNAMIC for Arc and DMLC plans. Possible values: STATIC, DYNAMIC.
>Radiation Type	{300A,00C6}	CS	2	Possible values: PHOTON, ELECTRON, NEUTRON.
>Primary Fluence Mode Sequence	{3002,0050}	SQ	3	
>>Fluence Mode	{3002,0051}	CS	1	If a non standard fluence mode is used for the beam quality this attribute will be set to NON_STANDARD. Otherwise this attribute is STANDARD Possible values: STANDARD, NON_STANDARD.
>>Fluence Mode ID	{3002,0052}	SH	1C	Identifier for the specific fluence-shaping mode. If "Fluence Mode" is set to NON_STANDARD this attribute must be set to either FFF or SRS. If "Fluence Mode" is set to STANDARD this attribute is null. Possible values: FFF, SRS.
>Treatment Machine Name	{300A,00B2}	SH	2	
>Primary Dosimeter Unit	{300A,00B3}	CS	3	Possible values: MU, MINUTE.
>Source-Axis Distance	{300A,00B4}	DS	3	
>Beam Limiting Device Sequence	{300A,00B6}	SQ	1	
>>RT Beam Limiting Device Type	{300A,00B8}	CS	1	Possible values: X, Y, ASYMX, ASYMY, MLCX, MLCY.
>>Number of Leaf/Jaw Pairs	{300A,00BC}	IS	1	
>>Leaf Position Boundaries	{300A,00BE}	DS	2C	
>Referenced Patient Setup Number	{300C,006A}	IS	3	
>Treatment Delivery Type	{300A,00CE}	CS	3	Possible values: TREATMENT, SETUP.
>Number of Wedges	{300A,00D0}	IS	1	Possible values: 0, 1.
>Wedge Sequence	{300A,00D1}	SQ	1C	
>>Wedge Number	{300A,00D2}	IS	1	Always 1.
>>Wedge Type	{300A,00D3}	CS	2	Possible values: STANDARD, DYNAMIC, MOTORIZED.
>>Wedge ID	{300A,00D4}	SH	3	
>>Wedge Angle	{300A,00D5}	IS	2	Always 60 for motorized wedge.

>>Wedge Factor	{300A,00D6}	DS	2	
>>Wedge Orientation	{300A,00D8}	DS	2	Possible values: 0.0, 90.0, 180.0, 270.0.
>>Effective Wedge Angle	{300A,00DE}	DS	3	
>Number of Compensators	{300A,00E0}	IS	1	Always 0.
>Number of Boli	{300A,00ED}	IS	1	
>Referenced Bolus Sequence	{300C,00B0}	SQ	1C	
>>Referenced ROI Number	{3006,0084}	IS	1	
>>Bolus ID	{300A,00DC}	SH	3	
>Number of Blocks	{300A,00F0}	IS	1	Electron plans support 0 to 1 blocks. Photon plans support 0 to 2 blocks. Possible values: 0, 1, 2.
>Total Block Tray Factor	{300A,00F2}	DS	3	
>Block Sequence	{300A,00F4}	SQ	1C	
>>Block Tray ID	{300A,00F5}	SH	3	
>>Source to Block Tray Distance	{300A,00F6}	DS	2	
>>Block Type	{300A,00F8}	CS	1	Possible values: SHIELDING, APERTURE.
>>Block Divergence	{300A,00FA}	CS	2	Possible values: PRESENT, ABSENT.
>>Block Mounting Position	{300A,00FB}	CS	3	Possible values: PATIENT_SIDE, SOURCE_SIDE.
>>Block Number	{300A,00FC}	IS	1	
>>Block Name	{300A,00FE}	LO	3	
>>Block Transmission	{300A,0102}	DS	2C	
>>Block Number of Points	{300A,0104}	IS	2	
>>Block Data	{300A,0106}	DS	2	
>Applicator Sequence	{300A,0107}	SQ	3	Used for Electron applicators or Photon cones.
>>Applicator ID	{300A,0108}	SH	1	
>>Applicator Type	{300A,0109}	CS	1	Possible values: ELECTRON_SQUARE, ELECTRON_RECT, ELECTRON_CIRC, ELECTRON_SHORT, ELECTRON_OPEN, PHOTON_CIRC, STEREOTACTIC, NEUTRON_CIRC.
>>Applicator Geometry Sequence	{300A,0431}	SQ	3	Only exported for photon cones.
>>>Applicator Aperture Shape	{300A,0432}	CS	1	Always SYM_CIRCULAR.
>>>Applicator Opening	{300A,0433}	FL	1C	
>>Applicator Description	{300A,010A}	LO	3	
>>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>>Applicator Thickness	{4001,1054}	DS	3	Physical thickness of applicator (in mm) parallel to radiation beam axis.
>Final Cumulative Meterset Weight	{300A,010E}	DS	1C	Always 1.0.
>Number of Control Points	{300A,0110}	IS	1	
>Control Point Sequence	{300A,0111}	SQ	1	

>>Control Point Index	(300A,0112)	IS	1	
>>Cumulative Meterset Weight	(300A,0134)	DS	2	Meterset weight at beginning of control point. Final control point will always have value 1.
>>Referenced Dose Reference Sequence	(300C,0050)	SQ	3	
>>>Referenced Dose Reference Number	(300C,0051)	IS	1	Always 1.
>>>Cumulative Dose Reference Coefficient	(300A,010C)	DS	2	Be aware these values are linearly interpolated from 0 to the final value and are not accurate.
>>Nominal Beam Energy	(300A,0114)	DS	3	Written as the nominal energy of the machine beam quality assigned to the beam. The energy used internally for dose calculation may differ from the nominal energy and will not be exported.
>>Dose Rate Set	(300A,0115)	DS	3	
>>Wedge Position Sequence	(300A,0116)	SQ	1C	If plan contains wedges.
>>>Referenced Wedge Number	(300C,00C0)	IS	1	Always 1.
>>>Wedge Position	(300A,0118)	CS	1	Possible values: IN, OUT.
>>Beam Limiting Device Position Sequence	(300A,011A)	SQ	1C	
>>>RT Beam Limiting Device Type	(300A,00B8)	CS	1	Possible values: X, Y, ASYMX, ASYMY, MLCX, MLCY.
>>>Leaf/Jaw Positions	(300A,011C)	DS	1	
>>Gantry Angle	(300A,011E)	DS	1C	
>>Gantry Rotation Direction	(300A,011F)	CS	1C	Possible values: NONE, CW, CC.
>>Gantry Pitch Rotation Direction	(300A,014C)	CS	3	Possible values: NONE, CW, CC.
>>Beam Limiting Device Angle	(300A,0120)	DS	1C	Same throughout entire beam.
>>Beam Limiting Device Rotation Direction	(300A,0121)	CS	1C	Possible values: NONE, CW, CC.
>>Patient Support Angle	(300A,0122)	DS	1C	Same throughout entire beam except for Wave Arc plans.
>>Patient Support Rotation Direction	(300A,0123)	CS	1C	Possible values: NONE, CW, CC.
>>Table Top Eccentric Axis Distance	(300A,0124)	DS	3	
>>Table Top Eccentric Angle	(300A,0125)	DS	1C	Always 0.0.
>>Table Top Eccentric Rotation Direction	(300A,0126)	CS	1C	Possible values: NONE, CW, CC.
>>Table Top Pitch Angle	(300A,0140)	FL	1C	
>>Table Top Pitch Rotation Direction	(300A,0142)	CS	1C	Only written in first control point for SMLC, DMLC and Static Arc. Possible values: NONE, CW, CC.
>>Table Top Roll Angle	(300A,0144)	FL	1C	
>>Table Top Roll Rotation Direction	(300A,0146)	CS	1C	Only written in first control point for SMLC and DMLC. Possible values: NONE, CW, CC.
>>Table Top Vertical Position	(300A,0128)	DS	2C	Always 0.0.

>>Table Top Longitudinal Position	{300A,0129}	DS	2C	Always 0.0.
>>Table Top Lateral Position	{300A,012A}	DS	2C	
>>Isocenter Position	{300A,012C}	DS	2C	Always written in first control point.
>>Source to Surface Distance	{300A,0130}	DS	3	Should be interpreted as the distance to the surface of the patient, i.e. including support, fixation and boli.
>>TomoTherapy HA Private Creator	{300D,0010}	LO	3	TOMO_HA_01
>>Tomo Projection Sinogram Data	{300D,10A7}	DS	3	TomoHelical plans only. Projection sinogram values for all 64 leaves of projection starting at the current control point. Empty if all 64 leaves have a sinogram value of zero. Empty for last control point.
>Brainlab Private Creator	{320B,0010}	LO	3	Brainlab - ONC - Multi-axial treatment machine
>Dynamic Tracking	{320B,1001}	CS	3	Dynamic Tracking for Vero. Possible values: ENABLED, DISABLED.
>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>Treatment Machine Commission Time	{4001,1001}	DT	3	The commission time of the treatment machine. Used together with Treatment Machine Name to determine the correct machine.
>Tomo IDMS Beam ID	{4001,1021}	LO	3	
>Tomo IDMS Machine ID	{4001,1020}	LO	3	
>Tomo Beam Isocenter	{4001,1023}	DS	3	
>Tomo machine revision	{4001,1025}	IS	3	
>Tomo beam revision	{4001,1026}	IS	3	
>Tomo intended back jaw position	{4001,1027}	DS	3	
>Tomo intended front jaw position	{4001,1028}	DS	3	
>Planned Blood Boron Concentration	{4001,1015}	DS	3	Populated for BNCT plans only.
>Tomo Projection Time FL	{4001,1068}	FL	3	Projection time for Tomo Helical and Tomo Direct plans.
>Varian Medical Systems	{3243,0010}	LO	3	Varian Medical Systems VISION 3243
>Beam Secondary Name	{3243,1009}	LO	3	Value copied from Beam Description.

### 7.2.5.11 Approval Module

Attribute name	Tag	Vr	Type	Comment
Approval Status	{300E,0002}	CS	1	Possible values: <ul style="list-style-type: none"> <li>APPROVED - Plan is approved.</li> <li>UNAPPROVED - Plan is not approved.</li> </ul>
Review Date	{300E,0004}	DA	2C	Approval date.
Review Time	{300E,0005}	TM	2C	Approval time.
Reviewer Name	{300E,0008}	PN	2C	User that approved data.

### 7.2.5.12 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	{0008,0016}	UI	1	Always 1.2.840.10008.5.1.4.1.1.481.5.
SOP Instance UID	{0008,0018}	UI	1	Generated RaySearch UID.
Specific Character Set	{0008,0005}	CS	1C	Possible values: ISO_IR 100, ISO_IR 192.
Instance Creation Date	{0008,0012}	DA	3	Last save date.
Instance Creation Time	{0008,0013}	TM	3	Last save time.
RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
Tomo Plan SOP Instance UID	{4001,102A}	UI	3	The SOP Instance UID that shall be used by RayGateway for creating plan instances in IDMS based on this RT Plan.

## 7.2.6 Spatial Registration IOD

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	Yes
Series	General Series Module	Yes
	Spatial Registration Series Module	Yes
Frame of Reference	Frame of Reference Module	Yes
Equipment	General Equipment Module	Yes
Spatial Registration	Spatial Registration Module	Yes
	Common Instance Reference Module	Yes
	SOP Common Module	Yes

### 7.2.6.1 Patient Module

Attribute values based on the patient demographics.

Attribute name	Tag	Vr	Type	Comment
Patient's Name	{0010,0010}	PN	2	
Patient ID	{0010,0020}	LO	2	
Patient's Birth Date	{0010,0030}	DA	2	
Patient's Sex	{0010,0040}	CS	2	Possible values: <ul style="list-style-type: none"> <li>• M - Male</li> <li>• F - Female</li> <li>• O - Other</li> </ul>

### 7.2.6.2 General Study Module

Attribute values based on the referenced image set.

Attribute name	Tag	Vr	Type	Comment
Study Instance UID	{0020,000D}	UI	1	

Study Date	(0008,0020)	DA	2	
Study Time	(0008,0030)	TM	2	
Referring Physician's Name	(0008,0090)	PN	2	
Study ID	(0020,0010)	SH	2	
Accession Number	(0008,0050)	SH	2	
Study Description	(0008,1030)	LO	3	

### 7.2.6.3 General Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	(0008,0060)	CS	1	Always REG.
Series Instance UID	(0020,000E)	UI	1	
Series Number	(0020,0011)	IS	2	
Series Date	(0008,0021)	DA	3	
Series Time	(0008,0031)	TM	3	
Protocol Name	(0018,1030)	LO	3	
Series Description	(0008,103E)	LO	3	
Patient Position	(0018,5100)	CS	2C	

### 7.2.6.4 Spatial Registration Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	(0008,0060)	CS	1	Always REG.

### 7.2.6.5 Frame of Reference Module

Attribute name	Tag	Vr	Type	Comment
Frame of Reference UID	(0020,0052)	UI	1	Always same as the referenced image series frame of reference.

### 7.2.6.6 General Equipment Module

Attribute name	Tag	Vr	Type	Comment
Manufacturer	(0008,0070)	LO	2	Always RaySearch Laboratories.
Institution Name	(0008,0080)	LO	3	As defined in Clinic Settings.
Manufacturer's Model Name	(0008,1090)	LO	3	Always RayStation.
Software Versions	(0018,1020)	LO	3	Always the current application version and the application version when object was last saved. May also contain the filter version applied when exporting the object.

### 7.2.6.7 Spatial Registration Module

Attribute name	Tag	Vr	Type	Comment
Content Date	(0008,0023)	DA	1	

Content Time	(0008,0033)	TM	1	
Instance Number	(0020,0013)	IS	1	
Content Label	(0070,0080)	CS	1	
Content Description	(0070,0081)	LO	2	
Content Creator's Name	(0070,0084)	PN	2	
Registration Sequence	(0070,0308)	SQ	1	
>Frame of Reference UID	(0020,0052)	UI	1C	
>Referenced Image Sequence	(0008,1140)	SQ	1C	
>>Referenced SOP Class UID	(0008,1150)	UI	1	
>>Referenced SOP Instance UID	(0008,1155)	UI	1	
>Matrix Registration Sequence	(0070,0309)	SQ	1	
>>Registration Type Code Sequence	(0070,030D)	SQ	2	
>>>Code Value	(0008,0100)	SH	1C	Possible values: 125021, 125023, 125024.
>>>Coding Scheme Designator	(0008,0102)	SH	1C	Always DCM.
>>>Code Meaning	(0008,0104)	LO	1	Possible values: Frame of Reference Identity, Acquisition Equipment Alignment, Image Content-based Alignment.
>>Matrix Sequence	(0070,030A)	SQ	1	
>>>Frame of Reference Transformation Matrix	(3006,00C6)	DS	1	
>>>Frame of Reference Transformation Matrix Type	(0070,030C)	CS	1	Always RIGID.

#### 7.2.6.8 Common Instance Reference Module

Attribute name	Tag	Vr	Type	Comment
Referenced Series Sequence	(0008,1115)	SQ	1C	
>Series Instance UID	(0020,000E)	UI	1	
>Referenced Instance Sequence	(0008,114A)	SQ	1	
>>Referenced SOP Class UID	(0008,1150)	UI	1	
>>Referenced SOP Instance UID	(0008,1155)	UI	1	
Studies Containing Other Referenced Instances Sequence	(0008,1200)	SQ	1C	If referenced images are located in a Study different from the REG object.
>Study Instance UID	(0020,000D)	UI	1	
>Referenced Series Sequence	(0008,1115)	SQ	1	
>>Series Instance UID	(0020,000E)	UI	1	
>>Referenced Instance Sequence	(0008,114A)	SQ	1	
>>>Referenced SOP Class UID	(0008,1150)	UI	1	
>>>Referenced SOP Instance UID	(0008,1155)	UI	1	

### 7.2.6.9 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	{0008,0016}	UI	1	Always 1.2.840.10008.5.1.4.1.1.66.1.
SOP Instance UID	{0008,0018}	UI	1	Generated RaySearch UID.
Specific Character Set	{0008,0005}	CS	1C	Possible values: ISO_IR 100, ISO_IR 192.
Instance Creation Date	{0008,0012}	DA	3	Last save date.
Instance Creation Time	{0008,0013}	TM	3	Last save time.

### 7.2.7 RT Ion Plan IOD

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	Yes
Series	RT Series Module	Yes
Frame of Reference	Frame of Reference Module	Yes
Equipment	General Equipment Module	Yes
Plan	RT General Plan Module	Yes
	RT Prescription Module	Yes
	RT Patient Setup Module	Yes
	RT Fraction Scheme Module	Yes
	RT Ion Beams Module	Yes
	Approval Module	Yes
	SOP Common Module	Yes

#### 7.2.7.1 Patient Module

Attribute values based on the patient demographics.

Attribute name	Tag	Vr	Type	Comment
Patient's Name	{0010,0010}	PN	2	
Patient ID	{0010,0020}	LO	2	
Patient's Birth Date	{0010,0030}	DA	2	
Patient's Sex	{0010,0040}	CS	2	Possible values: <ul style="list-style-type: none"> <li>• M - Male</li> <li>• F - Female</li> <li>• 0 - Other</li> </ul>

#### 7.2.7.2 General Study Module

Attribute values based on the referenced image set.

Attribute name	Tag	Vr	Type	Comment
----------------	-----	----	------	---------



Study Instance UID	(0020,000D)	UI	1	
Study Date	(0008,0020)	DA	2	
Study Time	(0008,0030)	TM	2	
Referring Physician's Name	(0008,0090)	PN	2	
Study ID	(0020,0010)	SH	2	
Accession Number	(0008,0050)	SH	2	
Study Description	(0008,1030)	LO	3	

#### 7.2.7.3 RT Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	(0008,0060)	CS	1	Always RTPLAN.
Series Instance UID	(0020,000E)	UI	1	Generated RaySearch UID.
Series Number	(0020,0011)	IS	2	
Series Date	(0008,0021)	DA	3	
Series Time	(0008,0031)	TM	3	
Series Description	(0008,103E)	LO	3	Exported from Beam Set Name.

#### 7.2.7.4 Frame of Reference Module

Attribute name	Tag	Vr	Type	Comment
Frame of Reference UID	(0020,0052)	UI	1	Always same as the referenced image series, except for plans with treatment setup Ocular Gaze.

#### 7.2.7.5 General Equipment Module

Attribute name	Tag	Vr	Type	Comment
Manufacturer	(0008,0070)	LO	2	Always RaySearch Laboratories.
Institution Name	(0008,0080)	LO	3	As defined in Clinic Settings.
Manufacturer's Model Name	(0008,1090)	LO	3	Always RayStation.
Software Versions	(0018,1020)	LO	3	Always the current application version and the application version when object was last saved. May also contain the filter version applied when exporting the object.

#### 7.2.7.6 RT General Plan Module

Attribute name	Tag	Vr	Type	Comment
RT Plan Label	(300A,0002)	SH	1	Exported from Beam Set Name.
RT Plan Name	(300A,0003)	LO	3	Exported from Treatment Plan Name.
RT Plan Description	(300A,0004)	ST	3	Exported from Beam Set Description.
RT Plan Date	(300A,0006)	DA	2	Last save time.
RT Plan Time	(300A,0007)	TM	2	Last save time.

Treatment Protocols	(300A,0009)	LO	3	Exported from plan delivery technique. Possible values: DoubleScattering, SingleScattering, PencilBeamScanning, LineScanning, UniformScanning, Wobbling, GridUniformScanning.
Plan Intent	(300A,000A)	CS	3	Treatment plans is always CURATIVE, QA plans always VERIFICATION, 3D Beam Commission plans always RESEARCH and plans with recalculated dose might have RECALCULATION. Possible values: CURATIVE, VERIFICATION, MACHINE_QA, RESEARCH, RECALCULATION.
Treatment Sites	(300A,000B)	LO	3	
RT Plan Geometry	(300A,000C)	CS	1	Always PATIENT.
Referenced Structure Set Sequence	(300C,0060)	SQ	1C	Based on the referenced Examination Structure Set.
>Referenced SOP Class UID	(0008,1150)	UI	1	Always 1.2.840.10008.5.1.4.1.1.481.3.
>Referenced SOP Instance UID	(0008,1155)	UI	1	
Referenced Dose Sequence	(300C,0080)	SQ	3	
>Referenced SOP Class UID	(0008,1150)	UI	1	
>Referenced SOP Instance UID	(0008,1155)	UI	1	
Referenced RT Plan Sequence	(300C,0002)	SQ	3	
>Referenced SOP Class UID	(0008,1150)	UI	1	
>Referenced SOP Instance UID	(0008,1155)	UI	1	
>RT Plan Relationship	(300A,0055)	CS	1	Possible values: VERIFIED_PLAN, EQUIVALENT.
>Brainlab Private Creator	(300B,0012)	LO	3	Brainlab - ONC - Beam Parameters
>Referenced Beam List	(300B,1210)	IS	3	Only used for Brainlab plans.

### 7.2.7.7 RT Prescription Module

Attribute name	Tag	Vr	Type	Comment
Prescription Description	(300A,000E)	ST	3	Exported from RT Plan Label.
Dose Reference Sequence	(300A,0010)	SQ	3	
>Dose Reference Number	(300A,0012)	IS	1	Always 1.
>Dose Reference UID	(300A,0013)	UI	3	
>Dose Reference Structure Type	(300A,0014)	CS	1	Possible values: POINT, VOLUME, SITE.
>Dose Reference Description	(300A,0016)	LO	3	Exported from Site Description or textualized prescription.
>Referenced ROI Number	(3006,0084)	IS	1C	
>Dose Reference Type	(300A,0020)	CS	1	Always TARGET.
>Constraint Weight	(300A,0021)	DS	3	
>Target Minimum Dose	(300A,0025)	DS	3	
>Target Prescription Dose	(300A,0026)	DS	3	The exported value is the RBE-weighted prescription dose.

>Target Underdose Volume Fraction	{300A,0028}	DS	3	Written for prescriptions of type Dose At Volume, Median Dose (always 50) and Near Minimum Dose (always 2).
>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>Target Prescription Effective Dose	{4001,1011}	DS	3	RaySearch Private Tag. Prescribed dose to Dose Reference if Dose Reference Type {300A,0020} is TARGET. The dose has been corrected for biological effect using user-defined modeling technique.

#### 7.2.7.8 RT Patient Setup Module

Attribute name	Tag	Vr	Type	Comment
Patient Setup Sequence	{300A,0180}	SQ	1	
>Patient Setup Number	{300A,0182}	IS	1	
>Patient Position	{0018,5100}	CS	1C	Not exported for plans with treatment setup Ocular Gaze. Possible values: HFS, HFP, FFS, FFP, HFDR, HFDL, FFDR, FFDL, SITTING.
>Patient Additional Position	{300A,0184}	LO	1C	Only exported for plans with treatment setup Ocular Gaze. Always OCULAR_GAZE.
>Referenced Setup Image Sequence	{300A,0401}	SQ	3	
>>Referenced SOP Class UID	{0008,1150}	UI	1	
>>Referenced SOP Instance UID	{0008,1155}	UI	1	
>Table Top Vertical Setup Displacement	{300A,01D2}	DS	3	Vertical displacement in IEC TABLE TOP coordinate system (in mm) relative to Localization Poi. Value is only written if the Localization Poi is defined on the planning CT. The displacement is defined as (Localization Poi - isocenter). Export of table top displacement can be disabled in the RayPhysics per machine.
>Table Top Longitudinal Setup Displacement	{300A,01D4}	DS	3	Longitudinal displacement in IEC TABLE TOP coordinate system (in mm) relative to Localization Poi. Value is only written if the Localization Poi is defined on the planning CT. The displacement is defined as (Localization Poi - isocenter). Export of table top displacement can be disabled in the RayPhysics per machine.
>Table Top Lateral Setup Displacement	{300A,01D6}	DS	3	Lateral displacement in IEC TABLE TOP coordinate system (in mm) relative to Localization Poi. Value is only written if the Localization Poi is defined on the planning CT. The displacement is defined as (Localization Poi - isocenter). Export of table top displacement can be disabled in the RayPhysics per machine.
>Motion Synchronization Sequence	{300A,0410}	SQ	3	
>>Respiratory Motion Compensation Technique	{0018,9170}	CS	1	Possible values: NONE, GATING.

>>Respiratory Signal Source	{0018,9171}	CS	1	Possible values: <ul style="list-style-type: none"> <li>• NONE</li> <li>• EXTERNAL_MARKER - Signal determined from external motion surrogate</li> <li>• SPIROMETER - Signal derived from flow sensor</li> <li>• BELT</li> <li>• INTERNAL_MARKER - Signal determined from internal motion surrogate</li> </ul>
>>Respiratory Interval Time	{0020,9254}	FD	3	Measured interval time in ms from maximum respiration peak to the next peak for the respiratory cycle in which this frame occurs.
>>IMPAC Private Creator	{300B,0010}	LO	3	IMPAC
>>Respiratory Phase Gating Duty Cycle	{300B,1020}	DS	3	Start and end phases for gating window for respiratory motion, given as percentages of the respiratory cycle.
medPhoton Private Creator	{30BB,0010}	LO	3	medPhoton 1.0
Patient Setup ID	{30BB,1000}	SH	3	Exported from Patient Setup IDCAS parameters.
Imaging Protocol ID	{30BB,1001}	SH	3	Exported from Patient Setup IDCAS parameters.

#### 7.2.7.9 RT Fraction Scheme Module

Attribute name	Tag	Vr	Type	Comment
Fraction Group Sequence	{300A,0070}	SQ	1	This sequence will always contain one item. Multiple fraction groups are not supported.
>Fraction Group Number	{300A,0071}	IS	1	Always 1.
>Referenced Dose Reference Sequence	{300C,0050}	SQ	3	
>>Referenced Dose Reference Number	{300C,0051}	IS	1	Always 1.
>Number of Fractions Planned	{300A,0078}	IS	2	
>Number of Fraction Pattern Digits Per Day	{300A,0079}	IS	3	Only written in plans with Mitsubishi Electric machines.
>Repeat Fraction Cycle Length	{300A,007A}	IS	3	Only written in plans with Mitsubishi Electric machines.
>Fraction Pattern	{300A,007B}	LT	3	Only written in plans with Mitsubishi Electric machines.
>Number of Beams	{300A,0080}	IS	1	
>Referenced Beam Sequence	{300C,0004}	SQ	1C	
>>Referenced Beam Number	{300C,0006}	IS	1	
>>Beam Dose	{300A,0084}	DS	3	Exported from Beam Dose at specified point, unless the user has chosen to "Export beam dose specification point dose as fraction of prescription dose" for the treatment machine.
>>Beam Dose Type	{300A,0090}	CS	1C	Possible values: PHYSICAL, EFFECTIVE.
>>Alternate Beam Dose	{300A,0091}	DS	3	

>>Alternate Beam Dose Type	{300A,0092}	CS	1C	Always EFFECTIVE.
>>Beam Meterset	{300A,0086}	DS	3	Only exported for non-passive plans.
>>Beam Dose Point Depth	{300A,0088}	FL	3	Should be interpreted as the distance from the surface of the patient to the beam dose specification point, i.e. including support, fixture and boli.
>>Beam Dose Point SSD	{300A,008A}	FL	3	Should be interpreted as the distance to the surface of the patient, including support, fixation and boli.
>>Beam Dose Specification Point	{300A,0082}	DS	3	Exported from Dose Specification Point.
>Number of Brachy Application Setups	{300A,00A0}	IS	1	Not supported. Always 0.

#### 7.2.7.10 RT Ion Beams Module

Attribute name	Tag	Vr	Type	Comment
Ion Beam Sequence	{300A,03A2}	SQ	1	
>Beam Number	{300A,00C0}	IS	1	Exported from Beam Number.
>Beam Name	{300A,00C2}	LO	1	Exported from Beam Name.
>Beam Description	{300A,00C3}	ST	3	Exported from Beam Description.
>Beam Type	{300A,00C4}	CS	1	Always STATIC.
>Radiation Type	{300A,00C6}	CS	1	Possible values: PROTON, ION.
>Radiation Mass Number	{300A,0302}	IS	1C	Exported if Radiation Type {300A,00C6} is ION. Possible values: 4, 12.
>Radiation Atomic Number	{300A,0304}	IS	1C	Exported if Radiation Type {300A,00C6} is ION. Possible values: 2, 6.
>Radiation Charge State	{300A,0306}	SS	1C	Exported if Radiation Type {300A,00C6} is ION. Possible values: 2, 6.
>Scan Mode	{300A,0308}	CS	1	If value is NONE, the private tag IBA Scattered Mode {300D, 1002} depicts if the beam delivery type is Single Scattering or Double Scattering. Possible values: NONE, UNIFORM, MODULATED, MODULATED_SPEC, LINE, WOBBLING.
>Modulated Scan Mode Type	{300A,0309}	CS	1C	Only written if Scan Mode {300A,0308} is MODULATED_SPEC. Possible values: STATIONARY, LEAPING, LINEAR, MIXED.
>Treatment Machine Name	{300A,00B2}	SH	2	Exported from Machine Reference.
>Manufacturer	{0008,0070}	LO	3	Exported from treatment machine design.
>Manufacturer's Model Name	{0008,1090}	LO	3	Exported from treatment machine.
>Primary Dosimeter Unit	{300A,00B3}	CS	1	Possible values: MU, NP.
>Virtual Source-Axis Distances	{300A,030A}	FL	1	Equal values in X and Y if Beam Limiting Device Angle is different than zero
>Ion Beam Limiting Device Sequence	{300A,03A4}	SQ	3	
>>RT Beam Limiting Device Type	{300A,00B8}	CS	1	Always MLCX.
>>Isocenter to Beam Limiting Device Distance	{300A,00BB}	FL	2	

>>Number of Leaf/Jaw Pairs	(300A,00BC)	IS	1	
>>Leaf Position Boundaries	(300A,00BE)	DS	1C	
>Referenced Patient Setup Number	(300C,006A)	IS	3	
>Referenced Reference Image Sequence	(300C,0042)	SQ	3	
>>Referenced SOP Class UID	(0008,1150)	UI	1	
>>Referenced SOP Instance UID	(0008,1155)	UI	1	
>>Reference Image Number	(300A,00C8)	IS	1	
>Treatment Delivery Type	(300A,00CE)	CS	1	Possible values: TREATMENT, SETUP.
>Referenced Dose Sequence	(300C,0080)	SQ	3	
>>Referenced SOP Class UID	(0008,1150)	UI	1	
>>Referenced SOP Instance UID	(0008,1155)	UI	1	
>Number of Wedges	(300A,00D0)	IS	1	
>Ion Wedge Sequence	(300A,03AA)	SQ	1C	
>>Wedge Number	(300A,00D2)	IS	1	
>>Wedge Type	(300A,00D3)	CS	2	Possible values: STANDARD, MOTORIZED, PARTIAL_STANDARD, PARTIAL_MOTORIZED.
>>Wedge ID	(300A,00D4)	SH	3	
>>Accessory Code	(300A,00F9)	LO	3	
>>Wedge Angle	(300A,00D5)	IS	2	
>>Wedge Orientation	(300A,00D8)	DS	2	
>>Isocenter to Wedge Tray Distance	(300A,00D9)	FL	1	
>Number of Compensators	(300A,00E0)	IS	1	Possible values: 0, 1.
>Ion Range Compensator Sequence	(300A,02EA)	SQ	1C	
>>Compensator Number	(300A,00E4)	IS	1	Always 0.
>>Material ID	(300A,00E1)	SH	2	
>>Compensator ID	(300A,00E5)	SH	3	
>>Isocenter to Compensator Tray Distance	(300A,02E4)	FL	1C	
>>Compensator Divergence	(300A,02E0)	CS	1	Possible values: ABSENT, PRESENT.
>>Compensator Mounting Position	(300A,02E1)	CS	1	Possible values: PATIENT_SIDE, SOURCE_SIDE.
>>Compensator Rows	(300A,00E7)	IS	1	
>>Compensator Columns	(300A,00E8)	IS	1	
>>Compensator Pixel Spacing	(300A,00E9)	DS	1	
>>Compensator Position	(300A,00EA)	DS	1	
>>Compensator Column Offset	(300A,02E5)	FL	1C	
>>Compensator Thickness Data	(300A,00EC)	DS	1	

>>Compensator Milling Tool Diameter	{300A,02E8}	FL	3	
>Number of Boli	{300A,00ED}	IS	1	Always 0.
>Number of Blocks	{300A,00F0}	IS	1	Possible values: 0, 1.
>Ion Block Sequence	{300A,03A6}	SQ	1C	
>>Isocenter to Block Tray Distance	{300A,00F7}	FL	1	
>>Block Type	{300A,00F8}	CS	1	Possible values: SHIELDING, APERTURE.
>>Block Divergence	{300A,00FA}	CS	1	Possible values: PRESENT, ABSENT.
>>Block Mounting Position	{300A,00FB}	CS	1	Possible values: PATIENT_SIDE, SOURCE_SIDE.
>>Block Number	{300A,00FC}	IS	1	Always 0.
>>Block Name	{300A,00FE}	LO	3	
>>Material ID	{300A,00E1}	SH	2	
>>Block Thickness	{300A,0100}	DS	1	
>>Block Number of Points	{300A,0104}	IS	1	
>>Block Data	{300A,0106}	DS	1	
>>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>>Collimation Mode	{4001,1057}	CS	3	Private tag. Only exported for plans with Multi-Layer Aperture (MLA). Always DYNAMIC.
>Snout Sequence	{300A,030C}	SQ	3	
>>Snout ID	{300A,030F}	SH	1	
>>Accessory Code	{300A,00F9}	LO	3	Only exported for Mevion Hyperscan. Exported value is the checksum meant for use in the Mevion SMC (Spot Map Converter).
>Applicator Sequence	{300A,0107}	SQ	3	Only for plans with Mitsubishi Electric machine.
>>Applicator ID	{300A,0108}	SH	1	
>>Applicator Type	{300A,0109}	CS	1	Possible values: ION_CIRC, ION_RECT.
>Number of Range Shifters	{300A,0312}	IS	1	Possible values: 0, 1.
>Range Shifter Sequence	{300A,0314}	SQ	1C	Possible values: 0, 1.
>>Range Shifter Number	{300A,0316}	IS	1	Always 0.
>>Range Shifter ID	{300A,0318}	SH	1	Always 'HYPER_RS' for Mevion Hyperscan.
>>Accessory Code	{300A,00F9}	LO	3	
>>Range Shifter Type	{300A,0320}	CS	1	Possible values: ANALOG, BINARY.
>Number of Lateral Spreading Devices	{300A,0330}	IS	1	0-2 Scatterers and 0 or 2 Magnets Possible values: 0, 1, 2, 3, 4.
>Lateral Spreading Device Sequence	{300A,0332}	SQ	1C	Scatterer(s) are exported for Mevion Double Scattering plans, Wobbling, Mitsubishi US, and for Pencil Beam Scanning if selected in referenced machine. Magnets are exported if selected in referenced machine.
>>Lateral Spreading Device Number	{300A,0334}	IS	1	

>>Lateral Spreading Device ID	{300A,0336}	SH	1	
>>Lateral Spreading Device Type	{300A,0338}	CS	1	Possible values: SCATTERER, MAGNET.
>Number of Range Modulators	{300A,0340}	IS	1	Possible values: 0, 1, 2.
>Range Modulator Sequence	{300A,0342}	SQ	1C	Not used for IBA and Mevion plans, since private SOBP Width and Depth tags are used to communicate nozzle settings.
>>Range Modulator Number	{300A,0344}	IS	1	
>>Range Modulator ID	{300A,0346}	SH	1	
>>Range Modulator Type	{300A,0348}	CS	1	Possible values: FIXED, WHL_FIXEDWEIGHTS, WHL_MODWEIGHTS.
>>Range Modulator Description	{300A,034A}	LO	3	
>Patient Support Type	{300A,0350}	CS	1	Possible values: TABLE, CHAIR.
>Patient Support ID	{300A,0352}	SH	3	Exported from treatment machine setting. Defaults to 'TABLE'.
>Fixation Light Azimuthal Angle	{300A,0356}	FL	3	
>Fixation Light Polar Angle	{300A,0358}	FL	3	
>Fixation Eye	{300A,0150}	CS	3	Possible values: <ul style="list-style-type: none"> <li>• L - Left eye</li> <li>• R - Right eye</li> </ul>
>Final Cumulative Meterset Weight	{300A,010E}	DS	1C	
>Number of Control Points	{300A,0110}	IS	1	
>Ion Control Point Sequence	{300A,03A8}	SQ	1	
>>Control Point Index	{300A,0112}	IS	1	
>>Cumulative Meterset Weight	{300A,0134}	DS	2	
>>Referenced Dose Reference Sequence	{300C,0050}	SQ	3	
>>>Referenced Dose Reference Number	{300C,0051}	IS	1	Always 1.
>>>Cumulative Dose Reference Coefficient	{300A,010C}	DS	2	
>>Nominal Beam Energy	{300A,0114}	DS	1C	Always 0 for IBA and Mevion plans, since private SOBP Width and Depth tags are used to communicate nozzle settings.
>>Meterset Rate	{300A,035A}	FL	3	
>>Ion Wedge Position Sequence	{300A,03AC}	SQ	1C	
>>>Referenced Wedge Number	{300C,00C0}	IS	1	
>>>Wedge Position	{300A,0118}	CS	1	Possible values: IN, OUT.
>>>Wedge Thin Edge Position	{300A,00DB}	FL	1C	
>>Range Shifter Settings Sequence	{300A,0360}	SQ	1C	



>>>Referenced Range Shifter Number	(300C,0100)	IS	1	Always 0.
>>>Range Shifter Setting	(300A,0362)	LO	1	For Mevion Hyperscan, the exported value will consist of a string '101011' representing energy selector plates starting with the most upstream plate (furthest away from isocenter), where each '1' and '0' indicates whether the specific plate is used or not. The exported value is always 'IN' for other PBS delivery techniques.
>>>Isocenter to Range Shifter Distance	(300A,0364)	FL	3	
>>>Range Shifter Water Equivalent Thickness	(300A,0366)	FL	3	
>>Lateral Spreading Device Settings Sequence	(300A,0370)	SQ	1C	Exported for Mevion Double Scattering plans, Sumitomo plans, Mitsubishi Uniform Scanning and Pencil Beam Scanning if scatterer(s)/magnets are set to be exported
>>>Referenced Lateral Spreading Device Number	(300C,0102)	IS	1	
>>>Lateral Spreading Device Setting	(300A,0372)	LO	1	For a referenced scatterer Sumitomo Wobbling machine specifies that Lateral Spreading Device Setting is set to a binary string indicating used Lollipops. Otherwise value IN is used.
>>>Isocenter to Lateral Spreading Device Distance	(300A,0374)	FL	3	
>>>Lateral Spreading Device Water Equivalent Thickness	(300A,033C)	FL	3	
>>Range Modulator Settings Sequence	(300A,0380)	SQ	1C	
>>>Referenced Range Modulator Number	(300C,0104)	IS	1	
>>>Range Modulator Gating Start Water Equivalent Thickness	(300A,0386)	FL	3	
>>>Range Modulator Gating Stop Water Equivalent Thickness	(300A,0388)	FL	3	
>>Beam Limiting Device Position Sequence	(300A,011A)	SQ	1C	
>>>RT Beam Limiting Device Type	(300A,00B8)	CS	1	Always MLCX.
>>>Leaf/Jaw Positions	(300A,011C)	DS	1	
>>Gantry Angle	(300A,011E)	DS	1C	
>>Gantry Rotation Direction	(300A,011F)	CS	1C	Always NONE.
>>Gantry Pitch Angle	(300A,014A)	FL	2C	
>>Gantry Pitch Rotation Direction	(300A,014C)	CS	2C	Always NONE.
>>Beam Limiting Device Angle	(300A,0120)	DS	1C	
>>Beam Limiting Device Rotation Direction	(300A,0121)	CS	1C	Always NONE.
>>Scan Spot Tune ID	(300A,0390)	SH	1C	

>>Number of Scan Spot Positions	(300A,0392)	IS	1C	
>>Scan Spot Position Map	(300A,0394)	FL	1C	
>>Scan Spot Meterset Weights	(300A,0396)	FL	1C	
>>Scanning Spot Size	(300A,0398)	FL	3	
>>Number of Paintings	(300A,039A)	IS	1C	
>>Patient Support Angle	(300A,0122)	DS	1C	
>>Patient Support Rotation Direction	(300A,0123)	CS	1C	Always NONE.
>>Table Top Pitch Angle	(300A,0140)	FL	2C	
>>Table Top Pitch Rotation Direction	(300A,0142)	CS	2C	Always NONE.
>>Table Top Roll Angle	(300A,0144)	FL	2C	
>>Table Top Roll Rotation Direction	(300A,0146)	CS	2C	Always NONE.
>>Table Top Vertical Position	(300A,0128)	DS	2C	Always 0.0.
>>Table Top Longitudinal Position	(300A,0129)	DS	2C	Always 0.0.
>>Table Top Lateral Position	(300A,012A)	DS	2C	Always 0.0.
>>Snout Position	(300A,030D)	FL	2C	
>>Isocenter Position	(300A,012C)	DS	2C	
>>Nominal Beam Energy Unit	(300A,0015)	CS	3	Only exported for proton plans. Note: This attribute is not part of the DICOM standard. Nominal beam energy is always exported in MEV/nucleon in accordance with the DICOM standard. Possible values: MEV, KV.
>>IMPAC Private Creator	(300B,0010)	LO	3	IMPAC
>>Line Spot Tune ID	(300B,1090)	SH	3	Private tag. Used for Sumitomo Line Scanning.
>>Number of Line Scan Spot Positions	(300B,1092)	IS	3	Private tag. Used for Sumitomo Line Scanning.
>>Line Scan Position Map	(300B,1094)	FL	3	Private tag. Used for Sumitomo Line Scanning.
>>Line Scan Meterset Weights	(300B,1096)	FL	3	Private tag. Used for Sumitomo Line Scanning.
>>Line Scanning Spot Size	(300B,1098)	FL	3	Private tag. Used for Sumitomo Line Scanning.
>>Number of Line Scan Spot Paintings	(300B,109A)	IS	3	Private tag. Used for Sumitomo Line Scanning.
>>RaySearch Private Creator	(4001,0010)	LO	3	RAYSEARCHLABS 2.0
>>Spill Length	(4001,1005)	FL	3	RaySearch Private tag. The synchrotron spill length in seconds.
>>Degradation	(4001,1006)	FL	3	RaySearch Private tag. The degradation applied by the degrader before the synchrotron.
>>Particles Per Spill	(4001,1007)	FL	3	RaySearch Private tag. Number of particles delivered during a spill for synchrotrons.
>>CNAO Nominal Beam Energy	(4001,1008)	DS	3	RaySearch Private tag. Nominal Beam Energy in MeV/nucleon. This attribute is only written for CNAO machines

>>CNAO Nominal Beam Energy Scale Factor	{4001,1009}	DS	3	RaySearch Private tag. Scale Factor that should be applied to control point's Nominal Beam Energy to convert it to MeV/u.. This attribute is only written for CNAO machines
>>Reference Beam Data ID	{4001,1030}	SH	3	RaySearch Private tag. Selected Reference Beam Data ID. Only exported for Mitsubishi Uniform Scanning
>>Reference Depth in Water	{4001,1031}	DS	3	RaySearch Private tag. Calculated as $ro\_range\_max - ro\_max\_sobp\_width/2 - (range\_shifter\_wet - ro\_range\_shifter\_wet\_min)$ [mm]. Here, ro is the Range Option used for the current beam (referenced in Reference Beam Data ID {4001,1030}), ro_range_max is the maximum Range for the ro, ro_max_sobp_width is the Max SOB Width for the ro, ro_range_shifter_wet_min is the minimum Range shifter WET for the ro, and range_shifter_wet is the Range shifter WET used for the current beam. Only exported for beams using a machine of type Mitsubishi Electric Co and treatment technique Uniform Scanning
>>Mitsubishi Reference Value	{4001,1032}	DS	3	RaySearch Private tag. Computed as $proton\_fluence * IDD\_SOBP(d0)$ , where proton_fluence is the total number of protons per area unit simulated for the current Beam for one treatment fraction [protons / cm <sup>2</sup> ], IDD_SOBP(d0) is the integrated depth dose of the entire SOB Width at d0 as used by the pencil beam dose engine in RayStation [Gy * cm <sup>2</sup> / proton], and d0 is the Reference Depth in Water {4001,1031} [mm]. Only exported for beams using a machine of type Mitsubishi Electric Co and treatment technique Uniform Scanning
>>Number Of Aperture Contour Points	{4001,1055}	IS	3	Private tag. Used for Multi-Layer Aperture (MLA). Number of contour points in the control point specific aperture contour.
>>Aperture Contour Point Map	{4001,1056}	FL	3	Private tag. Used for Multi-Layer Aperture (MLA). Contour point positions defined at isocenter.
>IMPAC Private Creator	{300B,0010}	LO	3	IMPAC
>Maximum Collimated Field Diameter	{300B,1002}	FL	3	Private tag. The maximum diameter (in mm) of a circle, centered about the beam axis, which contains the collimated field.
>Nominal SOB Width	{300B,100E}	FL	3	Private tag. Distance (in mm) between maximal Water-Equivalent distance to distal border of target and minimal Water-Equivalent Distance to proximal border of target. Required for passive plans. Supported values: [0,320]
>Planned Distal Target Distance	{300B,1004}	FL	3	Private tag. Maximal Water-Equivalent distance (in mm) to distal border of target. Required for passive plans. Supported values: [35,320].
>IBA Private Creator	{300D,0010}	LO	3	IBA
>IBA Scattered Mode	{300D,1002}	SH	3	Private tag. Used to depict the scattering mode. Only used when Scan Mode {300A, 0308} is NONE. Possible values: SINGLE, DOUBLE.
>RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
>Treatment Machine Commission Time	{4001,1001}	DT	3	RaySearch Private tag. The commission time of the treatment machine. Used together with Treatment Machine Name to determine the correct machine.

>RBE Model Name	{4001,1002}	ST	3	RaySearch Private tag. The RBE Model Name. Used to determine the correct RBE Model.
>RBE Model Commission Time	{4001,1003}	DT	3	RaySearch Private tag. The commission time of the RBE Model. Used together with RBE Model Name to determine the correct RBE Model.
>Block Milling Tool Diameter	{4001,1004}	FL	3	RaySearch Private tag. The block milling tool diameter in mm.
>Internal Treatment Machine Name	{4001,1012}	SH	3	RaySearch Private tag. The internal treatment machine name. This value will differ from Treatment Machine Name {300A,00B2} if a treatment machine name alias have been specified on the ion beam quality.
>Treatment Machine RBE Mode	{4001,1033}	CS	3	Private tag. Indicates RBE mode of treatment machine. Possible values: RBE_INCLUDED, RBE_NOT_INCLUDED.
>Mevion SMC Output	{4001,1064}	OB	3	Private tag. Binary output data from Mevion SMC.
>Mevion SMC Version	{4001,1065}	LO	3	Private tag. Version of Mevion SMC used to produce "Mevion SMC Output".

### 7.2.7.11 Approval Module

Attribute name	Tag	Vr	Type	Comment
Approval Status	{300E,0002}	CS	1	Possible values: <ul style="list-style-type: none"> <li>APPROVED - Plan is approved.</li> <li>UNAPPROVED - Plan is not approved.</li> </ul>
Review Date	{300E,0004}	DA	2C	Approval date.
Review Time	{300E,0005}	TM	2C	Approval time.
Reviewer Name	{300E,0008}	PN	2C	User that approved data.

### 7.2.7.12 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	{0008,0016}	UI	1	Always 1.2.840.10008.5.1.4.1.1.481.8.
SOP Instance UID	{0008,0018}	UI	1	Generated RaySearch UID.
Specific Character Set	{0008,0005}	CS	1C	Possible values: ISO_IR 100, ISO_IR 192.
Instance Creation Date	{0008,0012}	DA	3	Last save date.
Instance Creation Time	{0008,0013}	TM	3	Last save time.

## 7.3 DATA DICTIONARY OF PRIVATE ATTRIBUTES

All used Private Creators and Private Attributes are listed in the table below. Usage of Private Attributes are listed in each module specification.

Attribute name	Tag	VR	VM	Value
Brainlab Private Creator	{300B,0012}	LO	1	Brainlab - ONC - Beam Parameters
Referenced Beam List		IS	1-n	

Brainlab Private Creator	{320B,0010}	LO	1	Brainlab - ONC - Multi-axial treatment machine
Dynamic Tracking		CS	1	
GE Private Creator	{0009,0010}	LO	1	GEMS_PETD_01
GE PET Scan Time		DT	1	
IBA Private Creator	{300D,0010}	LO	1	IBA
IBA Scattered Mode		SH	1	
IMPAC Private Creator	{300B,0010}	LO	1	IMPAC
Line Spot Tune ID		SH	1	
Number of Line Scan Spot Positions		IS	1	
Line Scan Position Map		FL	1-n	
Line Scan Meterset Weights		FL	1-n	
Line Scanning Spot Size		FL	2	
Number of Line Scan Spot Paintings		IS	1	
Maximum Collimated Field Diameter		FL	1	
Nominal SOBPs Width		FL	1	
Planned Distal Target Distance		FL	1	
Respiratory Phase Gating Duty Cycle		DS	2	
medPhoton Private Creator	{30BB,0010}	LO	1	medPhoton 1.0
Patient Setup ID		SH	1	
Imaging Protocol ID		SH	1	
Patient Support Angle Offset Interval		DS	2	
Table Top Pitch Angle Offset Interval		DS	2	
Table Top Roll Angle Offset Interval		DS	2	
Isocenter Position Offset Interval		DS	6	
Philips Private Creator	{7053,0010}	LO	1	Philips PET Private Group
Philips SUV Scale Factor		DS	1	
Philips Activity Concentration Scale Factor		DS	1	
RaySearch Private Creator	{3001,0011}	LO	1	RAYSEARCHLABS 1.0
RaySearch Private Creator	{4001,0010}	LO	1	RAYSEARCHLABS 2.0
Tissue Name		ST	1	
Target Prescription Effective Dose		DS	1	
Treatment Machine Commission Time		DT	1	
Spill Length		FL	1	
Degrader		FL	1	
Particles Per Spill		FL	1	

CNAO Nominal Beam Energy		DS	1	
CNAO Nominal Beam Energy Scale Factor		DS	1	
RBE Model Name		ST	1	
RBE Model Commission Time		DT	1	
Block Milling Tool Diameter		FL	1	
Internal Treatment Machine Name		SH	1	
Tomo IDMS Beam ID		LO	1	
RBE Cell Type Name		ST	1	
Tomo IDMS Machine ID		LO	1	
Tomo Localization point		DS	3	
Tomo Beam Isocenter		DS	3	
Tomo intended back jaw position		DS	1	
Tomo beam revision		IS	1	
Tomo machine revision		IS	1	
Tomo intended front jaw position		DS	1	
Beam Dose Specification Point Name		ST	1	
ROI Contour Sequence		SQ	1	
Contour Name		LO	1	
Contour Description		ST	1	
Contour Generation Algorithm		CS	1	
Contour Interpreted Type		CS	1	
Contour Display Color		IS	3	
Number of Contour Points		IS	1	
Contour Data		OF	2-2n	
Number of Pixel Contour Points		IS	1	
Pixel Contour Data		OF	2-2n	
Referenced Structure Set SOP Instance UID		UI	1	
Referenced ROI Number		IS	1	
Reference Beam Data ID		SH	1	
Reference Depth in Water		DS	1	
Mitsubishi Reference Value		DS	1	
Contour Data Sequence		SQ	1	
Tomo Plan SOP Instance UID		UI	1	
Tomo Projection Time		DS	1	

Planned Blood Boron Concentration		DS	1	
ROI Index Pixel Data		OB	1	
Applicator Thickness		DS	1	
ROI Index Bits Allocated		US	1	
Number Of Aperture Contour Points		IS	1	
Aperture Contour Point Map		FL	2-2n	
Collimation Mode		CS	1	
Treatment Machine RBE Mode		CS	1	
RaySearch Checksum Algorithm Version		LO	1	
RaySearch Checksum Data		OB	1	
Tracking ROI Sequence		SQ	1	
Imaging Angles		DS	1-6	
Mevion SMC Output		OB	1	
Mevion SMC Version		LO	1	
Tomo Couch Insertion Position Y		DS	1	
Consent Expiry Date		DA	1	
Tomo Projection Time FL		FL	1	
Siemens IONTRIS Private Creator	(300B,0010)	LO	1	SIEMENS MED SYNGO RT
Plan Type		CS	1	
Checksum Encryption Code		OB	1	
Configuration Baseline		OB	1	
Fraction Group Code		CS	1	
Fraction Dose to Dose Reference		DS	1	
Imager Organ Program		LO	1	
Imaging Technique		CS	1	
Imager Angular Angle		DS	1	
Imager Iso-centric Angle		DS	1	
Imager Orbital Angle		DS	1	
Imager Vertical Position		DS	1	
Imager Longitudinal Position		DS	1	
Imager Lateral Position		DS	1	
Imager Angular Rotation Direction		DS	1	
Imager Isocentric Rotation Direction		DS	1	
Imager Orbital Rotation Direction		DS	1	
Referenced Treatment Beam Number		IS	1	

Alternative Treatment Machine Name Sequence		SQ	1	
TomoTherapy DD Private Creator	{300D,0020}	LO	1	TOMO_DD_01
Procedure Number		SH	1	
Procedure Type		CS	1	
Procedure Purpose		CS	1	
Detector channel count		IS	1	
Dataset length		IS	1	
Detector data scaling factor		DS	1	
Compression Type		CS	1	
Compression Factor		IS	1	
Active Projection Beam Sequence		SQ	1	
Detector unit of measure		CS	1	
Element width		DS	1	
Curvature radius		DS	1	
Surface to Axis Distance		DS	1	
Surface to Center Distance		DS	1	
Channel range		IS	2	
TomoTherapy HA Private Creator	{300D,0010}	LO	1	TOMO_HA_01
Tomo Structure Blocking		CS	1	
Tomo Overlap Precedence		IS	1	
Tomo Modulation Factor		DS	1	
Tomo Target Minimum Dose Penalty		IS	1	
Tomo Target Maximum Dose Penalty		IS	1	
Tomo Organ At Risk Maximum Dose Penalty		IS	1	
Tomo Organ At Risk Full Volume Dose Penalty		IS	1	
Tomo Primary Prescription Type		CS	1	
Tomo Gantry Period		DS	1	
Tomo Machine UID		UI	1	
Tomo Treatment Pitch		DS	1	
Tomo Couch Speed		DS	1	
Tomo Performed Work Items		CS	3	
Tomo Plan Geometry		CS	1	
Tomo Projection Sinogram Data		DS	64	
Tomo Reference Isocenter		DS	3	



Tomo Registration Translations		DS	3	
Tomo Registration Rotations		DS	3	
Tomo Delivered Secondary Meterset		DS	1	
Tomo Plan IVDT Sequence		SQ	1	
Tomo IVDT HU Entry		DS	1	
Tomo IVDT Density Entry		DS	1	
Varian Medical Systems	{3243,0010}	LO	1	Varian Medical Systems VISION 3243
Beam Secondary Name		LO	1	
Related Referenced Beam Number		IS	1	



## CONTACT INFORMATION



**RaySearch Laboratories AB (publ)**  
Sveavägen 44  
SE-111 34 Stockholm  
Sweden

### Contact details head office

P.O. Box 3297  
SE-103 65 Stockholm, Sweden  
Phone: +46 8 510 530 00  
Fax: +46 8 510 530 30  
[info@raysearchlabs.com](mailto:info@raysearchlabs.com)  
[www.raysearchlabs.com](http://www.raysearchlabs.com)

### RaySearch Americas

Phone: +1 877 778 3849

### RaySearch France

Phone: +33 975 433 632

### RaySearch Korea

Phone: +82 10 2230 2046

### RaySearch Belgium

Phone: +32 2 213 83 65

### RaySearch Germany

Phone: +49 30 89 36 06 90

### RaySearch Singapore

Phone: +65 81 28 59 80

### RaySearch China

Phone: +86 137 0111 5932

### RaySearch Japan

Phone: +81 3 4405 6902

### RaySearch UK

Phone: +44 7508 426 563