

# RAYSTATION 11A

DICOM Conformance Statement Accuray Driver



RayStation

11A

### *Declaration of conformity*



Complies with 93/42/EEC Medical Device Directive as amended by M1 to M5. A copy of the corresponding Declaration of Conformity is available on request.

### *Copyright*

This document contains proprietary information that is protected by copyright. No part of this document may be photocopied, reproduced or translated to another language without prior written consent of RaySearch Laboratories AB (publ).

All Rights Reserved. © 2021, RaySearch Laboratories AB (publ).

### *Printed material*

Hard copies of Instructions for Use and Release Notes related documents are available upon request.

### *Trademarks*

RayAdaptive, RayAnalytics, RayBiology, RayCare, RayCloud, RayCommand, RayData, RayIntelligence, RayMachine, RayOptimizer, RayPACS, RayPlan, RaySearch, RaySearch Laboratories, RayStation, RayStore, RayTreat, RayWorld and the RaySearch Laboratories logotype are trademarks of RaySearch Laboratories AB (publ)\*.

Third-party trademarks as used herein are the property of their respective owners, which are not affiliated with RaySearch Laboratories AB (publ).

RaySearch Laboratories AB (publ) including its subsidiaries is hereafter referred to as RaySearch.

\* Subject to registration in some markets.

# 1 OVERVIEW

This document specifies the DICOM interface for the treatment management system (TMS) RayTreat Accuray driver with Accuray TomoTherapy and Radixact Treatment Delivery Devices(TDD).

## 1.1 NETWORK SERVICES

SOP Class Name	SOP Class UID	Provider of Service (SCP)	User of Service (SCU)
Transfer			
CT Image Storage	1.2.840.10008.5.1.4.1.1.2	Yes	No
RT Beams Treatment Record Storage	1.2.840.10008.5.1.4.1.1.481.4	Yes	No
Spatial Registration (REG) Storage	1.2.840.10008.5.1.4.1.1.66.1	Yes	No
Query/Retrieve			
Study Root Query/Retrieve Information Model – FIND	1.2.840.10008.5.1.4.1.2.2.1	Yes	No
Workflow Management			
Unified Procedure Step - Push SOP Class-Trial (Retired)	1.2.840.10008.5.1.4.34.4.1	Yes	No
Unified Procedure Step - Pull SOP Class-Trial (Retired)	1.2.840.10008.5.1.4.34.4.3	Yes	No
Verification			
Verification SOP Class	1.2.840.10008.1.1	Yes	No

# TABLE OF CONTENTS

OVERVIEW	2
NETWORK SERVICES	2
TABLE OF CONTENTS	3
INTRODUCTION	4
REVISION HISTORY	4
AUDIENCE	5
REMARKS	5
TERMS AND DEFINITIONS	5
BASICS OF DICOM COMMUNICATION	6
ABBREVIATIONS	7
REFERENCES	7
NETWORKING	8
IMPLEMENTATION MODEL	8
AE SPECIFICATIONS:	9
MEDIA INTERCHANGE	13
TRANSFORMATION OF DICOM TO CDA	14
SUPPORT OF CHARACTER SETS	15
SECURITY	16
SECURITY PROFILES	16
ASSOCIATION LEVEL SECURITY	16
APPLICATION LEVEL SECURITY	16
ANNEXES	17
IOD CONTENTS	17
DATA DICTIONARY OF PRIVATE ATTRIBUTES	33
CODE TERMINOLOGY AND TEMPLATES	34
GRAYSCALE IMAGE CONSISTENCY	34
STANDARD EXTENDED/SPECIALIZED/PRIVATE SOP CLASSES	34
PRIVATE TRANSFER SYNTAXES	34

## 3 INTRODUCTION

### 3.1 REVISION HISTORY

Date	Version	Comment
2018-12-13	1.0	Accuray Driver DCS for RayStation Release 8B
2019-02-14	2.0	Accuray Driver DCS for RayStation Release 8B SP1
2019-05-24	1.0	Accuray Driver DCS for RayStation Release 9A
2019-12-03	1.0	Accuray Driver DCS for RayStation Release 9B
2020-03-04	2.0	Accuray Driver DCS for RayStation Release 9B SP1
2020-05-11	1.0	Accuray Driver DCS for RayStation Release 10A
2020-11-27	1.0	Accuray Driver DCS for RayStation Release 10B
2021-05-12	1.0	Accuray Driver DCS for RayStation Release 11A

#### 3.1.1 TOMO Driver changelog

##### 3.1.1.1 Updates between 2.2.0.0 (RayStation 8B) - 2.3.0.0 (RayStation 8B SP1)

- If an exception is thrown when processing the request queue, the driver application will now be taken down. It should then be automatically restarted by the Windows service manager. Previously, any further requests would instead fail and the driver would have to be restarted manually.
- Validation has been added to ensure that for plans with multiple isocenters, the differences in the delivery positions for each beam are as expected. This is done by comparing to the relative positions of the isocenters for each beam in the plan to the table top displacements in the patient setup sequence and the absolute table top positions in the control point sequence.

##### 3.1.1.2 Updates between 2.3.0.0 (RayStation 8B SP1) - 3.0.0.0 (RayStation 9A)

- More detail in logging of plan validation results.
- Adding validation to make it impossible to have a second drivers connected to the same database as the first driver if the drivers has somehow been severely misconfigured.

##### 3.1.1.3 Updates between 3.0.0.0 (RayStation 9A) - 9.1.0 (RayStation 9B)

- Uses the same version numbering as the corresponding RayStation release.
- Validates the checksum that are produced by RayStation DICOM export when plan is created by a RayStation version above 9.1.0.0 (RayStation 9B) and when the plan has DICOM attribute `Manufacturer` set to `RaySearch Laboratories`.
- No longer validates that the plan has been asked for when setting session in progress.
- No longer requires machine to exist in MachineDB to be able to create tolerance tables.
- Automatic configuration of driver features, synchronizable through Clinic Settings.
- Logs more session information when first receiving it on the driver.
- Displays progress percentage in RayTreat while session is in progress.
- UPS and BDI is no longer sent as `CONTINUATION` if no meterset has been delivered, even if a treatment record has been received in a previous session for that fraction.
- Enables parsing of treatment records outside of a session context for offline recording. Validating the `Treatment Machine Name`, `Current Fraction Number`, `Patient ID`, `Patient's Name`, `Patient's Sex`, `Patient's Birth Date` and the `Referenced SOP Instance UID` of the `Referenced RT Plan`, in the Treatment Record.
- Now prefers the `Study Instance UID` of the Treatment Record to be the same as the plan that was delivered, also prefers the `Referenced Series Sequence` of the treatment record to contain a reference to the plan that was delivered. If set, this will enable better possibilities for offline treatment recording.
- Improved readability of DCS. Shows the correct indentation of attributes inside a sequence. Removes all attributes from the "Created SOP Instance(s)" chapter where the value is just read but never written. Type 1 values that are not actually read by us now has the comment "Value not read".
- No longer sends any BDI as part of the UPS (Even though one is created and sent for storage in the PACS).
- Immediately refreshes the plans `Series Instance UID` once reference is received from RayTreat. Since the `Series Instance UID` will be unknown to us because it was regenerated by iDMS. This is done by querying the PACS.

##### 3.1.1.4 Updates between 9.1.0 (RayStation 9B) - 9.2.0 (RayStation 9B SP1)

- UPS for a session with a previous treatment record of 0 delivered MU now sends its treatment records in the `Input Information Sequence` as part of the UPS.
- Table top positions are now included when delivering QA session.

##### 3.1.1.5 Updates between 9.2.0 (RayStation 9B SP1) - 10.0.0 (RayStation 10A)

- No changes affecting driver in this release.

##### 3.1.1.6 Updates between 10.0.0 (RayStation 10A) - 10.1.0 (RayStation 10B)

- Parses out the `CurrentFractionNumber` for all beams in the treatment records and always expects all to be the same.
- Introduced `IsClinical` concept to the driver. A Clinical driver can never communicate with non-clinical RayTreat. Or vice versa.
- Made connection between driver and RayTreat more secure by forcing the usage of HTTPS for a clinical driver.
- New RaaS service now handles the data synchronization with the driver (all drivers will be connected to the same RaaS service). Driver will also always send the PACS data to this service instead of to RayTreat.

### 3.1.1.7 Updates between 10.1.0 (RayStation 10B) - 11.0.0 (RayStation 11A)

- Supports creating a TDW-1 Trial (Retired) BDI and sends it as part of the UPS.
- Added feature to cache and send additional DICOM references in the UPS if such are created by RayStation. (IE: 2nd gen DICOM objects created for CyberKnife).
- Now supports receiving CMOVE requests. Supported SopClasses: RTPlan and BDI (Retired).
- CMOVE requests are still configurable and defaulted to not allowed because of late testing.
- Changes made to how a Registration is paired to the RTStruct received since the driver does not know of the Planning image used for Registration by Accuray.
- Supports receiving empty table top positions in treatment record.

## 3.2 AUDIENCE

This document is written for users that need to understand how Accuray will integrate into their healthcare facility. This includes both those responsible for overall imaging network policy and architecture, as well as integrators who need to have a detailed understanding of the DICOM features of the product. This document contains some basic DICOM definitions so that any reader may understand how this product implements DICOM features. However, integrators are expected to fully understand all the DICOM terminology, how the tables in this document relate to the product's functionality, and how that functionality integrates with other devices that support compatible DICOM features.

## 3.3 REMARKS

This document is written for users that need to understand how RayTreat Accuray Driver will integrate into their healthcare facility. This includes both those responsible for overall imaging network policy and architecture, as well as integrators who need to have a detailed understanding of the DICOM features of the product. This document contains some basic DICOM definitions so that any reader may understand how this product implements DICOM features. However, integrators are expected to fully understand all the DICOM terminology, how the tables in this document relate to the product's functionality, and how that functionality integrates with other devices that support compatible DICOM features.

This Conformance Statement is not supposed to replace validation with other DICOM equipment to ensure proper exchange of intended information. In fact, the user should be aware of the following important issues:

- The comparison of different Conformance Statements is just the first step towards assessing interconnectivity and interoperability between the product and other DICOM conformant equipment.
- Test procedures should be defined and executed to validate the required level of interoperability with specific compatible DICOM equipment, as established by the healthcare facility.

### 3.3.1 Interoperability validation needed

When using RayTreat Accuray Driver together with other software, the DICOM conformance statements must be compared and relevant validation tests run. The DICOM standard by itself does not guarantee interoperability. The Conformance Statement does, however, facilitate a first-level comparison for interoperability between different applications supporting compatible DICOM functionality. RaySearch is also active within the IHE-RO. Contact RaySearch for more info regarding adherence to IHE-RO profiles.

### 3.3.2 DICOM revision

The module tables listed in the last two chapters are based on part 3 of the DICOM-standard edition 2020a. For extra clarity all attributes in the referenced modules have been listed, even the ones that are not used by Accuray.

## 3.4 TERMS AND DEFINITIONS

Informal definitions are provided for the following terms used in this Conformance Statement. The DICOM Standard is the authoritative source for formal definitions of these terms.

**Abstract Syntax** – the information agreed to be exchanged between applications, generally equivalent to a Service/Object Pair (SOP) Class. Examples : Verification SOP Class, Modality Worklist Information Model Find SOP Class, Computed Radiography Image Storage SOP Class.

**Application Entity (AE)** – an end point of a DICOM information exchange, including the DICOM network or media interface software; i.e., the software that sends or receives DICOM information objects or messages. A single device may have multiple Application Entities.

**Application Entity Title** – the externally known name of an Application Entity, used to identify a DICOM application to other DICOM applications on the network.

**Application Context** – the specification of the type of communication used between Application Entities. Example: DICOM network protocol.

**Association** – a network communication channel set up between Application Entities.

**Attribute** – a unit of information in an object definition; a data element identified by a tag. The information may be a complex data structure (Sequence), itself composed of lower level data elements. Examples: Patient ID (0010,0020), Accession Number (0008,0050), Photometric Interpretation (0028,0004), Procedure Code Sequence (0008,1032).

**Information Object Definition (IOD)** – the specified set of Attributes that comprise a type of data object; does not represent a specific instance of the data object, but rather a class of similar data objects that have the same properties. The Attributes may be specified as Mandatory (Type 1), Required but possibly unknown (Type 2), or Optional (Type 3), and there may be conditions associated with the use of an Attribute (Types 1C and 2C). Examples: MR Image IOD, CT Image IOD, Print Job IOD.

**Joint Photographic Experts Group (JPEG)** – a set of standardized image compression techniques, available for use by DICOM applications.

**Module** – a set of Attributes within an Information Object Definition that are logically related to each other. Example: Patient Module includes Patient Name, Patient ID, Patient Birth Date, and Patient Sex.

**Negotiation** – first phase of Association establishment that allows Application Entities to agree on the types of data to be exchanged and how that data will be encoded.

**Protocol Data Unit (PDU)** – a packet (piece) of a DICOM message sent across the network. Devices must specify the maximum size packet they can receive for DICOM messages.

**Service Class Provider (SCP)** – role of an Application Entity that provides a DICOM network service; typically, a server that performs operations requested by another

Application Entity (Service Class User). Examples: Picture Archiving and Communication System (image storage SCP, and image query/retrieve SCP), Radiology Information System (modality worklist SCP).

**Service Class User (SCU)** – role of an Application Entity that uses a DICOM network service; typically, a client. Examples: imaging modality (image storage SCU, and modality worklist SCU), imaging workstation (image query/retrieve SCU)

**Service/Object Pair (SOP) Class** – the specification of the network or media transfer (service) of a particular type of data (object); the fundamental unit of DICOM interoperability specification. Examples: Ultrasound Image Storage Service, Basic Grayscale Print Management.

**Service/Object Pair (SOP) Instance** – an information object; a specific occurrence of information exchanged in a SOP Class. Examples: a specific x-ray image.

**Tag** – a 32-bit identifier for a data element, represented as a pair of four digit hexadecimal numbers, the “group” and the “element”. If the “group” number is odd, the tag is for a private (manufacturer-specific) data element. Examples: (0010,0020) [Patient ID], (07FE,0010) [Pixel Data], (0019,0210) [private data element]

**Transfer Syntax** – the encoding used for exchange of DICOM information objects and messages. Examples: JPEG compressed (images), little endian explicit value representation.

**Unique Identifier (UID)** – a globally unique “dotted decimal” string that identifies a specific object or a class of objects; an ISO-8824 Object Identifier. Examples: Study Instance UID, SOP Class UID, SOP Instance UID.

**Value Representation (VR)** – the format type of an individual DICOM data element, such as text, an integer, a person’s name, or a code. DICOM information objects can be transmitted with either explicit identification of the type of each data element (Explicit VR), or without explicit identification (Implicit VR); with Implicit VR, the receiving application must use a DICOM data dictionary to look up the format of each data element.

### 3.5 BASICS OF DICOM COMMUNICATION

This section describes terminology used in this Conformance Statement for the non-specialist. This section is not a substitute for training about DICOM, and it makes many simplifications about the meanings of DICOM terms.

Two Application Entities (devices) that want to communicate with each other over a network using DICOM protocol must first agree on several things during an initial network “handshake”. One of the two devices must initiate an Association (a connection to the other device), and ask if specific services, information, and encoding can be supported by the other device (Negotiation).

DICOM specifies a number of network services and types of information objects, each of which is called an Abstract Syntax for the Negotiation. DICOM also specifies a variety of methods for encoding data, denoted Transfer Syntaxes. The Negotiation allows the initiating Application Entity to propose combinations of Abstract Syntax and Transfer Syntax to be used on the Association; these combinations are called Presentation Contexts. The receiving Application Entity accepts the Presentation Contexts it supports.

For each Presentation Context, the Association Negotiation also allows the devices to agree on Roles – which one is the Service Class User (SCU - client) and which is the Service Class Provider (SCP - server). Normally the device initiating the connection is the SCU, i.e., the client system calls the server, but not always.

The Association Negotiation finally enables exchange of maximum network packet (PDU) size, security information, and network service options (called Extended Negotiation information). The Application Entities, having negotiated the Association parameters, may now commence exchanging data. Common data exchanges include queries for worklists and lists of stored images, transfer of image objects and analyses (structured reports), and sending images to film printers. Each exchangeable unit of data is formatted by the sender in accordance with the appropriate Information Object Definition, and sent using the negotiated Transfer Syntax. There is a Default Transfer Syntax that all systems must accept, but it may not be the most efficient for some use cases. Each transfer is explicitly acknowledged by the receiver with a Response Status indicating success, failure, or that query or retrieve operations are still in process.

Two Application Entities may also communicate with each other by exchanging media (such as a CD-R). Since there is no Association Negotiation possible, they both use a Media Application Profile that specifies “pre-negotiated” exchange media format, Abstract Syntax, and Transfer Syntax.

### 3.6 ABBREVIATIONS

Name	Meaning
AE	Application Entity
CT	Computed Tomography
DICOM	Digital Imaging and Communications in Medicine
IHE / IHE-RO	Integrating the Healthcare Enterprise. IHE-RO deals with integrating Radiation Oncology.
IOD	Information Object Definition
JPEG	Joint Photographic Experts Group
MR	Magnetic Resonance Imaging
PACS	Picture Archiving and Communication System
PET	Positron Emission Tomography
PTS	Proton Planning System (used by IBA)
RT	Radiotherapy
SCP	Service Class Provider
SCU	Service Class User
SOP	Service-Object Pair
TDD	Treatment Delivery Device
TMS	Treatment Management System
TPS	Treatment Planning System

### 3.7 REFERENCES

- NEMA PS3 Digital Imaging and Communications in Medicine (DICOM) Standard, available free at <http://medical.nema.org/>




## 4 NETWORKING

The following diagram illustrates the application dataflow between RayTreat Accuray Driver and the specific TDD.

### 4.1 IMPLEMENTATION MODEL

#### 4.1.1 Application data flow

The following diagram illustrates the application data flow between RayTreat Accuray Driver and the specific TDD

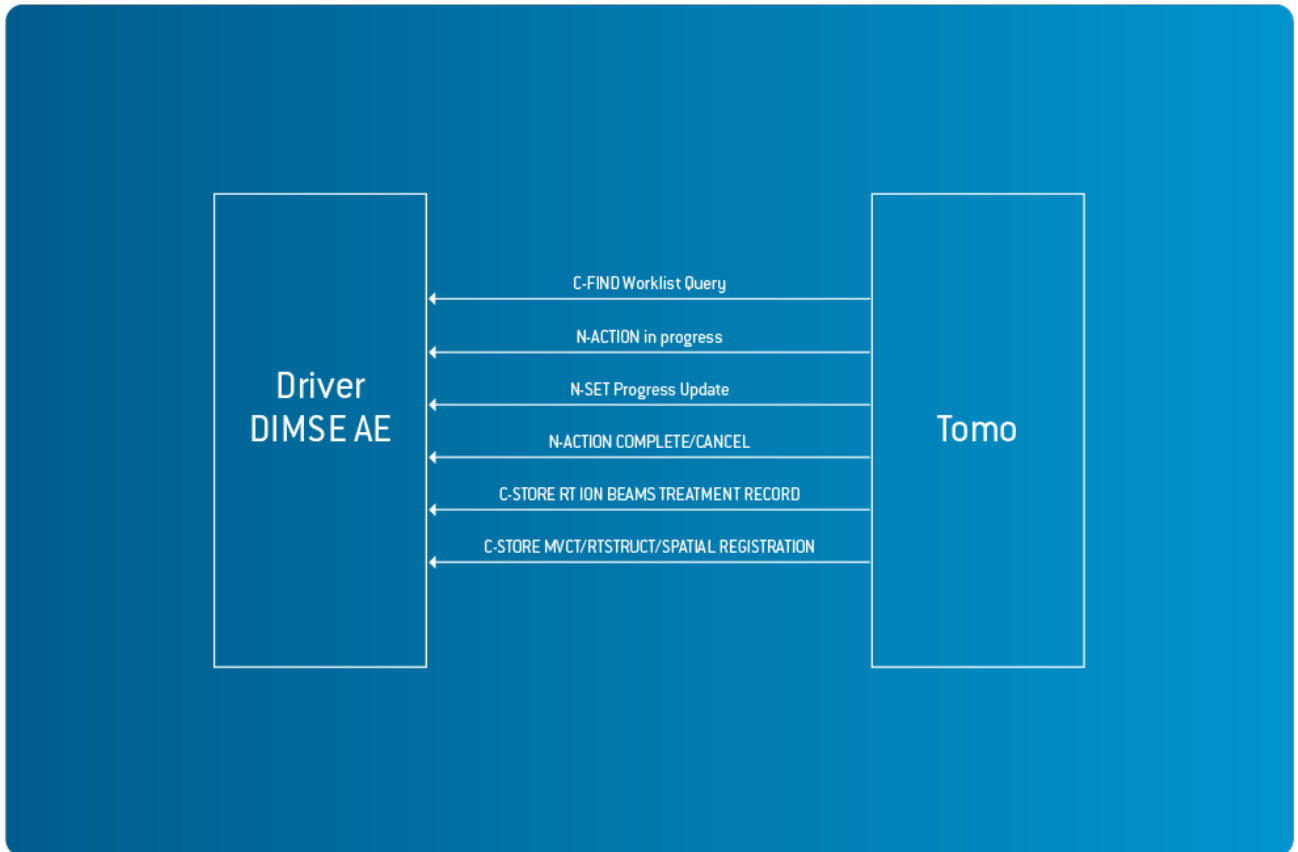

*The scenario starts with a C-FIND query for Unified Procedure Steps from the remote client. The client can then take responsibility for the UPS by setting it to IN PROGRESS. Once the UPS is IN PROGRESS the RayTreat Tomo Driver only allows requests corresponding to the current session until the session has been completed in the application.*

#### 4.1.2 Functional Definition of AEs

##### 4.1.2.1 Functional Definition of "Accuray Application Entity"

The following operations are supported:

###### CT Image

- C-STORE for setup CT images.
- C-MOVE for planning CT images

###### Spatial Registration (SR0)

- C-STORE for registration between setup and planning images.

###### RT Image

- C-STORE for setup RT Images.

###### RT Structure Set

- C-MOVE for planning RT Structure Set.

###### RT Beams Treatment Record

- C-STORE for delivery result.
- C-MOVE for previous delivered results.

###### Beams Delivery Instructions

- C-MOVE for BDIs related to the Unified Procedure Step.

**Modality Performed Procedure Step - PULL**

- C-FIND for worklist query.
- N-ACTION for UPS status changes.
- N-SET for progress update.

**Verification**

- C-ECHO for connection verification

**4.1.3 Sequence of Real World Activities**

**4.1.3.1 Prepare session**

Once the patient is checked in to the session, Unified Procedure Steps will be created and available for Worklist queries.

**4.1.3.2 Manual cancellation**

The procedure step can be canceled by the user in the application. Further requests relation to the session will be rejected.

**4.1.3.3 Complete session**

All sessions, including canceled sessions, needs to be completed by the user in the application before another session can be started.

**4.2 AE SPECIFICATIONS:**

**4.2.1 RayTreat Accuray Driver Application Entity**

**4.2.1.1 SOP Classes**

SOP Class Name	SOP Class UID	Provider of Service (SCP)	User of Service (SCU)
Transfer			
CT Image Storage	1.2.840.10008.5.1.4.1.1.2	Yes	No
RT Beams Treatment Record Storage	1.2.840.10008.5.1.4.1.1.481.4	Yes	No
Spatial Registration (REG) Storage	1.2.840.10008.5.1.4.1.1.66.1	Yes	No
Query/Retrieve			
Study Root Query/Retrieve Information Model – FIND	1.2.840.10008.5.1.4.1.2.2.1	Yes	No
Workflow Management			
Unified Procedure Step - Push SOP Class-Trial (Retired)	1.2.840.10008.5.1.4.34.4.1	Yes	No
Unified Procedure Step - Pull SOP Class-Trial (Retired)	1.2.840.10008.5.1.4.34.4.3	Yes	No
Verification			
Verification SOP Class	1.2.840.10008.1.1	Yes	No

**4.2.1.2 Association Policies**

Not applicable

**4.2.1.3 General**

The DICOM standard Application context shall be specified.

Application Context Name	1.2.840.10008.3.1.1.1
--------------------------	-----------------------

The maximum PDU size in not configurable and is set to 16384 for SCU and unlimited for SCP.

**4.2.1.4 Number of Associations**

Any number of incoming concurrent associations are accepted.

**4.2.1.5 Asynchronous Nature**

RayTreat Accuray Driver does not support asynchronous communication (multiple outstanding transactions over a single Association).

**4.2.1.6 Implementation Identity Information**

Not applicable

4.2.1.7 Association Initiation Policy

The implementation for this Application Entity is:

Implementation Class UID	1.3.6.1.4.1.30071.8
Implementation Version Name	fo-dicom-raysearch 4.0.4 (based on official fo-dicom 4.0.2)

4.2.1.8 Activity C-ECHO

4.2.1.8.1 Description and Sequencing of Activities

A C-ECHO request can always be sent to the Accuray driver.

4.2.1.8.2 Accepted Presentation Context

Abstract Syntax Name	Abstract Syntax UID	Transfer Syntax name	Transfer Syntax UID	Role	Extended Negotiation
Verification	1.2.840.10008.1.1	Implicit VR Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR Little Endian	1.2.840.10008.1.2.1		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		

4.2.1.9 Activity C-FIND

4.2.1.9.1 Description and Sequencing of Activities

Not applicable

4.2.1.9.2 Accepted Presentation Context

Abstract Syntax Name	Abstract Syntax UID	Transfer Syntax name	Transfer Syntax UID	Role	Extended Negotiation
Unified Procedure Step - Pull SOP Class - Trial (Retired)	1.2.840.10008.5.1.4.34.4.3	Implicit VR Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR Little Endian	1.2.840.10008.1.2.1		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		

4.2.1.10 Activity C-MOVE

4.2.1.10.1 Description and Sequencing of Activities

Not applicable

4.2.1.10.2 Accepted Presentation Context

RayTreat Accuray Driver does not support C-MOVE requests.

4.2.1.11 Activity C-STORE

4.2.1.11.1 Description and Sequencing of Activities

Not applicable

4.2.1.11.2 Accepted Presentation Context

Abstract Syntax Name	Abstract Syntax UID	Transfer Syntax name	Transfer Syntax UID	Role	Extended Negotiation
CT Image Storage	1.2.840.10008.5.1.4.1.1.2	Implicit VR Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR Little Endian	1.2.840.10008.1.2.1		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
RT Structure Set Storage	1.2.840.10008.5.1.4.1.1.481.3	Implicit VR Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR Little Endian	1.2.840.10008.1.2.1		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
Spatial Registration Storage	1.2.840.10008.5.1.4.1.1.66.1	Implicit VR Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR Little Endian	1.2.840.10008.1.2.1		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		

RT Beams Treatment Record Storage	1.2.840.10008.5.1.4.1.1.481.4	Implicit VR Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR Little Endian	1.2.840.10008.1.2.1		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		

4.2.1.11.3 Status Response

Service Status	Further meaning	Error Code	Reason
Failure	Storage Cannot Understand	Cxxx	Cannot find session or validation failed.
	SOP class not supported	0122	SOP class not supported.
Success	Success	0000	

4.2.1.12 Activity N-ACTION

4.2.1.12.1 Description and Sequencing of Activities

Not applicable

4.2.1.12.2 Accepted Presentation Context

Abstract Syntax Name	Abstract Syntax UID	Transfer Syntax name	Transfer Syntax UID	Role	Extended Negotiation
Unified Procedure Step - Pull SOP Class - Trial (Retired)	1.2.840.10008.5.1.4.34.4.3	Implicit VR Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR Little Endian	1.2.840.10008.1.2.1		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		

4.2.1.12.3 Status Response

Service Status	Further meaning	Error Code	Reason
Refused	QueryRetrieveUnableToProcess	C000	Unknown SOP Instance UID.
	NoLongerUpdateUps	C300	The UPS may no longer be updated.
	IncorrectTransactionUid	C301	The correct Transaction UID was not provided.
	AlreadyInProgress	C302	The UPS is already IN PROGRESS
	IncorrectTransactionUid	C301	The correct Transaction UID was not provided.
	SopInstanceUidDoesNotExist	C307	Specified SOP Instance UID does not exist or is nota UPS Instance managed by this SCP
Failure	AlreadyCanceled	B304	The UPS is already in the requested state of CANCELED.
	AlreadyCompleted	B306	The UPS is already in the requested state of COMPLETED.
Success	Success	0000	

4.2.1.13 Activity N-SET

4.2.1.13.1 Description and Sequencing of Activities

Not applicable

4.2.1.13.2 Accepted Presentation Context

Abstract Syntax Name	Abstract Syntax UID	Transfer Syntax name	Transfer Syntax UID	Role	Extended Negotiation
Unified Procedure Step - Pull SOP Class - Trial (Retired)	1.2.840.10008.5.1.4.34.4.3	Implicit VR Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR Little Endian	1.2.840.10008.1.2.1		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		

4.2.1.13.3 Status Response

Service Status	Further meaning	Error Code	Reason
Refused	QueryRetrieveUnableToProcess	C000	Unknown SOP Instance UID.

NETWORKING

	NoLongerUpdateUps	C300	The UPS may no longer be updated.
	IncorrectTransactionUid	C301	The correct Transaction UID was not provided.
	SopInstanceUidDoesNotExists	C307	Specified SOP Instance UID does not exist or is nota UPS Instance managed by this SCP
Success	Success	0000	

## 5 MEDIA INTERCHANGE

Not applicable

## 6 TRANSFORMATION OF DICOM TO CDA

Not applicable

## 7 SUPPORT OF CHARACTER SETS

RayTreat Accuray Driver support the following charactersets in addition to the default

- ISO\_IR 192



## 8 SECURITY

### 8.1 SECURITY PROFILES

No Security Profiles are supported.

### 8.2 ASSOCIATION LEVEL SECURITY

RayTreat Accuray Driver checks the following values for validation of received Association Open Requests:

- Called AE Title.

### 8.3 APPLICATION LEVEL SECURITY

None supported.

## 9 ANNEXES

### 9.1 IOD CONTENTS

#### 9.1.1 Created SOP Instance(s)

##### 9.1.1.1 RT Beams Treatment Record IOD

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	Yes
Series	RT Series Module	Yes
Equipment	General Equipment Module	Yes
Treatment Record	RT General Treatment Record Module	Yes
	RT Patient Setup Module	Yes
	RT Treatment Machine Record Module	Yes
	RT Beams Session Record Module	Yes
	SOP Common Module	Yes

##### 9.1.1.1.1 Patient Module

Attribute name	Tag	Vr	Type	Comment
Patient's Name	(0010,0010)	PN	2	
Patient ID	(0010,0020)	LO	2	
Patient's Birth Date	(0010,0030)	DA	2	
Patient's Sex	(0010,0040)	CS	2	

##### 9.1.1.1.2 General Study Module

Attribute name	Tag	Vr	Type	Comment
Study Instance UID	(0020,000D)	UI	1	
Study Date	(0008,0020)	DA	2	
Study Time	(0008,0030)	TM	2	
Referring Physician's Name	(0008,0090)	PN	2	
Study ID	(0020,0010)	SH	2	
Study Description	(0008,1030)	LO	3	

##### 9.1.1.1.3 RT Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	(0008,0060)	CS	1	Can only be updated by the driver to: RTRECORD.
Series Instance UID	(0020,000E)	UI	1	
Series Number	(0020,0011)	IS	2	
Operators' Name	(0008,1070)	PN	2	

##### 9.1.1.1.4 General Equipment Module

Attribute name	Tag	Vr	Type	Comment
Manufacturer	(0008,0070)	LO	2	
Manufacturer's Model Name	(0008,1090)	LO	3	
Software Versions	(0018,1020)	LO	3	

##### 9.1.1.1.5 RT General Treatment Record Module

Attribute name	Tag	Vr	Type	Comment
Instance Number	(0020,0013)	IS	1	
Treatment Date	(3008,0250)	DA	2	
Treatment Time	(3008,0251)	TM	2	
Referenced RT Plan Sequence	(300C,0002)	SQ	2	
>Referenced SOP Class UID	(0008,1150)	UI	1	
>Referenced SOP Instance UID	(0008,1155)	UI	1	
Referenced Treatment Record Sequence	(3008,0030)	SQ	3	
>Referenced SOP Class UID	(0008,1150)	UI	1	
>Referenced SOP Instance UID	(0008,1155)	UI	1	

## 9.1.1.1.6 RT Patient Setup Module

Attribute name	Tag	Vr	Type	Comment
Patient Setup Sequence	(300A,0180)	SQ	1	
>Patient Setup Number	(300A,0182)	IS	1	
>Patient Position	(0018,5100)	CS	1C	
>Setup Technique	(300A,01B0)	CS	3	

## 9.1.1.1.7 RT Treatment Machine Record Module

Attribute name	Tag	Vr	Type	Comment
Treatment Machine Sequence	(300A,0206)	SQ	1	
>Manufacturer	(0008,0070)	LO	2	
>Institution Name	(0008,0080)	LO	2	
>Manufacturer's Model Name	(0008,1090)	LO	2	
>Device Serial Number	(0018,1000)	LO	2	

## 9.1.1.1.8 RT Beams Session Record Module

Attribute name	Tag	Vr	Type	Comment
Number of Fractions Planned	(300A,0078)	IS	2	
Primary Dosimeter Unit	(300A,00B3)	CS	1	
Treatment Session Beam Sequence	(3008,0020)	SQ	1	
>Referenced Beam Number	(300C,0006)	IS	3	
>Beam Name	(300A,00C2)	LO	3	
>Beam Type	(300A,00C4)	CS	1	
>Radiation Type	(300A,00C6)	CS	1	
>Beam Limiting Device Leaf Pairs Sequence	(3008,00A0)	SQ	1	
>>RT Beam Limiting Device Type	(300A,00B8)	CS	1	
>>Number of Leaf/Jaw Pairs	(300A,00BC)	IS	1	
>Referenced Patient Setup Number	(300C,006A)	IS	3	
>Number of Wedges	(300A,00D0)	IS	1	
>Number of Compensators	(300A,00E0)	IS	2	
>Number of Boli	(300A,00ED)	IS	2	
>Number of Blocks	(300A,00F0)	IS	2	
>Current Fraction Number	(3008,0022)	IS	2	
>Treatment Delivery Type	(300A,00CE)	CS	2	
>Treatment Termination Status	(3008,002A)	CS	1	

>Treatment Verification Status	(3008,002C)	CS	2	
>Specified Primary Meterset	(3008,0032)	DS	3	
>Delivered Primary Meterset	(3008,0036)	DS	3	
>Number of Control Points	(300A,0110)	IS	1	
>Control Point Delivery Sequence	(3008,0040)	SQ	1	
>>Referenced Control Point Index	(300C,00F0)	IS	3	
>>Treatment Control Point Date	(3008,0024)	DA	1	
>>Treatment Control Point Time	(3008,0025)	TM	1	
>>Specified Meterset	(3008,0042)	DS	2	
>>Delivered Meterset	(3008,0044)	DS	1	
>>Nominal Beam Energy	(300A,0114)	DS	3	
>>Wedge Position Sequence	(300A,0116)	SQ	3	
>>>Referenced Wedge Number	(300C,00C0)	IS	1	
>>>Wedge Position	(300A,0118)	CS	1	
>>Gantry Angle	(300A,011E)	DS	1C	
>>Gantry Rotation Direction	(300A,011F)	CS	1C	
>>Beam Limiting Device Angle	(300A,0120)	DS	1C	
>>Beam Limiting Device Rotation Direction	(300A,0121)	CS	1C	
>>Patient Support Angle	(300A,0122)	DS	1C	
>>Patient Support Rotation Direction	(300A,0123)	CS	1C	
>>Table Top Pitch Angle	(300A,0140)	FL	1C	
>>Table Top Pitch Rotation Direction	(300A,0142)	CS	1C	
>>Table Top Roll Angle	(300A,0144)	FL	1C	
>>Table Top Roll Rotation Direction	(300A,0146)	CS	1C	
>>Table Top Vertical Position	(300A,0128)	DS	2C	
>>Table Top Longitudinal Position	(300A,0129)	DS	2C	
>>Table Top Lateral Position	(300A,012A)	DS	2C	

## 9.1.1.1.9 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	(0008,0016)	UI	1	
SOP Instance UID	(0008,0018)	UI	1	

## 9.1.1.2 RT Beams Delivery Instruction RETIRED IOD

IE	Module	Used
Patient	Patient Module	Yes
	Clinical Trial Subject Module	No
Study	General Study Module	Yes
	Patient Study Module	No
	Clinical Trial Study Module	No
Series	General Series Module	Yes
	Clinical Trial Series Module	No
Equipment	General Equipment Module	Yes
Plan	RT Beams Delivery Instruction Module	Yes
	Common Instance Reference Module	No

General Reference Module	No
SOP Common Module	Yes

## 9.1.1.2.1 Patient Module

Attribute name	Tag	Vr	Type	Comment
Patient's Name	(0010,0010)	PN	2	
Patient ID	(0010,0020)	LO	2	
Patient's Birth Date	(0010,0030)	DA	2	
Patient's Sex	(0010,0040)	CS	2	Can only be updated by the driver to: M, F, O.

## 9.1.1.2.2 General Study Module

Attribute name	Tag	Vr	Type	Comment
Study Instance UID	(0020,0000)	UI	1	
Study Date	(0008,0020)	DA	2	
Study Time	(0008,0030)	TM	2	
Referring Physician's Name	(0008,0090)	PN	2	
Study ID	(0020,0010)	SH	2	
Accession Number	(0008,0050)	SH	2	
Study Description	(0008,1030)	LO	3	

## 9.1.1.2.3 General Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	(0008,0060)	CS	1	Can only be updated by the driver to: PLAN.
Series Instance UID	(0020,000E)	UI	1	
Series Number	(0020,0011)	IS	2	
Series Date	(0008,0021)	DA	3	
Series Time	(0008,0031)	TM	3	
Series Description	(0008,103E)	LO	3	

## 9.1.1.2.4 General Equipment Module

Attribute name	Tag	Vr	Type	Comment
Manufacturer	(0008,0070)	LO	2	
Station Name	(0008,1010)	SH	3	

## 9.1.1.2.5 RT Beams Delivery Instruction Module

Attribute name	Tag	Vr	Type	Comment
Beam Task Sequence	(0074,1020)	SQ	1	
>Beam Task Type	(0074,1022)	CS	1	Can only be updated by the driver to: <ul style="list-style-type: none"> <li>TREAT - Treat</li> </ul>
>Treatment Delivery Type	(300A,00CE)	CS	1	Can only be updated by the driver to: TREATMENT, CONTINUATION.
>Continuation Start Meterset	(0074,0120)	FD	1C	
>Continuation End Meterset	(0074,0121)	FD	1C	
>Current Fraction Number	(3008,0022)	IS	1	
>Referenced Beam Number	(300C,0006)	IS	1	
>Table Top Vertical Setup Displacement	(300A,01D2)	DS	2	
>Table Top Longitudinal Setup Displacement	(300A,01D4)	DS	2	
>Table Top Lateral Setup Displacement	(300A,01D6)	DS	2	

>Referenced RT Plan Sequence	(300C,0002)	SQ	3	
>>Referenced Series Sequence	(0008,1115)	SQ	1C	
>>>Series Instance UID	(0020,000E)	UI	1	
>>>Retrieve AE Title	(0008,0054)	AE		
>>>Referenced SOP Sequence	(0008,1199)	SQ	1	
>>>>Referenced SOP Class UID	(0008,1150)	UI	1	
>>>>Referenced SOP Instance UID	(0008,1155)	UI	1	
>>Study Instance UID	(0020,000D)	UI	1	
Omitted Beam Task Sequence	(300C,0111)	SQ	3	
>Referenced Beam Number	(300C,0006)	IS	1	
>Reason for Omission	(300C,0112)	CS	1	Can only be updated by the driver to: ALREADY_TREATED.

## 9.1.1.2.6 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	(0008,0016)	UI	1	
SOP Instance UID	(0008,0018)	UI	1	
Specific Character Set	(0008,0005)	CS	1C	
Instance Creation Date	(0008,0012)	DA	3	
Instance Creation Time	(0008,0013)	TM	3	

## 9.1.1.3 Unified Procedure Step RETIRED IOD

IE	Module	Used
Unified Procedure Step	SOP Common Module	Yes
	Unified Procedure Step Relationship Module	Yes
	Unified Procedure Step Scheduled Procedure Information Module	Yes
	Unified Procedure Step Progress Information Module	Yes
	Unified Procedure Step Performed Procedure Information Module	No
	Patient Demographic Module	No
	Patient Medical Module	No
	Visit Identification Module	No
	Visit Status Module	No
	Visit Admission Module	No
	Transaction Module	Yes

## 9.1.1.3.1 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	(0008,0016)	UI	1	
SOP Instance UID	(0008,0018)	UI	1	
Timezone Offset From UTC	(0008,0201)	SH	3	

## 9.1.1.3.2 Unified Procedure Step Relationship Module

Attribute name	Tag	Vr	Type	Comment
Patient's Name	(0010,0010)	PN		
Patient ID	(0010,0020)	LO		
Patient's Birth Date	(0010,0030)	DA		

Patient's Sex	(0010,0040)	CS		Can only be updated by the driver to:
				<ul style="list-style-type: none"> <li>• M - Male</li> <li>• F - Female</li> <li>• 0 - Other</li> </ul>
Referenced Request Sequence	(0040,A370)	SQ		
>Requested Procedure Code Sequence	(0032,1064)	SQ		
>>Code Value	(0008,0100)	SH	1C	
>>Coding Scheme Designator	(0008,0102)	SH	1C	
>>Coding Scheme Version	(0008,0103)	SH	1C	
>>Code Meaning	(0008,0104)	LO	1	
>Requested Procedure ID	(0040,1001)	SH		

## 9.1.1.3.3 Unified Procedure Step Scheduled Procedure Information Module

Attribute name	Tag	Vr	Type	Comment
Scheduled Procedure Step Priority	(0074,1200)	CS		Can only be updated by the driver to:
				<ul style="list-style-type: none"> <li>• HIGH - High</li> <li>• MEDIUM - Medium</li> <li>• LOW - Low</li> </ul>
Procedure Step Label	(0074,1204)	LO		
Scheduled Station Name Code Sequence	(0040,4025)	SQ		
>Code Value	(0008,0100)	SH	1C	
>Coding Scheme Designator	(0008,0102)	SH	1C	
>Code Meaning	(0008,0104)	LO	1	
Scheduled Procedure Step Start DateTime	(0040,4005)	DT		
Expected Completion DateTime	(0040,4011)	DT		
Scheduled Workitem Code Sequence	(0040,4018)	SQ		
>Code Value	(0008,0100)	SH	1C	
>Coding Scheme Designator	(0008,0102)	SH	1C	
>Code Meaning	(0008,0104)	LO	1	
Scheduled Processing Parameters Sequence	(0074,1210)	SQ		
>Value Type	(0040,A040)	CS	1	Can only be updated by the driver to:
				<ul style="list-style-type: none"> <li>• TEXT - Text</li> </ul>
>Concept Name Code Sequence	(0040,A043)	SQ	1	
>>Code Value	(0008,0100)	SH	1C	
>>Coding Scheme Designator	(0008,0102)	SH	1C	
>>Code Meaning	(0008,0104)	LO	1	
>Text Value	(0040,A160)	UT	1C	
Input Information Sequence	(0040,4021)	SQ		
>Study Instance UID	(0020,0000)	UI	1	
>Referenced Series Sequence	(0008,1115)	SQ	1C	
>>Series Instance UID	(0020,000E)	UI	1	
>>Retrieve AE Title	(0008,0054)	AE		
>>Referenced SOP Sequence	(0008,1199)	SQ	1	
>>>Referenced SOP Class UID	(0008,1150)	UI	1	

>>>Referenced SOP Instance UID	(0008,1155)	UI	1	
Study Instance UID	(0020,0000)	UI		
Input Availability Flag	(0040,4020)	CS	1	Can only be updated by the driver to: <ul style="list-style-type: none"> <li>COMPLETE - Complete</li> <li>INCOMPLETE - Incomplete</li> </ul>

## 9.1.1.3.4 Unified Procedure Step Progress Information Module

Attribute name	Tag	Vr	Type	Comment
Procedure Step State	(0074,1000)	CS		Can only be updated by the driver to: <ul style="list-style-type: none"> <li>SCHEDULED - Scheduled</li> <li>IN PROGRESS - In Progress</li> <li>CANCELED - Canceled</li> <li>COMPLETED - Completed</li> </ul>

## 9.1.1.3.5 Transaction Module

Attribute name	Tag	Vr	Type	Comment
Transaction UID	(0008,1195)	UI	3	

## 9.1.2 Usage of Attributes From Received IODs

## 9.1.2.1 CT Image IOD

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	Yes
Series	General Series Module	Yes
Frame of Reference	Frame of Reference Module	Yes
Equipment	General Equipment Module	No
Image	General Image Module	No
	Image Plane Module	No
	Image Pixel Module	No
	Contrast/Bolus Module	No
	CT Image Module	Yes
	Multi-energy CT Image Module	No
	SOP Common Module	Yes

## 9.1.2.1.1 Patient Module

Attribute name	Tag	Vr	Type	Comment
Patient's Name	(0010,0010)	PN	2	
Patient ID	(0010,0020)	LO	2	
Patient's Birth Date	(0010,0030)	DA	2	
Patient's Sex	(0010,0040)	CS	2	

## 9.1.2.1.2 General Study Module

Attribute name	Tag	Vr	Type	Comment
Study Instance UID	(0020,0000)	UI	1	

## 9.1.2.1.3 General Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	(0008,0060)	CS	1	Value not read



Series Instance UID	(0020,000E)	UI	1	
Series Date	(0008,0021)	DA	3	
Series Time	(0008,0031)	TM	3	
Series Description	(0008,103E)	LO	3	
Patient Position	(0018,5100)	CS	2C	

## 9.1.2.1.4 Frame of Reference Module

Attribute name	Tag	Vr	Type	Comment
Frame of Reference UID	(0020,0052)	UI	1	

## 9.1.2.1.5 CT Image Module

Attribute name	Tag	Vr	Type	Comment
Image Type	(0008,0008)	CS	1	Value not read
Samples per Pixel	(0028,0002)	US	1	Value not read
Photometric Interpretation	(0028,0004)	CS	1	Value not read
Bits Allocated	(0028,0100)	US	1	Value not read
Bits Stored	(0028,0101)	US	1	Value not read
High Bit	(0028,0102)	US	1	Value not read
Rescale Intercept	(0028,1052)	DS	1	Value not read
Rescale Slope	(0028,1053)	DS	1	Value not read
KVP	(0018,0060)	DS	2	If larger than 1000 the image received will be registered as an MVCT instead of as an CBCT.
Table Height	(0018,1130)	DS	3	
Patient Support Angle	(300A,0122)	DS	3	
Table Top Pitch Angle	(300A,0140)	FL	3	
Table Top Roll Angle	(300A,0144)	FL	3	
Table Top Longitudinal Position	(300A,0129)	DS	3	
Table Top Lateral Position	(300A,012A)	DS	3	

## 9.1.2.1.6 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	(0008,0016)	UI	1	
SOP Instance UID	(0008,0018)	UI	1	

## 9.1.2.2 RT Structure Set IOD

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	Yes
Series	RT Series Module	Yes
Equipment	General Equipment Module	No
Structure Set	Structure Set Module	Yes
	ROI Contour Module	Yes
	RT ROI Observations Module	Yes
	SOP Common Module	Yes

## 9.1.2.2.1 Patient Module

Attribute name	Tag	Vr	Type	Comment
----------------	-----	----	------	---------

Patient's Name	(0010,0010)	PN	2	
Patient ID	(0010,0020)	LO	2	
Patient's Birth Date	(0010,0030)	DA	2	
Patient's Sex	(0010,0040)	CS	2	

## 9.1.2.2.2 General Study Module

Attribute name	Tag	Vr	Type	Comment
Study Instance UID	(0020,000D)	UI	1	

## 9.1.2.2.3 RT Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	(0008,0060)	CS	1	Value not read
Series Instance UID	(0020,000E)	UI	1	

## 9.1.2.2.4 Structure Set Module

Attribute name	Tag	Vr	Type	Comment
Structure Set Label	(3006,0002)	SH	1	Value not read
Referenced Frame of Reference Sequence	(3006,0010)	SQ	3	
>Frame of Reference UID	(0020,0052)	UI	1	Value not read
>RT Referenced Study Sequence	(3006,0012)	SQ	3	
>>Referenced SOP Class UID	(0008,1150)	UI	1	Value not read
>>Referenced SOP Instance UID	(0008,1155)	UI	1	Value not read
>>RT Referenced Series Sequence	(3006,0014)	SQ	1	
>>>Series Instance UID	(0020,000E)	UI	1	
>>>Contour Image Sequence	(3006,0016)	SQ	1	
>>>>Referenced SOP Class UID	(0008,1150)	UI	1	
>>>>Referenced SOP Instance UID	(0008,1155)	UI	1	
Structure Set ROI Sequence	(3006,0020)	SQ	1	Value not read

## 9.1.2.2.5 ROI Contour Module

Attribute name	Tag	Vr	Type	Comment
ROI Contour Sequence	(3006,0039)	SQ	1	
>Referenced ROI Number	(3006,0084)	IS	1	
>Contour Sequence	(3006,0040)	SQ	3	
>>Contour Geometric Type	(3006,0042)	CS	1	Value not read
>>Number of Contour Points	(3006,0046)	IS	1	Value not read
>>Contour Data	(3006,0050)	DS	1	

## 9.1.2.2.6 RT ROI Observations Module

Attribute name	Tag	Vr	Type	Comment
RT ROI Observations Sequence	(3006,0080)	SQ	1	
>Observation Number	(3006,0082)	IS	1	Value not read
>Referenced ROI Number	(3006,0084)	IS	1	
>RT ROI Interpreted Type	(3006,00A4)	CS	2	
>ROI Physical Properties Sequence	(3006,00B0)	SQ	3	
>>ROI Physical Property	(3006,00B2)	CS	1	
>>ROI Physical Property Value	(3006,00B4)	DS	1	

## 9.1.2.2.7 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	(0008,0016)	UI	1	
SOP Instance UID	(0008,0018)	UI	1	

## 9.1.2.3 RT Plan IOD

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	Yes
Series	RT Series Module	Yes
Frame of Reference	Frame of Reference Module	Yes
Equipment	General Equipment Module	Yes
Plan	RT General Plan Module	Yes
	RT Patient Setup Module	Yes
	RT Fraction Scheme Module	Yes
	RT Beams Module	Yes
	RT Brachy Application Setups Module	No
	SOP Common Module	Yes

## 9.1.2.3.1 Patient Module

Attribute name	Tag	Vr	Type	Comment
Patient's Name	(0010,0010)	PN	2	
Patient ID	(0010,0020)	LO	2	
Patient's Birth Date	(0010,0030)	DA	2	
Patient's Sex	(0010,0040)	CS	2	

## 9.1.2.3.2 General Study Module

Attribute name	Tag	Vr	Type	Comment
Study Instance UID	(0020,0000)	UI	1	
Study Date	(0008,0020)	DA	2	
Study Time	(0008,0030)	TM	2	
Referring Physician's Name	(0008,0090)	PN	2	
Study ID	(0020,0010)	SH	2	
Accession Number	(0008,0050)	SH	2	
Study Description	(0008,1030)	LO	3	

## 9.1.2.3.3 RT Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	(0008,0060)	CS	1	Value not read
Series Instance UID	(0020,000E)	UI	1	

## 9.1.2.3.4 Frame of Reference Module

Attribute name	Tag	Vr	Type	Comment
Frame of Reference UID	(0020,0052)	UI	1	

## 9.1.2.3.5 General Equipment Module

Attribute name	Tag	Vr	Type	Comment
Software Versions	(0018,1020)	LO	3	

## 9.1.2.3.6 RT General Plan Module

Attribute name	Tag	Vr	Type	Comment
RT Plan Label	(300A,0002)	SH	1	
RT Plan Geometry	(300A,000C)	CS	1	Value not read

## 9.1.2.3.7 RT Patient Setup Module

Attribute name	Tag	Vr	Type	Comment
Patient Setup Sequence	(300A,0180)	SQ	1	
>Patient Setup Number	(300A,0182)	IS	1	
>Patient Position	(0018,5100)	CS	1C	
>Setup Technique	(300A,01B0)	CS	3	

## 9.1.2.3.8 RT Fraction Scheme Module

Attribute name	Tag	Vr	Type	Comment
Fraction Group Sequence	(300A,0070)	SQ	1	
>Fraction Group Number	(300A,0071)	IS	1	Value not read
>Number of Fractions Planned	(300A,0078)	IS	2	
>Number of Beams	(300A,0080)	IS	1	Value not read
>Number of Brachy Application Setups	(300A,00A0)	IS	1	Value not read

## 9.1.2.3.9 RT Beams Module

Attribute name	Tag	Vr	Type	Comment
Beam Sequence	(300A,00B0)	SQ	1	
>Beam Number	(300A,00C0)	IS	1	
>Beam Name	(300A,00C2)	LO	3	
>Beam Type	(300A,00C4)	CS	1	
>Radiation Type	(300A,00C6)	CS	2	
>Beam Limiting Device Sequence	(300A,00B6)	SQ	1	Value not read
>Treatment Delivery Type	(300A,00CE)	CS	3	
>Number of Wedges	(300A,00D0)	IS	1	
>Number of Compensators	(300A,00E0)	IS	1	
>Number of Boli	(300A,00ED)	IS	1	
>Number of Blocks	(300A,00F0)	IS	1	
>Number of Control Points	(300A,0110)	IS	1	Value not read
>Control Point Sequence	(300A,0111)	SQ	1	
>>Control Point Index	(300A,0112)	IS	1	
>>Nominal Beam Energy	(300A,0114)	DS	3	
>>Gantry Angle	(300A,011E)	DS	1C	
>>Gantry Rotation Direction	(300A,011F)	CS	1C	
>>Beam Limiting Device Angle	(300A,0120)	DS	1C	
>>Beam Limiting Device Rotation Direction	(300A,0121)	CS	1C	
>>Patient Support Angle	(300A,0122)	DS	1C	
>>Patient Support Rotation Direction	(300A,0123)	CS	1C	
>>Table Top Pitch Angle	(300A,0140)	FL	1C	
>>Table Top Pitch Rotation Direction	(300A,0142)	CS	1C	
>>Table Top Roll Angle	(300A,0144)	FL	1C	

>>Table Top Roll Rotation Direction	{300A,0146}	CS	1C	
>>Table Top Vertical Position	{300A,0128}	DS	2C	
>>Table Top Longitudinal Position	{300A,0129}	DS	2C	
>>Table Top Lateral Position	{300A,012A}	DS	2C	
>>Isocenter Position	{300A,012C}	DS	2C	

## 9.1.2.3.10 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	{0008,0016}	UI	1	
SOP Instance UID	{0008,0018}	UI	1	
RaySearch Private Creator	{4001,0010}	LO	3	RAYSEARCHLABS 2.0
RaySearch Checksum Algorithm Version	{4001,1060}	LO	3	RaySearch checksum algorithm version used to calculate the checksum of the plan.
RaySearch Checksum Data	{4001,1061}	OB	3	RaySearch custom checksum calculation specific for the current checksum algorithm version.

## 9.1.2.4 RT Beams Treatment Record IOD

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	Yes
Series	RT Series Module	Yes
Equipment	General Equipment Module	No
Treatment Record	RT General Treatment Record Module	Yes
	RT Patient Setup Module	Yes
	RT Treatment Machine Record Module	Yes
	RT Beams Session Record Module	Yes
	SOP Common Module	Yes
	Common Instance Reference Module	Yes

## 9.1.2.4.1 Patient Module

Attribute name	Tag	Vr	Type	Comment
Patient's Name	{0010,0010}	PN	2	
Patient ID	{0010,0020}	LO	2	
Patient's Birth Date	{0010,0030}	DA	2	
Patient's Sex	{0010,0040}	CS	2	

## 9.1.2.4.2 General Study Module

Attribute name	Tag	Vr	Type	Comment
Study Instance UID	{0020,000D}	UI	1	

## 9.1.2.4.3 RT Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	{0008,0060}	CS	1	Value not read
Series Instance UID	{0020,000E}	UI	1	

## 9.1.2.4.4 RT General Treatment Record Module

Attribute name	Tag	Vr	Type	Comment
Instance Number	{0020,0013}	IS	1	Value not read

Treatment Date	{3008,0250}	DA	2	
Treatment Time	{3008,0251}	TM	2	
Referenced RT Plan Sequence	{300C,0002}	SQ	2	
>Referenced SOP Class UID	{0008,1150}	UI	1	
>Referenced SOP Instance UID	{0008,1155}	UI	1	

## 9.1.2.4.5 RT Patient Setup Module

Attribute name	Tag	Vr	Type	Comment
Patient Setup Sequence	{300A,0180}	SQ	1	
>Patient Setup Number	{300A,0182}	IS	1	Value not read

## 9.1.2.4.6 RT Treatment Machine Record Module

Attribute name	Tag	Vr	Type	Comment
Treatment Machine Sequence	{300A,0206}	SQ	1	
>Treatment Machine Name	{300A,00B2}	SH	2	Used for validation purposes when reading treatment record from disk.

## 9.1.2.4.7 RT Beams Session Record Module

Attribute name	Tag	Vr	Type	Comment
Primary Dosimeter Unit	{300A,00B3}	CS	1	
Treatment Session Beam Sequence	{3008,0020}	SQ	1	
>Referenced Beam Number	{300C,0006}	IS	3	
>Beam Name	{300A,00C2}	LO	3	
>Beam Type	{300A,00C4}	CS	1	Value not read
>Radiation Type	{300A,00C6}	CS	1	Value not read
>Beam Limiting Device Leaf Pairs Sequence	{3008,00A0}	SQ	1	Value not read
>Number of Wedges	{300A,00D0}	IS	1	Value not read
>Current Fraction Number	{3008,0022}	IS	2	
>Treatment Delivery Type	{300A,00CE}	CS	2	
>Treatment Termination Status	{3008,002A}	CS	1	
>Specified Primary Meterset	{3008,0032}	DS	3	
>Specified Secondary Meterset	{3008,0033}	DS	3	
>Delivered Primary Meterset	{3008,0036}	DS	3	
>Delivered Secondary Meterset	{3008,0037}	DS	3	
>Specified Treatment Time	{3008,003A}	DS	3	
>Delivered Treatment Time	{3008,003B}	DS	3	
>Number of Control Points	{300A,0110}	IS	1	Value not read
>Control Point Delivery Sequence	{3008,0040}	SQ	1	
>>Treatment Control Point Date	{3008,0024}	DA	1	Value not read
>>Treatment Control Point Time	{3008,0025}	TM	1	Value not read
>>Delivered Meterset	{3008,0044}	DS	1	
>>Gantry Angle	{300A,011E}	DS	1C	
>>Patient Support Angle	{300A,0122}	DS	1C	
>>Table Top Pitch Angle	{300A,0140}	FL	1C	
>>Table Top Roll Angle	{300A,0144}	FL	1C	
>>Table Top Vertical Position	{300A,0128}	DS	2C	
>>Table Top Longitudinal Position	{300A,0129}	DS	2C	

>>Table Top Lateral Position	(300A,012A)	DS	2C	
>Tomo Registration	(300D,0010)	LO	3	TOMO_HA_01
>Tomo Registration Translation	(300D,10B0)	DS	3	
>Tomo Registration Rotation	(300D,10B1)	DS	3	
>Tomo bug Hi art 3.6.1	(300D,0010)	LO	3	LO
>Tomo Registration Translation Bug	(300D,10B0)	DS	3	
>Tomo Registration Rotation Bug	(300D,10B1)	SH	3	

## 9.1.2.4.8 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	(0008,0016)	UI	1	
SOP Instance UID	(0008,0018)	UI	1	

## 9.1.2.4.9 Common Instance Reference Module

Attribute name	Tag	Vr	Type	Comment
Referenced Series Sequence	(0008,1115)	SQ	1C	
>Series Instance UID	(0020,000E)	UI	1	Used to find the series of the Referenced RT Plan [300C,0002].
>Referenced Instance Sequence	(0008,114A)	SQ	1	
>>Referenced SOP Class UID	(0008,1150)	UI	1	Used to find the series of the Referenced RT Plan [300C,0002].
>>Referenced SOP Instance UID	(0008,1155)	UI	1	Used to find the series of the Referenced RT Plan [300C,0002].

## 9.1.2.5 Spatial Registration IOB

IE	Module	Used
Patient	Patient Module	Yes
Study	General Study Module	Yes
Series	General Series Module	Yes
	Spatial Registration Series Module	No
Frame of Reference	Frame of Reference Module	Yes
Equipment	General Equipment Module	No
Spatial Registration	Spatial Registration Module	Yes
	Common Instance Reference Module	Yes
	SOP Common Module	Yes

## 9.1.2.5.1 Patient Module

Attribute name	Tag	Vr	Type	Comment
Patient's Name	(0010,0010)	PN	2	
Patient ID	(0010,0020)	LO	2	
Patient's Birth Date	(0010,0030)	DA	2	
Patient's Sex	(0010,0040)	CS	2	

## 9.1.2.5.2 General Study Module

Attribute name	Tag	Vr	Type	Comment
Study Instance UID	(0020,000D)	UI	1	

## 9.1.2.5.3 General Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	(0008,0060)	CS	1	Value not read
Series Instance UID	(0020,000E)	UI	1	

Patient Position	(0018,5100)	CS	2C	
------------------	-------------	----	----	--

## 9.1.2.5.4 Frame of Reference Module

Attribute name	Tag	Vr	Type	Comment
Frame of Reference UID	(0020,0052)	UI	1	

## 9.1.2.5.5 Spatial Registration Module

Attribute name	Tag	Vr	Type	Comment
Content Date	(0008,0023)	DA	1	Value not read
Content Time	(0008,0033)	TM	1	Value not read
Instance Number	(0020,0013)	IS	1	Value not read
Content Label	(0070,0080)	CS	1	Value not read
Registration Sequence	(0070,0308)	SQ	1	
>Frame of Reference UID	(0020,0052)	UI	1C	
>Referenced Image Sequence	(0008,1140)	SQ	1C	
>>Referenced SOP Class UID	(0008,1150)	UI	1	
>>Referenced SOP Instance UID	(0008,1155)	UI	1	
>Matrix Registration Sequence	(0070,0309)	SQ	1	
>>Matrix Sequence	(0070,030A)	SQ	1	
>>>Frame of Reference Transformation Matrix	(3006,00C6)	DS	1	
>>>Frame of Reference Transformation Matrix Type	(0070,030C)	CS	1	Value not read

## 9.1.2.5.6 Common Instance Reference Module

Attribute name	Tag	Vr	Type	Comment
Referenced Series Sequence	(0008,1115)	SQ	1C	
>Series Instance UID	(0020,000E)	UI	1	Used to match RT Struct and Image against Spatial registration
>Referenced Instance Sequence	(0008,114A)	SQ	1	
>>Referenced SOP Class UID	(0008,1150)	UI	1	
>>Referenced SOP Instance UID	(0008,1155)	UI	1	

## 9.1.2.5.7 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	(0008,0016)	UI	1	
SOP Instance UID	(0008,0018)	UI	1	

## 9.1.2.6 RT Beams Delivery Instruction RETIRED IOD

IE	Module	Used
Patient	Patient Module	No
	Clinical Trial Subject Module	No
Study	General Study Module	Yes
	Patient Study Module	No
	Clinical Trial Study Module	No
Series	General Series Module	Yes
	Clinical Trial Series Module	No
Equipment	General Equipment Module	No
Plan	RT Beams Delivery Instruction Module	Yes



	Common Instance Reference Module	No
	General Reference Module	No
	SOP Common Module	Yes

## 9.1.2.6.1 General Study Module

Attribute name	Tag	Vr	Type	Comment
Study Instance UID	(0020,000D)	UI	1	

## 9.1.2.6.2 General Series Module

Attribute name	Tag	Vr	Type	Comment
Modality	(0008,0060)	CS	1	Value not read
Series Instance UID	(0020,000E)	UI	1	

## 9.1.2.6.3 RT Beams Delivery Instruction Module

Attribute name	Tag	Vr	Type	Comment
Beam Task Sequence	(0074,1020)	SQ	1	
>Beam Task Type	(0074,1022)	CS	1	Value not read
>Treatment Delivery Type	(300A,00CE)	CS	1	Supported values: TREATMENT, CONTINUATION.
>Continuation Start Meterset	(0074,0120)	FD	1C	
>Continuation End Meterset	(0074,0121)	FD	1C	
>Current Fraction Number	(3008,0022)	IS	1	Value not read
>Referenced Beam Number	(300C,0006)	IS	1	
Omitted Beam Task Sequence	(300C,0111)	SQ	3	
>Referenced Beam Number	(300C,0006)	IS	1	
>Reason for Omission	(300C,0112)	CS	1	Value not read

## 9.1.2.6.4 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	(0008,0016)	UI	1	
SOP Instance UID	(0008,0018)	UI	1	

## 9.1.2.7 Unified Procedure Step RETIRED IOD

IE	Module	Used
Unified Procedure Step	SOP Common Module	Yes
	Unified Procedure Step Relationship Module	Yes
	Unified Procedure Step Scheduled Procedure Information Module	Yes
	Unified Procedure Step Progress Information Module	Yes
	Unified Procedure Step Performed Procedure Information Module	No
	Patient Demographic Module	No
	Patient Medical Module	No
	Visit Identification Module	No
	Visit Status Module	No
	Visit Admission Module	No
	Transaction Module	Yes

## 9.1.2.7.1 SOP Common Module

Attribute name	Tag	Vr	Type	Comment
SOP Class UID	(0008,0016)	UI	1	Value not read
SOP Instance UID	(0008,0018)	UI	1	
Timezone Offset From UTC	(0008,0201)	SH	3	

## 9.1.2.7.2 Unified Procedure Step Relationship Module

Attribute name	Tag	Vr	Type	Comment
Patient's Name	(0010,0010)	PN		
Patient ID	(0010,0020)	LO		

## 9.1.2.7.3 Unified Procedure Step Scheduled Procedure Information Module

Attribute name	Tag	Vr	Type	Comment
Scheduled Station Name Code Sequence	(0040,4025)	SQ		
>Code Value	(0008,0100)	SH	1C	
>Coding Scheme Designator	(0008,0102)	SH	1C	
>Code Meaning	(0008,0104)	LO	1	Value not read
Scheduled Procedure Step Start DateTime	(0040,4005)	DT		The date used to filter the UPSs to return in the worklist query. Supports either a range in the format 'yyyyMMddHHmmss-yyyyMMddHHmmss' or a single date in the format 'yyyyMMddHHmmss'.
Scheduled Workitem Code Sequence	(0040,4018)	SQ		
>Code Value	(0008,0100)	SH	1C	
>Coding Scheme Designator	(0008,0102)	SH	1C	
>Code Meaning	(0008,0104)	LO	1	
Study Instance UID	(0020,0000)	UI		
Input Availability Flag	(0040,4020)	CS	1	Value not read

## 9.1.2.7.4 Unified Procedure Step Progress Information Module

Attribute name	Tag	Vr	Type	Comment
Procedure Step State	(0074,1000)	CS		
Procedure Step Progress Information Sequence	(0074,1002)	SQ		
>Procedure Step Progress	(0074,1004)	DS		
>Procedure Step Progress Description	(0074,1006)	ST		

## 9.1.2.7.5 Transaction Module

Attribute name	Tag	Vr	Type	Comment
Transaction UID	(0008,1195)	UI	3	

## 9.1.3 Attribute Mapping

Not applicable

## 9.1.4 Coerced/Modified Fields

Not applicable

## 9.2 DATA DICTIONARY OF PRIVATE ATTRIBUTES

All used Private Creators are listed in the table below. Usage of Private Attributes are listed in each module specification.

Attribute name	Tag	VR	VM	Value
RaySearch Private Creator	(4001,0010)	LO	1	RAYSEARCHLABS 2.0
Tomo bug Hi art 3.6.1	(300D,0010)	LO	1	LO

Tomo Registration	(300D,0010)	L0	1	TOMO_HA_01
-------------------	-------------	----	---	------------

9.3 CODE TERMINOLOGY AND TEMPLATES

Not applicable

9.4 GRAYSCALE IMAGE CONSISTENCY

Not applicable

9.5 STANDARD EXTENDED/SPECIALIZED/PRIVATE SOP CLASSES

9.5.1 Standard extended SOP Class

9.5.2 Specialized SOP Class

Not applicable

9.5.3 Private SOP Class

Not applicable

9.6 PRIVATE TRANSFER SYNTAXES

Not applicable





## CONTACT INFORMATION



**RaySearch Laboratories AB (publ)**  
Sveavägen 44  
SE-111 34 Stockholm  
Sweden

### Contact details head office

P.O. Box 3297  
SE-103 65 Stockholm, Sweden  
Phone: +46 8 510 530 00  
Fax: +46 8 510 530 30  
info@raysearchlabs.com  
www.raysearchlabs.com

### RaySearch Americas

Phone: +1 877 778 3849

### RaySearch France

Phone: +33 1 76 53 72 02

### RaySearch Korea

Phone: +82 10 2230 2046

### RaySearch Australia

Phone: +61 411 534 316

### RaySearch Belgium

Phone: +32 475 36 80 07

### RaySearch Germany

Phone: +49 30 893 606 90

### RaySearch Singapore

Phone: +65 81 28 59 80

### RaySearch China

Phone: +86 137 0111 5932

### RaySearch Japan

Phone: +81 3 44 05 69 02

### RaySearch UK

Phone: +44 2039 076791